


User Manual

H006A

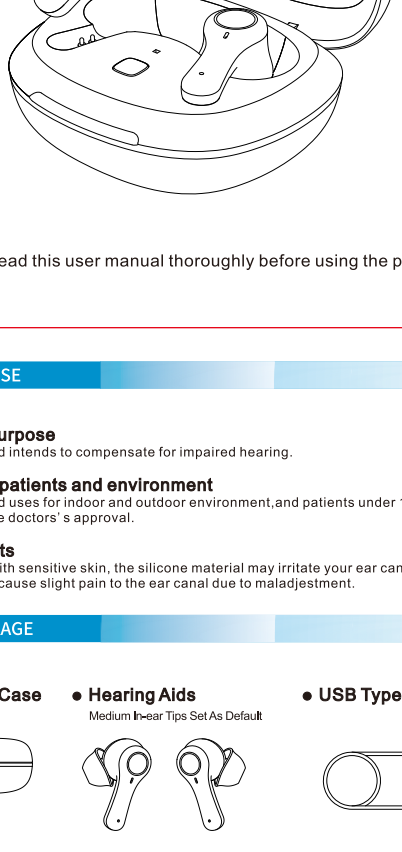

Please read this user manual thoroughly before using the product.

INTENDED USE

1. Medical purpose

-The hearing aid intends to compensate for impaired hearing.

2. Intended patients and environment


-The hearing aid uses for indoor and outdoor environment, and patients under 18 years old must require the doctors' approval.

3. Side effects

- For patients with sensitive skin, the silicone material may irritate your ear canals .
- First use may cause slight pain to the ear canal due to maladjustment.


IN THE PACKAGE

● Charging Case

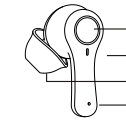


● Hearing Aids


Medium In-Ear Tips Set As Default

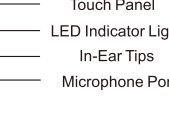


● USB Type-C Cable

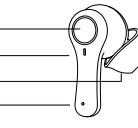

● In-Ear Tips

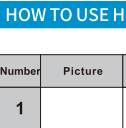


● User Manual

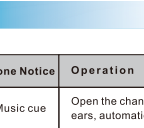

ASSEMBLE THE EARPLUGS


Please choose the appropriate size earplugs, and assemble the ear tips as below picture.

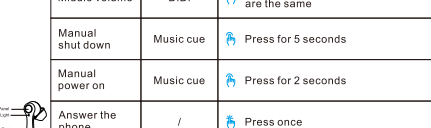


In-Ear tips

HOW TO WEAR YOUR HEARING AID



-Hold the hearing aid with your index finger and thumb, with the little pore facing up .
-Gently insert the hearing aid into the ear canal, and then rotate it backwards until the product is stable in the ear and will not fall . Insert the Left and Right hearing aid into your ear canals as shown.

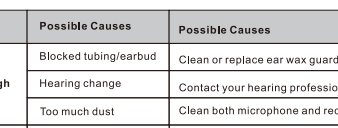

KNOW YOUR HEARING AIDS

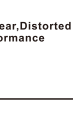

Left Hearing Aid

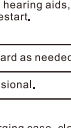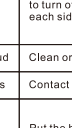

Right Hearing Aid

HOW TO USE HEARING AIDS

Number	Picture	Function	Tone Notice	Operation
1		Turn on	Music cue	Open the charging case, wear the hearing aids on your ears, automatically enter hearing aid mode after 10 seconds
2		Turn off	/	Put the hearing aids back to charging case, close the case cover, then Hearing aids automatically turn off
3		Volume adjustment	Music cue	Press once a total of seven volume levels, which will be increased and cyclically adjusted by the volume buttons
4		Minimum volume	DI	Press once
5		Maximum volume	DiDiDi	Press once
6		Middle volume	DiDi	Press once the beeps of the middle five gears are the same
7		Manual shut down	Music cue	Press for 5 seconds
8		Manual power on	Music cue	Press for 2 seconds
9		Answer the phone	/	Press once
10		Hang up the phone	/	Press for 2 seconds
11		Reject the call	/	Press for 2 seconds
12		Low battery in charging compartment		Blinking blue LED light
13		Hearing aid fully charged		Red LED light is off
14		Factory reset the charging case		Press and hold the reset button of the charging compartment for 5 seconds

Notice: If the battery is too low after use, put the hearing aid into the charging case to charging

TROUBLESHOOTING GUIDE

Symptom	Possible Causes	Possible Causes
Not Loud Enough	Blocked tubing/earbud	Clean or replace ear wax guard as needed.
	Hearing change	Contact your hearing professional.
	Too much dust	Clean both microphone and receiver with brush.
Inconsistent Performance	Restart required	Put the hearing aid into charging case, close the top cover to turn off. Then take out the hearing aids, press and hold each side for 2 seconds to restart.
	Blocked tubing/earbud	Clean or replace ear wax guard as needed.
Unclear, Distorted Performance	Restart required	Put the hearing aid into charging case, close the top cover to turn off. Then take out the hearing aids, press and hold each side for 2 seconds to restart.
	Blocked tubing/ earbud compartment	Clean or replace ear wax guard as needed.
Dead Hearing aids	Defective hearing aids	Contact your hearing professional.
	Charge required	Put the hearing aid into charging case, close the top cover and wait for charging.
	Blocked earbud/ tubing	Clean or replace ear wax guard as needed.
No LED Indicator when hearing aids are inserted into the charging port	Ear canal filled with oil wax	Contact your hearing professional.
	Low battery	Connect the Type-c cable to your charging case and plug it into the wall. The on-board battery LED's will flash for a few seconds to indicate the connection to the power source.
Blinking blue LED light	Low battery	Insert your type-c cord to charge the case
Hearing aids whistling in the charging case	Dead battery	Put the hearing aid into charging case, close the top cover, and insert your type -c cord to charge the case.

FAULT SELF-CHECK

Project	Possible Causes	Solution
Hearing aid has no sound	The button time is too short to start	Press and hold touch panel for 3 seconds to restart
	Low battery	Put the hearing aid into the charging case to charging
	Blocked up the earplugs' sound port or microphone port	Clean the sound port or microphone port
The voices is very weak and unclear	Low battery	Put the hearing aid into the charging case to charging
	The volume is too low	Double click to adjust the volume level
	Earplugs loose	Re-wearing the hearing aid

Hearing aid failure self-checking table

TECHNICAL SPECIFICATION

Model	H006A
Output, OSPL90 Max.	≤105dB
FOG50 MAX Gain	≤33dB
EQ Input MAX	±33dB
Frequency Range(Not narrower than)	200Hz-4500Hz W
Working Current	≤8mA/
Total Harmonic Distortion	<5%
Working Voltage	D.C 3.7V
Play Time	8hours

Notice: The above parameter index testing is on basis of FONIX8000Hearing Aid Test System; Inspection standard: IEC 60118-7 2005, 20CC coupling cavity HA-1

HOW TO CONNECT H006A HEARING AID USING BLUETOOTH ON

IPHONE AND ANDROID PHONE

Step 1: Turn on the phone's Bluetooth search

Step 2: Then open the charging box of the H006A hearing aid (do not take out the hearing aid), wait for a few seconds, and the search result will automatically pop up in the Bluetooth search of your phone

Step 3: You will see the hearing aid name "Sengdong" or "coniler" in the Bluetooth search, click to complete the connection

Step 4: After the connection is completed, you can use the phone volume keys to control the volume of the hearing aid. Of course, if you use the phone to listen to music or watch movies, it will be conveyed to your ears through the hearing aid

PS: The Bluetooth function only needs to be set once. After that, every time you use the hearing aid, it will automatically connect to the Bluetooth of the mobile phone. Please turn on the Bluetooth function of the phone when using the hearing aid

HOW TO USE THE APP (IPHONE USER)


Step 1: Open the APP Store, click the search box and enter "coniler", download and install the app (can also scan the QR code below to download the app)

Step 2: First, please turn on the Bluetooth function of your iphone, then open the "coniler" APP, click "Search device", and then open the charging box of your H006A hearing aid (please do not take out the hearing aid), then your APP will Automatically find hearing aids

Step 3: After waiting for a few seconds, you will see the name of the hearing aid "Sengdong Right" and "Sengdong Left" displayed in the APP, click on them to complete the connection (if you only need to use a hearing aid for one ear, only Click on the hearing aid for the ear you want to use)

Step 4: After completing the connection, you will be able to use the APP to control the volume, Frequency band and usage mode of the hearing aid

PS: There is no need to use the APP to control the hearing aid every time. In order to improve the battery life of the hearing aid, we have canceled the setting of the automatic connection to the APP. Therefore, every time you use the APP, you need to follow the above steps to connect the hearing aid. When connecting the APP, please pay attention to whether the hearing aid has power (blue light flashes on the earphone), please charge first if there is no power

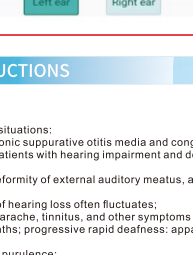

HOW TO USE THE APP (ANDROID USER)

Step 1: Please be sure to use your mobile phone to log in to the official website (<http://www.srmcm.com/Coniler/>). After into the website, click "Official" to automatically download the application.

Step 2: First, please turn on the Bluetooth function of your iphone, then open the "coniler" APP, click "Search device", and then open the charging box of your H006A hearing aid (please do not take out the hearing aid), then your APP will Automatically find hearing aids

Step 3: After waiting for a few seconds, you will see the name of the hearing aid "Sengdong Right" and "Sengdong Left" displayed in the APP, click on them to complete the connection (if you only need to use a hearing aid for one ear, only Click on the hearing aid for the ear you want to use)



Step 4: After completing the connection, you will be able to use the APP to control the volume, Frequency band and usage mode of the hearing aid

PS: There is no need to use the APP to control the hearing aid every time. In order to improve the battery life of the hearing aid, we have canceled the setting of the automatic connection to the APP. Therefore, every time you use the APP, you need to follow the above steps to connect the hearing aid. When connecting the APP, please pay attention to whether the hearing aid has power (blue light flashes on the earphone), please charge first if there is no power

HOW TO USE THE APP

1 Adjust the volume.

Move the volume slider, you will get different volume, select comfortable volume.

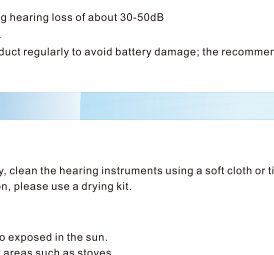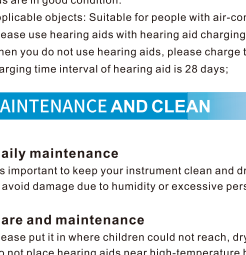

*These programs are preset according to the needs of most users in common scenarios. We suggest that you try each program and choose the most appropriate one.

*If your left and right ears are connected, changing program in the app will change the two devices at the same time.


HOW TO USE THE APP

2 Adjust the program.

Press the square for changing the program. Select the mode close to the audiometric curve of yourself according to the icon or Select the comfortable program for you. If you don't have an audiogram, you can choose the option with the best sound quality

3 Adjust the Equalizer band.

Move the slider for changing the gain for different frequency, select the comfortable sound.

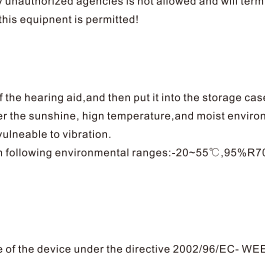

IMPORTANT SAFE INSTRUCTIONS

Δ 1.CONTRAINDICATION

The device is unsuitable in the following situations:
-Usage is prohibited for patients with chronic suppurative otitis media and congenital ear atresia.
-This hearing aid can only compensate patients with hearing impairment and do not have hearing therapeutic value.
1) Congenital atresia of external ear or deformity of external auditory meatus, a large number of foreign bodies in external acoustic meatus;
2) Deafness is unstable, and the degree of hearing loss often fluctuates;
3) Complaint with headache, dizziness, earache, tinnitus, and other symptoms that are not suitable for adaptation;
4) Sudden deafness in the past three months; progressive rapid deafness; apparent unilateral deafness or other acute ear diseases;
5) Effusion in the middle ear and frequent purulence;
6) Central deafness; non-organic deafness;
7) Patients with acute otitis externa and tympanitis;
8) Patients with acute and chronic suppurative otitis media (in the period of purulent infection);
9) Patients who are allergic to the material of this product.

Δ 2.WARNINGS

-Keep small swallow-able accessories out of children's reach.
-Dropping, immersing in liquid, and accessing strong electromagnetic fields or over heat will damage hearing aids.
-If a hearing aid or accessories remain in the ear, please go to see a doctor as soon as possible, do not take it out by yourself.
-Warning not to use in explosive or oxygen-enriched atmospheres; it may damage the device.
-Warning that only the intended person can use it; using a hearing aid designed for someone else may result in hearing damage.
-Statement about common conditions that could damage the device
-Warning about common special needs of particular patient groups
-Warning that results from the risk management file permissible environmental conditions of transporting storage.
-Each warning/safety contains: nature of Hazard, likely consequences, precautions

Δ 3.NOTICE

-Please Do not use it at high volume for a long time. Otherwise, it will injure your ear. If you feel uncomfortable, please stop using it.
-Please check the product every time before using it, and avoid wearing the earplug with foreign objects in ear canal. It will affect the sound amplifier.
-Avoid shocking and dropping the product. Do not keep or use it in high humidity environments (such as the sauna room).
-Do not put the scrapped battery as standard rubbish.
-Do not use the product in a potential explosive or dangerous environment .
-Do not share your hearing aids with others. Your hearing aids setting is personal and individual, share with someone else will result in others' hearing damage.
-Be careful with any unusual skin reactions. In some instances, wearing a hearing aid can cause an allergic reaction leading to it. Chiness or rashes. If you notice any uncommon skin condition, stop wearing the hearing aids, and check with a medical specialist.
-Do not place the hearing aids near the magnetic objects.
-Be sure to take off the hearing aids when getting a diagnosis or treatment.
-Never pry the inside of the hearing aids with a needle or similar object; otherwise, it will damage the device.
-Do not put the hearing aid near heating appliances or in direct sunlight, and you should also avoid high moisture areas such as bathrooms or near a humidifier.
-Warning that results from the risk management file permissible environmental conditions of transporting storage.
-Do not use the hearing aids close to a cell phone, otherwise it may increase the noise and affect the sound quality.
-When you do not use hearing aids, please charge the product regularly to avoid battery damage; the recommended charging time interval of hearing aid is 28 days;

SERVICE AND REPAIR

1. If for any reason, your hearing aids and charger do not operate properly, do NOT attempt to fix it yourself.
Not only are you likely to violate any applicable warranties or insurance, you may also cause further damage.
2. If your hearing aids and charger fail or perform poorly, check the guide on the next page for possible solutions.
If problems continue, contact your hearing professional for advice and assistance. Most of the problems can find in your hearing professional's office or clinic.

MATTERS NEEDING ATTENTION, WARNINGS, AND TIPS

1. Before fitting the hearing aid, it should go through professional inspection and hearing test, and use it under the guidance of a doctor or an experienced fitting technician of hearing aid.
2. Do not immerse the hearing aid in any liquid. Remove the hearing aid before swimming and shampooing.
After swimming or washing, wait for the ear canal to dry before wearing the hearing aid.
3. Please use the function and volume control switch directly, open and close as lightly as possible to avoid excessive force
4. Please send the hearing aids regularly to the fitting center for maintenance and inspection to ensure that the hearing aids are in good condition.
5. Applicable objects: Suitable for people with air-conducting hearing loss of about 30-50dB
6. Please use hearing aids with hearing aid charging case:
7. When you do not use hearing aids, please charge the product regularly to avoid battery damage; the recommended charging time interval of hearing aid is 28 days;

MAINTENANCE AND CLEAN

1.Daily maintenance

It is important to keep your instrument clean and dry. Daily, clean the hearing instruments using a soft cloth or tissue. To avoid damage due to humidity or excessive perspiration, please use a drying kit.

2.Care and maintenance

Please put it in where children could not reach, dry and no exposed in the sun.
Do not place hearing aids near high-temperature heating areas such as stoves.
Do not leave the hearing aid in a humid place for a long time.
Do not scratch the microphone with sharp objects to avoid damage to the microphone
Do not disassemble or attempt any repairs by yourself. Ask the dispenser for repair.

3.Cleaning

1) Cleaning the earplug
Remove the wax plugged in the earpiece by using a stick and brush.
2) Cleaning hook
Use a soft cloth or tissue, wipe of any moisture, sweat, or debris from the surface of the hook .
Note: Do not use alcohol or other cleaning solvents to clean the hook as this could damage the protective covering.

HELPFUL HINTS

-Do not take apart your hearing aids or insert the cleaning tools inside them.
(When not in use, place your hearing aids in the charger and store them.) - In a dry, safe place .
-Away from direct sunlight or heat to avoid extreme temperatures.
-Where you can easily find them.
-Safely out of reach from children and pets. Synergy Charger Care
-Keep your charger clean. Heat, moisture, and bright materials can result in poor performance.
-Using the provided cleaning brush, keep the charging ports clear of debris.
-Do not use water, solvents, or cleaning fluids to clean the charging ports
-Keep lid closed as much as possible to avoid dust and debris build- up.
-Store your charger in a clean and dry location
-Fully charge hearing aid batteries every night.

MAINTENANCE STORAGE AND DISPOSAL

1.Maintenance

1) We do not authorize any institution or individual to maintain and repair the product. If you suspect that the products have any issues please contact the manufacturer or distributor to handle them.
2) The user must not attempt any repairs to the device or any of its accessories. Please contact the retailer for repair.
3) Opening of the equipment by unauthorized agencies is not allowed and will terminate any warranty claim.
Warning: No modification of this equipment is permitted!

2.Storage

1) If you Do not use it, turn off the hearing aid, and then put it into the storage case
2) Do not put the hearing aid under the sunshine, high temperature, and moist environment or someplace, which may get in touch with fire or is vulnerable to vibration.
3) Keep the device storage in following environmental ranges: -20~55℃, 95%RH/00-1060hPa

3.Disposal

Please dispose of the device under the directive 2002/96/EC- WEEE (Waste Electrical and Electronic Equipment) and Don't dispose with domestic waste. If you have any queries, please refer to the local authorities responsible for waste disposal.

HEARING AIDS CARE

Keep your hearing aids clean. Heat, moisture, and foreign substances can result in poor performance.
-Use a cleaning brush or soft cloth to clean debris from around the user control and microphone; inspect the receiver, earbud, or ear tip and wax guard regularly.
-Never use water, solvents, cleaning fluids, or oil to clean your hearing aids.
-Charging contacts should always be kept clean and dry. Your hearing professional can provide further information on additional maintenance procedures for your hearing aids, if needed

IMPORTANT NOTICE FOR PROSPECTIVE HEARING AID USERS

Good health practice requires that a person with a hearing loss have a medical evaluation by a licensed physician (preferably a specialist who specializes in diseases of the ear) before purchasing a hearing aid. Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists or otorhinolaryngologists. The purpose of the medical evaluation is to assure that all medically treatable conditions that may affect hearing are identified and treated before the hearing aid is purchased .

Warning to Hearing Aid Dispensers.
A hearing aid dispenser should advise a prospective hearing aid user to consult promptly with a licensed physician (preferably an ear specialist) before dispensing a hearing aid if the hearing aid dispenser determines through inquiry, actual observation, or review of any other available information concerning the prospective user, that the prospective user has any of the following conditions:

(i) Visible congenital or traumatic deformity of the ear.
(ii) History of active drainage from the ear within the previous 90 days.
(iii) History of sudden or rapidly progressive hearing loss within the previous 90 days.
(iv) Acute or chronic dizziness.
(v) Unilateral hearing loss of sudden or recent onset within the previous 90 days.
(vi) Audiometric air - bone gap equal to or greater than 15 decibels at 500 hertz (Hz), 1,000 Hz, and 2,000 Hz
(vii) Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.
(viii) Pain or discomfort in the ear special care should be exercised in selecting and fitting a hearing aid whose maximum sound pressure level exceeds 132 decibels because there may be risk of impairing the remaining hearing of the hearing aid user.
(This provision is required only for those hearing aids with a maximum sound pressure capability greater than 132 decibels (dB).)

NORMALIZED SYMBOLS

	Consult instructions for use
	B type applied part
	3 Disposal Please dispose of the device under the directive 2002/96/EC- WEEE (Waste Electrical and Electronic Equipment) and Don't dispose with domestic waste.
	Manufacturer
	Date of manufacture
	Serial number
	Warning/Danger: Improper use might cause serious injuries.
	IP41 IP code of the device: this device's grade of against ingress of solid foreign objects - the actual observation, or review of any other available information concerning the prospective user, that the prospective user has any of the following conditions:

AFTER SALES & FEEDBACK

Please contact us via the message system of the marketplace if you have any questions about the usage, trouble shooting and feedback/comment of our products.
Company Name: Shenzhen Yazhouyafang Management Co., Ltd.
Email Address: coniler20@gmail.com

Enjoy life's conversation