



VISCOJECT-BIO 2.2 Injector Set

Handling instructions for CT SPHERIS 204,
CT ASPHINA 404 and high diopter ZEISS toric IOLs




VISCOJECT-BIO 2.2 Injector Set

Preparing the compatible ZEISS IOLs for use.



Place a small amount of OVD* in the cartridge tip.
Place a line of OVD in each channel of the cartridge.
Do not touch the tip of the cartridge.

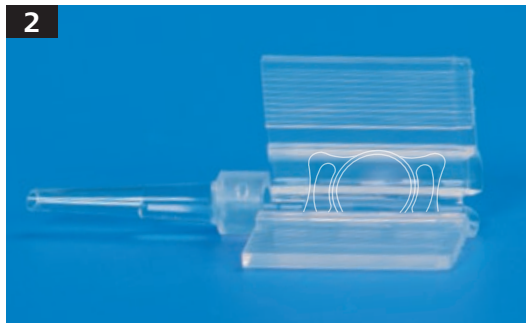

Place the lens in the center of the open cartridge. Push down on the edges of the optic to ensure that the optic haptics are positioned under the ledges of the cartridge.



For toric and multifocal IOLs, mind the orientation holes and the meridian marks.



Maintain the IOL in place with the folding forceps while folding the cartridge. Close the cartridge until you hear a click.



Place the cartridge into the injector and move it forward as far as it will go. The lens is ready for implantation.



Place the bevelled edge downward and gently push the plunger to inject the IOL.

**Carl Zeiss Meditec recommends the use of ZEISS OVDs.*

EN_32_022_001811 Printed in Germany. CZ-X/2018 International edition: Only for sale in selected countries.
The contents of the brochure may differ from the current status of approval of the product or service offering in your country. Please contact our regional representatives for more information. Subject to changes in design and scope of delivery and due to ongoing technical development. VISCOJECT is a trademark of Medcel (CE 0482).
© Carl Zeiss Meditec AG, 2018. All rights reserved.



Carl Zeiss Meditec AG
Goeschwitzer Str. 51–52
07745 Jena
Germany
www.meditec.zeiss.com/iol
www.meditec.zeiss.com/contacts



0297

CT ASPHINA 404
CT SPHERIS 204