
## EN Chlamydia Test

The cerascreen® Chlamydia Test is a send-in sampling kit for the determination of Chlamydia trachomatis and Neisseria gonorrhoeae. The test is used to detect the bacteria Chlamydia trachomatis and Neisseria gonorrhoeae in your native urine, which can transmit sexually transmitted diseases (chlamydia and gonorrhoea). The sample is analyzed in vitro in our specialised cerascreen® partner diagnostics laboratory.

## Your test kit contains

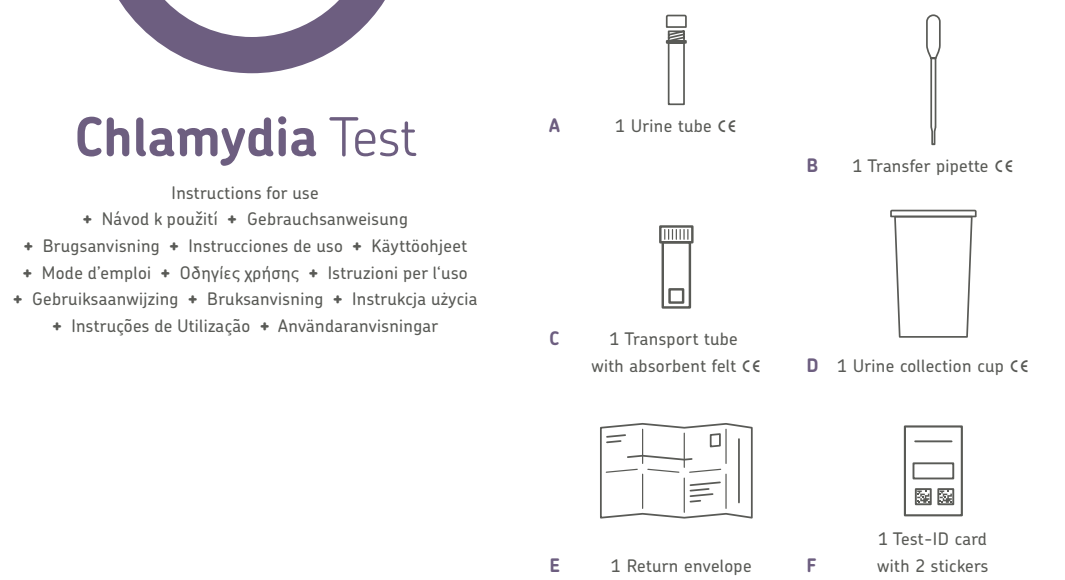
## Information and tips

 Read the Instructions for use carefully before you begin and follow the listed steps.

**Manufacturer:**  
cerascreen GmbH  
Güterbahnhofstrasse 16  
19059 Schwerin  
Germany

**SKU:** 99090856  
CRS3-Chlam\_INT\_000  
Valid from April 2021

**CE** **Supplier | Dodavatel | Lieferant | Leverander | Proveedor | Toimittaja | Fournisseur | Προμηθευτής | Fornitore | Leverancier | Leverador | Dostawca | Fornecedor | Leverantör:**  
**A, B, C, D:** Sarstedt AG & Co.KG, Sarstedtstraße 1,51588 Nümbrecht, Germany
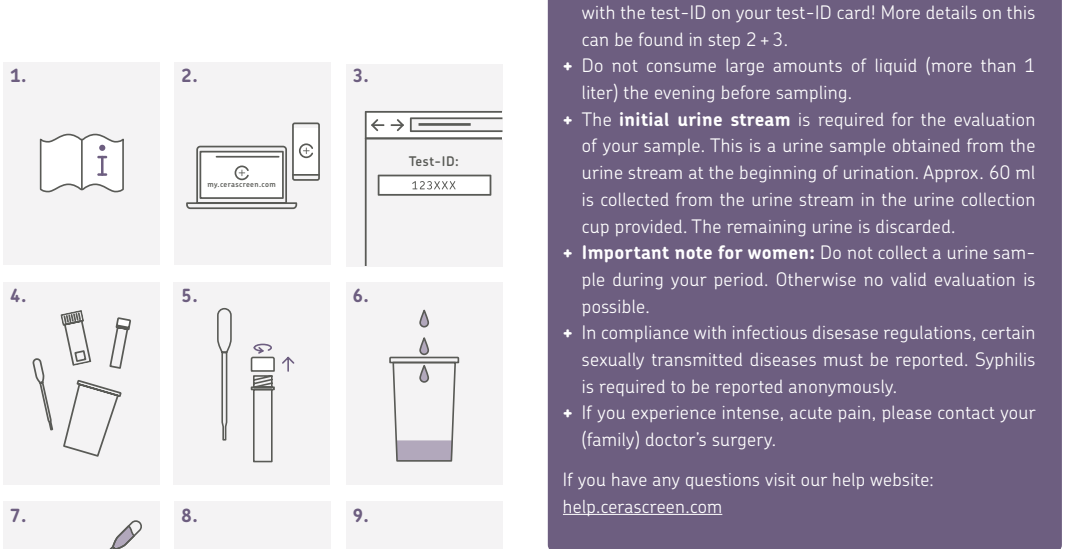
## Directions for use

- Check that all conditions are fulfilled for sample taking. Please read the section "Information and tips" carefully.
- If you are a new customer, please create an account at [my.cerascreen.com](https://my.cerascreen.com) or with the **my cerascreen®** app. See our FAQs on [help.cerascreen.com](https://help.cerascreen.com) for further instructions.
- To receive your result report, log in via the **my cerascreen®** app (Android, iOS) or on [my.cerascreen.com](https://my.cerascreen.com). Insert your six-digit ID-number according to your test-ID card under "activate test" and follow the instructions.
- Place all test kit items on a clean surface and check that none of the components are missing.
- Open the urine tube and place it on a flat surface. Make sure that it cannot fall over. Lay the transfer pipette next to it.
- Take the urine collection cup and collect an initial urine sample in it (see "IMPORTANT").
- Use the transfer pipette to transfer your urine into the smaller urine tube. Fill the tube completely.
- Close the urine tube by screwing on the cap.
- Remove the two stickers from the test ID card. Stick one sticker on the filled test tube and the other sticker on the transport tube.
- Slide the filled collection tube into the larger transport tube without removing the absorbent felt. Seal the transport tube well.
- Place the transport tube containing the urine collection tube into the return envelope. Close the envelope carefully and send it immediately by post free of charge.
- Check that you have activated the test and dispose the used components in the clear plastic bag in the household waste.

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