


H5020 5 IN 1 Multi-functional beauty instrument




This product is not a medical equipment, Pls do not used as disease treatment!

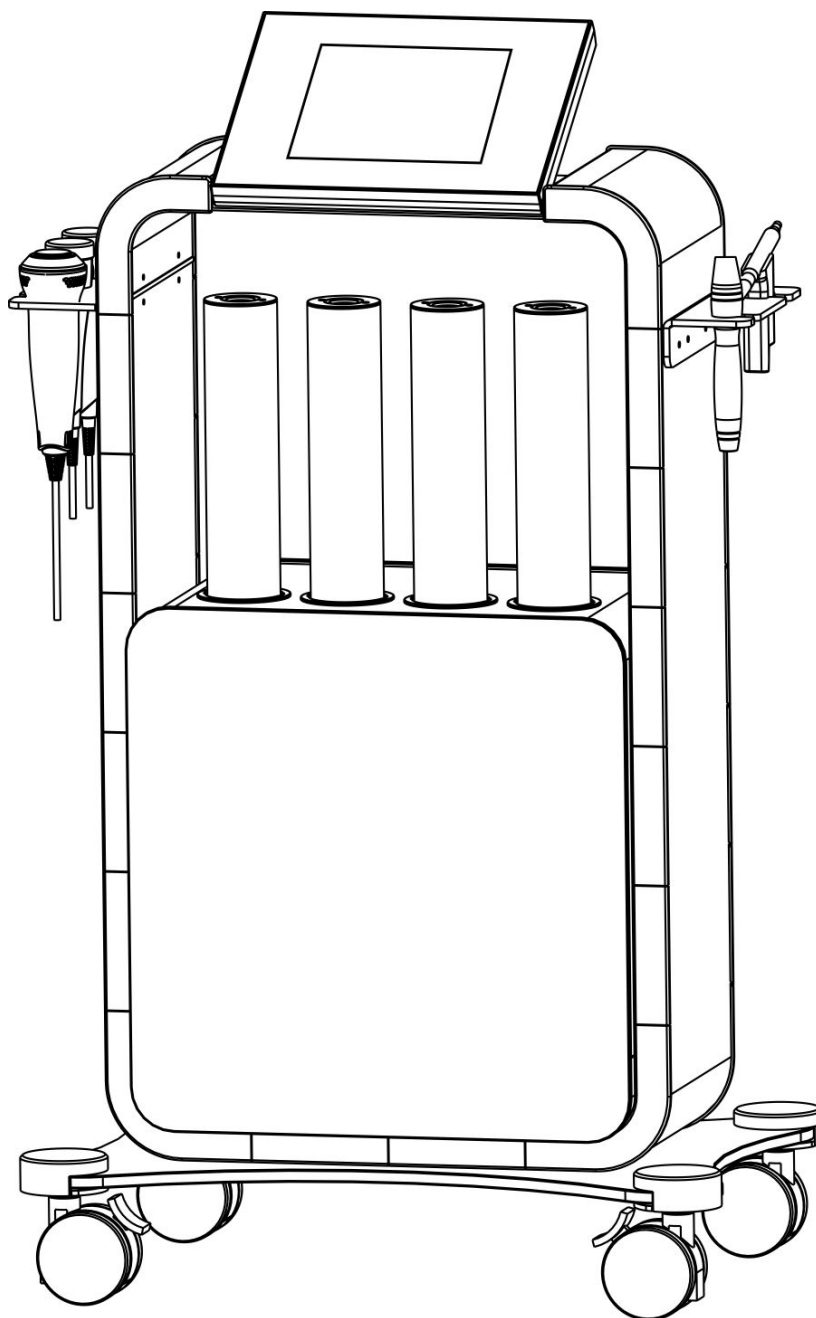

Specifications and Technical Data:

Product Name	5 IN 1 Multi-functional beauty instrument
Item No.	H5020
Function	Water& oxygen jet facial clean、Ultrasonic、RF Lifting、Spray Water Oxygen、Cool and Hot Hammer
Rated Voltage/Power	AC 100-120V~
	AC 220-240V~
Working Frequency	50/60Hz
Rated Power	190W
Dimension(LXWXH)	61*41*90.1cm
N.W.	42Kg


Packing List:

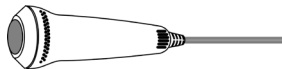Item	Product Description	QTY
Host Machine		1 set
Power cable		1 pc
RF handle		1 pc
Ultrasonic handle		1 pc
Cool and hot hammer handle		1 pc
Spray Water Oxygen handle		1 pc
Water & oxygen jet facial clean handle		1 pc
Grinding head		1 set
Hard plastic head		1 set
Silicone head		1 set
Instruction Manual		1copy



Host Machine Instruction


1. Display screen
2. Accessories
3. Accessories holder
4. Bottle(A,B,C,D)
5. Machine frame
6. Front panel
7. Power switch
8. Fuse holder
9. Power cord socket
10. Back panel
11. Waste water bottle

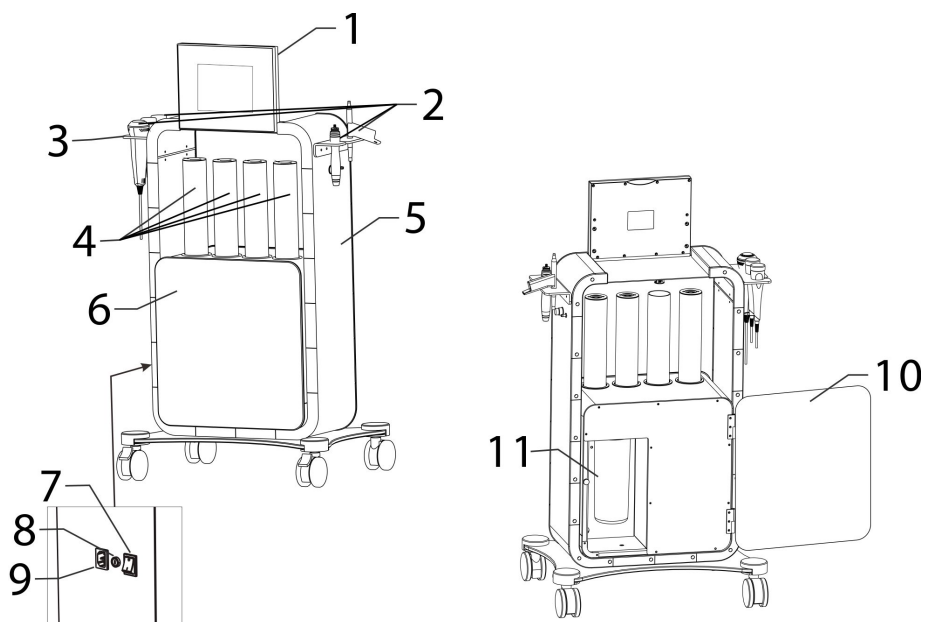

Instruction and Attentions:

1. Follow the labels to install the accessories.
2. Fill the mineral water or primary liquid in the bottle, Then open the top button.
3. Turn on the power switch on the side of the main frame. The power switch indicator light is on.
4. Lift the display screen gently
5. Screen display and operation: Power on-Function selection-various functional
6. Touch the display screen to choose the functions as your needed

Function Selection:

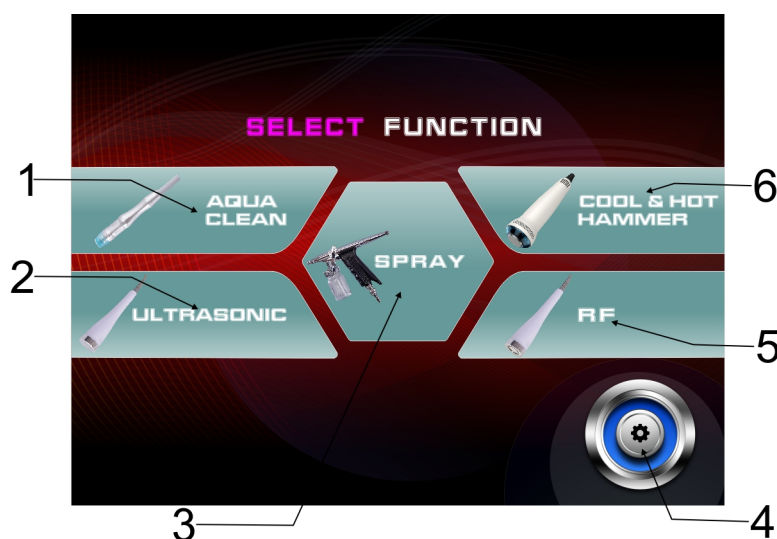

- 1.Aqua Clean
2. Ultrasonic
- 3.Spray
- 4.Setting
- 5.RF
- 6.Cool & Hot Hammer



Skin management process

1. Aqua Clean (Deep cleansing)
2. Spray water oxygen (Inject essence into the skin surface)
3. RF (Lifting the skin.)
4. Ultrasonic (Introduction the essence)
5. Cool & Hot Hammer (Shrinking pores)

Operational steps:


Aqua clean function

Through the negative pressure forming circuit, the ultra-small bubbles and nutrient fluids are fully combined. The specially designed small spiral suction heads directly act on the skin, maintain the ultra-small bubbles for long periods of contact with the skin, promote the peeling effect, can clean the surface deep, and remove the aging keratin cells., removing sebum, Thoroughly remove all impurities, mites and grease residues from the head of the follicle, and supplement the skin with nutrients to make the skin moist, delicate and shiny.

Prep for use


Add the primary liquid or mineral water to the water bottle

1. Press the button in the middle of the bottle.
2. Unscrew the top cover in the direction of the arrow.
3. Add the primary liquid or mineral water to the water bottle
4. Tighten the top cover in the opposite direction and press the middle button (see right image)
5. Insert the bottles into the corresponding sockets of the host. (see below image)

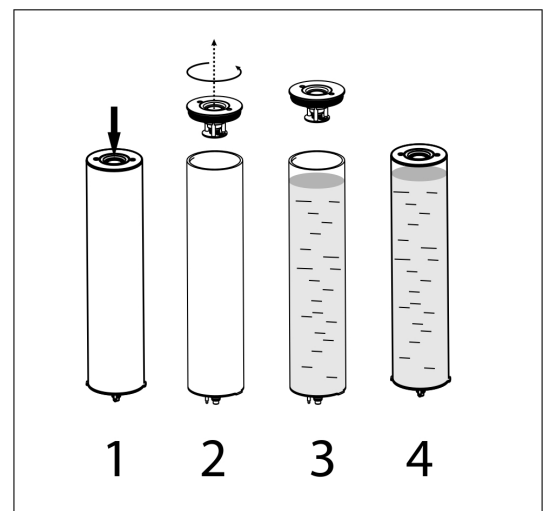


Attention: Check that the bottom sealing ring is intact before inserting(see below image). If the sealing ring is missing or damaged, it will cause leakage! When inserting, pls check a large and small bayonet on two sides, pls match them, otherwise you will not be able to insert it. Insert 4 water bottles in the sockets as you need



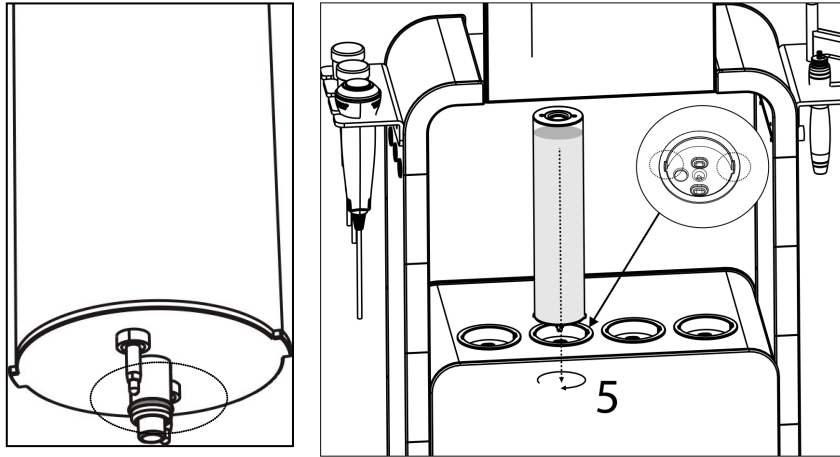

How to use

(1) Touch the Aqua clean on screen and choice the Aqua clean function.

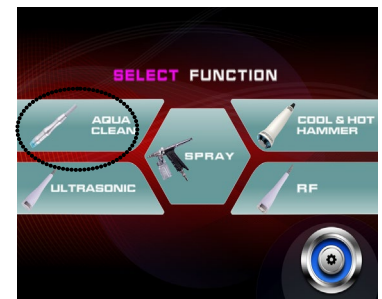

Functional introduction

1. Reduction of time

2. Use time

3. Increase of time

4. Back to Main Menu

5. Energy Reduction

6. Energy bar

7. Percentage of water

bottles used

8. Bottle A


(can turn on hydrogen rich
water function)

9. Bottle B,C,D

10. Bottle working condition

11. Energy Increase

12. Start Key

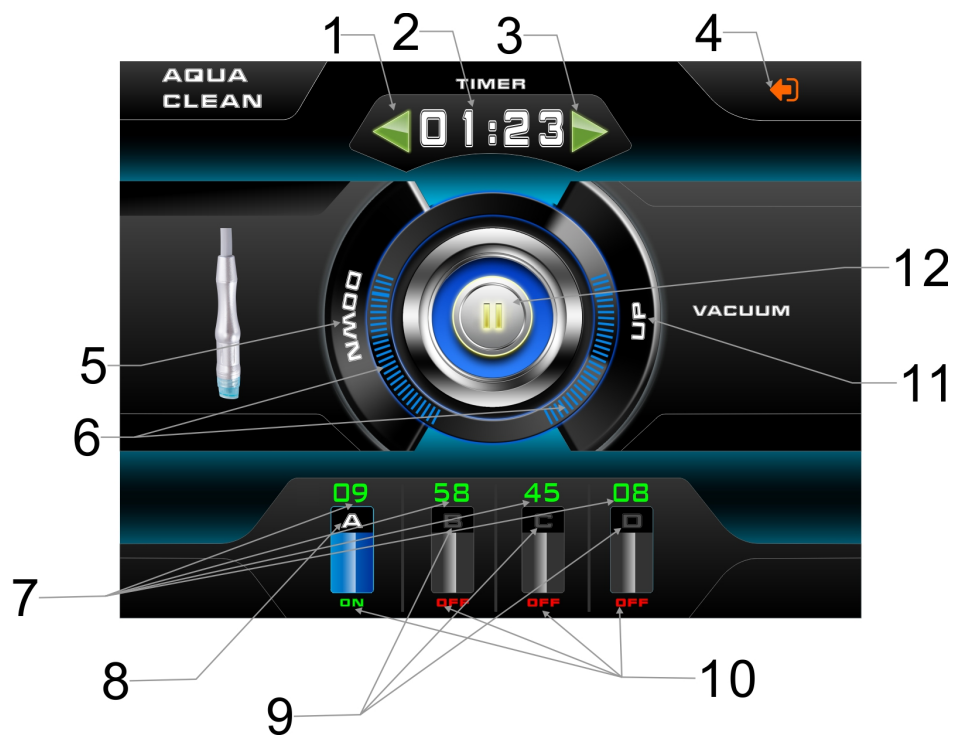



(2). press the button in the middle of the cap (otherwise water won't come out), then press

Reduction of time  or Increase of time  adjust usage time as what you need. (suggest 15-20min), press UP and DOWN adjust suck power according needed (suggest starting from level 2), choice the bottle you need, you can choice one or all, bottle ABCD as show on screen, when open,

the bottom of the bottle is display ON , close display OFF , Bottle A can also turn on the

hydrogen water function, In this state  click the icon again to appear  means that the

hydrogen water has been turned on., when you want to close, just click on the icon to become  can be close. (Noted: Only Bottle A with this function, B,C,D do not have)

Explanation:(Hydrogen Water) solved the scientific problem that hydrogen is insoluble in water with nanometer gas liquid mixing technology, and used physical methods to evenly wrap hydrogen molecules in water to promote a stable combination of hydrogen and water. It has the characteristics of high hydrogen concentration and good stability. Hydrogen is a colorless, odorless, non-toxic, and odorless gas. The unique nature of hydrogen determines that hydrogen has many advantages in biology. A more obvious feature is strong penetration, which can easily enter any part of the cell such as the nucleus and mitochondria, so hydrogen rich water has a very powerful skin function.

Adjust the percentage of water per bottle, Just click the number on the screen , the number will appear on the top



1. Show input number window
2. Keyboard area for input number
3. Cancel an input number
4. Confirm the input No. and Exit



In the number input interface, click the percentage required to enter the number above, such as entering 23, which means that 23 % of the water in the water bottle is used. press  confirmation, system saves the value and exits the interface. If the input is wrong, press  enter again.


When you have adjusted the number what you need, press  start working.

(3) operation

Choice the grinding head you need and fix on the handle



, three grinding head for you choice

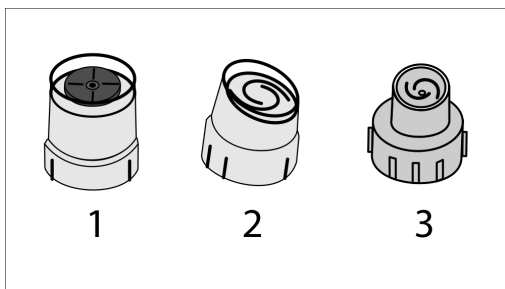

1. Grinding head
2. Hard plastic head
3. Silicone head


Grinding head is suitable for exfoliating dead skin, women's skin is thin, so only do the T area, only do it once, men can operate on the face according to the situation. Hard plastic heads are suitable for regular skin use. Different models are used for cleaning the whole face for different parts. If you need to do nutrient fluid introduction, you need to use the point-pull method.

Silicone head is suitable for sensitive skin, no irritation, and a good sense of comfort.

The whole face cleaning includes sucking black head time in 15-20 minutes. To go black head must cooperate with black head export liquid to soften black head, so that it can be twice as effective




Attentions for Operation:

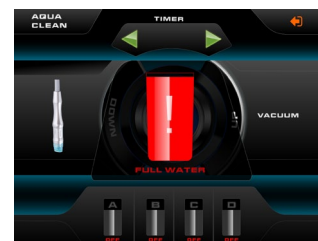

1、 When the waste water bottle is full, there will be a warning logo on the screen and stop working, Please open the door on the back of the machine, remove the waste water bottle according to the arrow direction, clean up the waste water and install it back, and then continue to use it.

2、 The water bottle's seal ring is missing or damaged, which can cause leakage. If it leaks accidentally, please immediately turn off the power , pull out the water bottle, dry the water on the host with dry cloth, and place it in a ventilated dry place until the host is completely dry. Try to boot again. Do not open the machine when the water is leaking! otherwise it will cause the machine to be scrapped! Please pay attention!

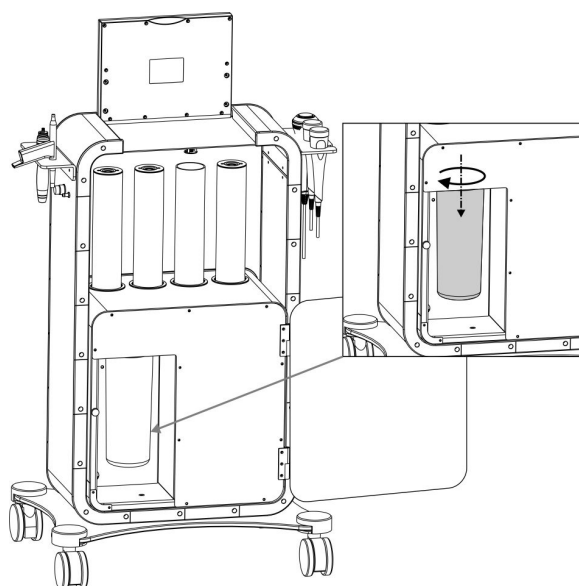

3. When moving the machine, please press the bottle key back to keep the water bottle in a sealed state, to prevent the water from spilling out.

4、 Please keep the grinding head clean, clean the probe

after each use and use alcohol or other disinfection tools to sterilize or replace the new grinding head.

5、Nursing time generally does not exceed 15 minutes. Those with a small area of care should appropriately reduce the nursing time. Too long will make the skin tired and poorly absorbed; The probe should not stay in a certain position of the skin for a long time, so as not to damage the skin tissue.

6、 During the nursing process, it is necessary to avoid beauty for a long time in the state of maximum energy output. This state will not only affect the absorption of the skin, but even damage the skin tissue.


Spray Water Oxygen function

Efficiency: water or nutrient solution into mist, make skin easier to absorb, oxygen-filled water, activated cell skin, nutrition introduction.



Handle introduction

1. Trigger (control the volume of the mist)
2. bottle
3. handle

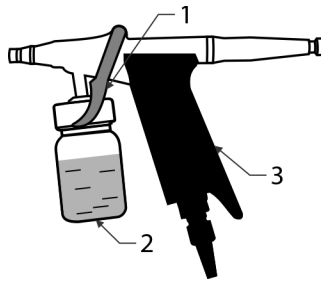


Prep before use



Add pure water or essence to the water bottle


Instruction:

- (1) Choice Spray water oxygen icon at the function selection enters to Spray water oxygen interface

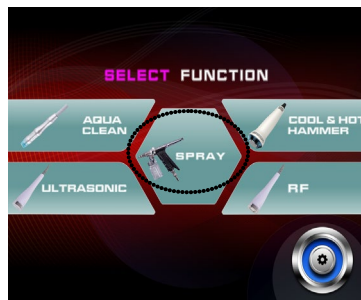

Interface introduction

1. adjusting time
2. start pause keys
3. exit the current function key

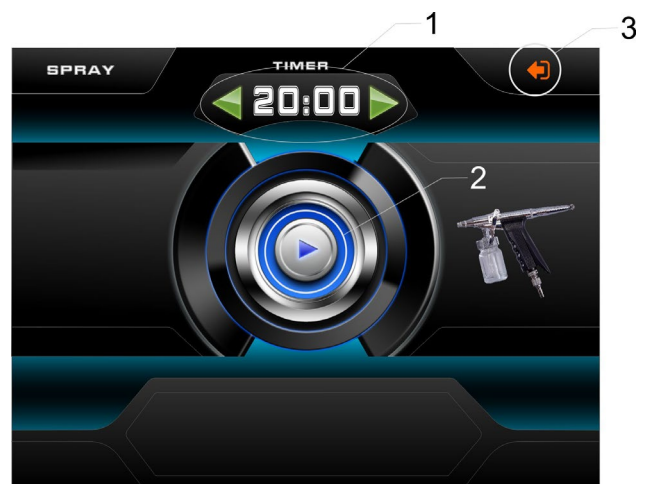

- (2) Increase  or decrease  to the proper time (recommended 15 minutes)



After finish all the setting, pls press  for start.



How to use

- Under start up of host output, Hold the handle and pull the trigger gently, The mist can come out of the nozzle, put different kinds of liquids in the bottle as needed; When change to put different liquid, Pls pause the output of the host firstly, Clean the bottle to put the different liquid inside. After spraying, Pls pause the output of the host, Then use a cotton pad or soft paper to blot up the excess liquid from the customers' face.

Attentions For Use:


It is better to pull the trigger to the middle position to control the mist volume, If pull to the trigger to the end, It is easy to cause drop water. Pay more attention to the important positions such as mouth, nose and eyes when giving facial care to customers, Ask the customers close their eyes and mouth, Spray the mist from the customer's forehead to avoid the liquid to flow into customers' nostrils, Do not forget to keep a proper distance from the customers' face.

RF(Radio frequency)

Functions: Positioning tissue heating, opening pores, dredging lymph nodes, improving fine lines and wrinkles, restoring skin elasticity and enhancing facial contour

Instruction:

- (1) Press RF on the function selection interface to enter

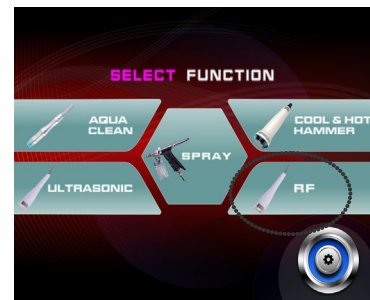

Interface introduction

1. Adjusting Time
2. Energy bar (Shows the energy output)
3. Energy Reduction
4. Energy Increase
5. Output waveform selection
6. Exit current function key
7. Start /pause key



- (2) Increase  or decrease  to the proper time (recommended 30 minutes) , Then press the energy reduction or energy increase to the proper energy, Generally start to select the second gear, To



test it on your hand first before doing treatment for the customers. The customers can feel the slight heat

The output wave state has three modes for choose ,  Continuous wave output state,  Mixed wave output state,  Wide pulse output state. Press the output wave state as your needed, Generally start to select .

Attentions For Use:

1. It must besmear essence liquid, emulsion, the lubrication liquid such as original liquid or gel before use, Do not use it directly on the dry skin, RF head must be close to the skin, Don't let it leave the skin.
2. If the temperature is overheated during use, it only needs to be moved to other parts and move in a circle faster, Then the temperature will drop down; If the temperature is not hot enough, To move RF head in a circle slower, The temperature will rise slowly.
3. It is not recommended for use during pregnancy or lactation. It is not recommended for use in sensitive areas, skin with trauma, redness, inflammation, allergies or other skin diseases.
4. Remove the ear and neck metal accessories before treatment.

Ultrasonic

Efficiency:

- 1、 Soften blood clots red face elimination: application of ultrasonic pulse therapy, cooperate with the appropriate ointment for activating blood circulation and as a result of ultrasonic mechanical action, heat effect and physical and chemical action, for the deformation curve of the blood vessels to appropriate massage, soften blood clots, dilate blood vessels, promote blood flow and metabolism function, can rectify deformation imperceptible blood-vessel, and gradually returned to normal, thereby eliminating blush effect.
- 2、 Soft firm to fat cure acne: there are many forms of acne, is one of the more common with purulent sex and acne, but there is a larger acne, redness is hard, touch very painful, improper handling tends to form hard scar knot. In these patients, it is best not to prick the early extrusion, the application of ultrasonic beauty machine with acne inflammation cream, gently massage the sore body surface, can be a single massage can also be a few together the whole massage, ultrasound can impact the lymph nodes, blood and lymph circulation, accelerate the backlog of inflammatory cells under the action of ultrasonic change shape, combined with ultrasonic antichloristic drug import, swollen form improving the congestion phenomenon of acne, soften the skin and indurations.
- 3、 Eliminate freckles and lighten spots: ultrasound can destroy the inner lining of color spot cells, prevent it from forming and decomposing pigment, so as to weaken and disappear with metabolism. Ultrasonic can enhance the absorption of nursing products (such as Chinese medicine face cream, high unit vitamin c, freckle cream, etc.), freckle removal effect is significant.
- 4、 Decontamination and desorption: the mechanical massage function of ultrasonic is used to regulate the subcutaneous cells and enhance the permeability of cell membrane on the one hand, on the other hand, strengthen the blood circulation and metabolism on the other hand. Combined with proper blood circulation and desalting nursing products, the effect is significant.



- 5、Eliminate pouch and black rim of the eye, due to the effect of mechanical massage of the ultrasonic beauty instrument can accelerate the blood and lymph circulation, and thermal effect, can refine and break down fat, accelerate skin absorption, make excessive accumulation of moisture and fat to dissipate, pouch and black rim of the eye will disappear.(use ultrasonic beauty instrument when doing eye care, probe should be along skin texture to massage, each pouch care not more than five minutes, efforts to gently, probe can't eye direction.)
- 6、Often do massage to improve the cortex, use ultrasonic beauty instrument to undertake proper treatment, can eliminate spot sore as soon as possible, and make the skin is improved or restored.

Instruction:

- (1) Press ultrasonic on the function selection interface to enter

Interface introduction

1. Adjusting time
2. Energy (Shows the energy output)
3. Energy reduction
4. Energy increase
5. Output wave state selection
6. Exit current function key
7. Start/pause key



(2) Increase  or decrease  to the proper

time (recommended 20 minutes), Then press the energy reduction or energy increase to the proper energy, Generally start to select the second gear

(3)The output wave state has three modes for choose,  Continuous wave output state,  Mixed wave output state,  Wide pulse wave output state。 Press the output wave state as your needed,

Generally start to select .

Put the lotion or essence on the skin before using ultrasonic probe, To move ultrasonic probe in a circle from the bottom to up。

Attentions:

- 1、Keep the ultrasonic probe clean, clean it after each use and sterilize it with alcohol or other disinfectant.
- 2、Do not put the ultrasonic probe on the eyeball。
- 3、Make sure to put enough care essence on the skin during treatment.



- 4、 Generally, the nursing time should not exceed 15 minutes. Those with a small area of the nursing area should reduce the nursing time appropriately. The probe should not stay in a certain position of the skin for a long time to avoid damaging the skin tissue.
- 5、 In the process of nursing care, it is necessary to avoid doing beauty treatment for a long time in the state of maximum energy output, which will not only affect the absorption of the skin, but even damage the skin tissue.

Cool and hot hammer

Efficiency:

Cold therapy


- 1、 It is good for constriction of blood vessels, tight capillary, thick capillary。
- 2、 Strengthen the raw fiber, make the skin become tight and elastic, smooth wrinkles, crow's feet。
- 3、 The harmful waste products of constricting blood vessels, effectively reducing the floating of red veins and recirculation of veins at any time are easier to be purified and detoxified。
- 4、 Lower metabolic rate, inhibit the activity of slow cells, reduce nerve sensitivity, have the effect of calm, offset relieve skin, swollen, bilge, pain, especially for allergic skin has good curative effect.
- 5、 Acne removal red blood silk, cold treatment + blue light can ease the skin evacuation red blood silk, anti-inflammatory acne

Hot therapy

- 1、 Slight heat, make a person feel warm, comfortable, can relax muscle certainly, balance self-discipline nerve, solve ache, ache。
 - 2、 Stimulate the blood circulation to accelerate, increase the amount of oxygen in the blood, open the pores, expand the blood vessels, and make the nutrients on the skin absorbed quickly, which is the most gentle, comfortable and effective way to import
 - 3、 Increase the skin metabolic rate, so that the tissue contains heat energy to increase, natural elimination of waste toxins, on the elimination of pouch, black eye circles have certain effects.
- Cold and hot therapy is used alternately, have the effect that adds, multiply to hairdressing, can activate cell completely naturally, improve, purify skin


Instruction:

- (1) Press cool and hot hammer on the function selection interface to enter

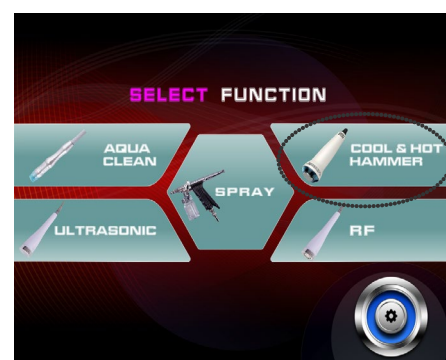

Interface introduction

1. Setting temperature
2. Start/pause
3. Cool hammer
4. Hot hammer
5. Displays the current probe temperature
6. Exit the current function key

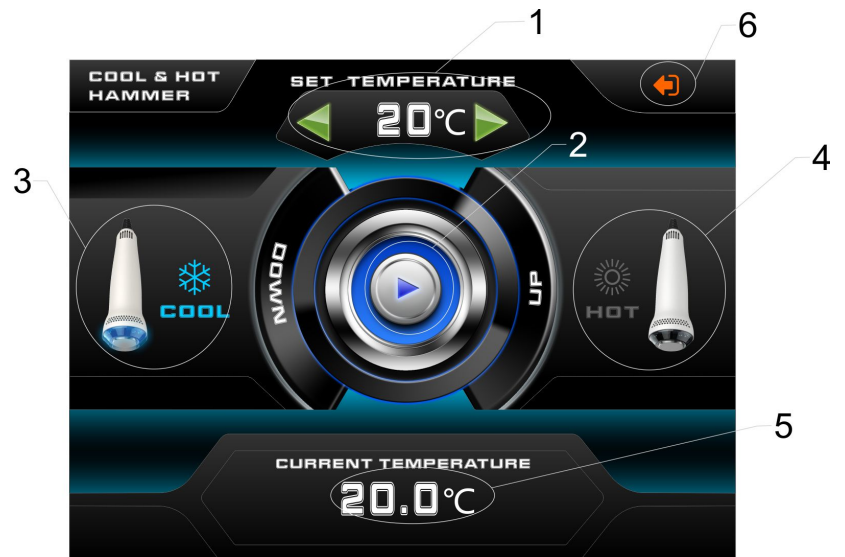

(2) Operating steps

- Increase  or decrease  to the proper temperature, Press  to start cool

therapy, Press  to start hot therapy. The temperature range of the hot hammer probe: 40~55 degree, The temperature range of the cool hammer probe: 15~20 degree, The current probe temperature

is displayed in the interface . If working in cool therapy, The probe of the hammer

will emit blue light . If working in hot therapy, The probe of the hammer will emit red light . We can distinguish the working state at this time from the light color of the probe, If you stop working, the lights go out. It also displays the working state on the screen.

Attentions:

1. Clean your face with makeup and dust to get the best effect.
2. Don't let the hot and cold hammer get wet. Please wipe it off with dry cloth.
3. Please keep the good heat emission, Do not block the heat emission holes of the handle



4. This device is not intended for use by persons with weak physical ability, slow reaction or mental disorders (including children) unless it is used safely under the direction or assistance of persons responsible for its safety.
5. Children should be supervised to ensure that they do not play with equipment.

Attentions should be paid in using for the whole machine:

1. When working, pay attention to the attachment and the host connection, not excessive pulling.
2. If the power cord is damaged, it must be replaced with special wire or special components purchased from the manufacturer or its maintenance department.
3. In the working process of the host have abnormal situation is found, please power off and unplug the power plug immediately, please contact your dealer near check repair or professional maintenance staff to check, please open the casing maintenance without authorization is prohibited.
4. Do not place or use in the bathroom or in wet places.
5. Do not place with inflammable or explosive objects.

In case of any discrepancy between the drawings on this manual and the real objects, the real objects shall prevail. The company reserves the right to continuously update and improve the products. The above contents are subject to change without prior notice.

