

Limited Warranty

True Innovations seating products are warranted against defects in material and workmanship while owned by the "Original Purchaser". The term "Original Purchaser" is defined as that party or entity which purchases True Innovations furniture from an authorized True Innovations retailer as shown by the **original sales receipt**. This warranty does not apply to "as is" or previously owned seating products.

True Innovations warrants its products to be free from defects in materials and workmanship (normal wear and tear, fading and stretching excepted), for one year from the date of purchase, except as provided below under the limited warranty.

The exceptions or exclusions to this warranty are as follows:

1. Any warranty claim that is submitted after the limited warranty period, or without proper proof of purchase will not include the cost of parts, labor, or delivery.

2. This warranty shall not apply to the True Innovations products which have been subjected to misuse, neglect, alteration, modification, or attachments, either caused by the original purchaser, shipping, storage, accident, fire, flood, or acts of God.

In the event of a defect in material or workmanship covered by this warranty, True Innovations will repair or replace, at their discretion, without charge the warranted furniture within a reasonable amount of time. We reserve the right to require damaged parts to be returned to us upon request. Natural markings such as scars, brands, grain variations, wrinkles, color variations, etc. will be considered normal characteristics and not construed as defects. We cannot warrant leather against scratching or scuffing, as all leather is subject to blemishes with use.

For claims or questions concerning this warranty, please contact True Innovations Customer Service department on the web at www.trueinnovations.com. Be prepared to provide your name, address, daytime phone number, the model and purchase order number of the product, date of purchase, and the store where it was purchased. True Innovations requires that the original sales receipt be submitted with all requests to confirm that you are the original purchaser according to our written warranty.


This warranty gives you specific rights, and you may also have other rights, which vary from state to state.

Should you prefer to speak with a service representative, you may call toll free: 1-800-379-9773, or fax 714-544-5118.

**Better Homes
& Gardens®**

Lot Number :
Date of purchase :

Power Lift Recliner Chocolate Brown Assembly Instructions

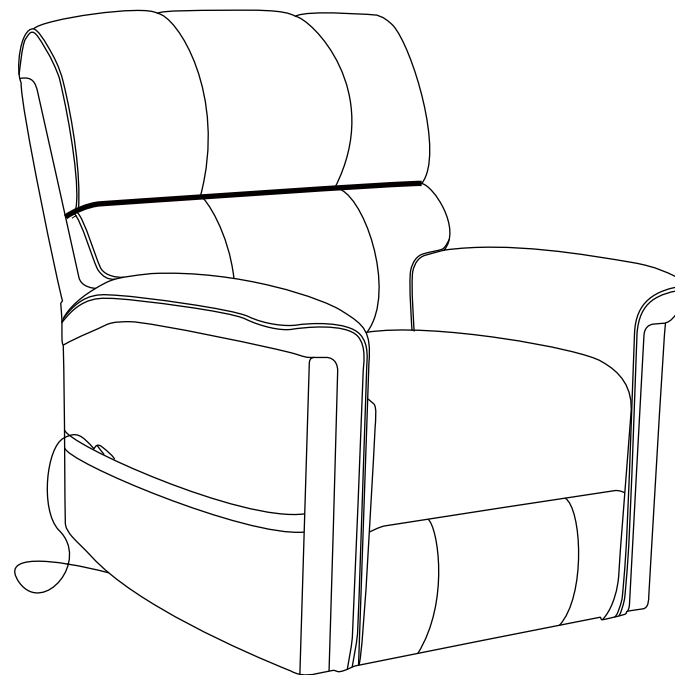

Model: BHT-10002

THIS INSTRUCTION BOOKLET CONTAINS **IMPORTANT** SAFETY INFORMATION.
PLEASE READ AND KEEP FOR FUTURE REFERENCE.

Customer Service Support: 1-800-379-9773





Thank you for selecting one of True Innovations seating solutions.

It's True...Fast, Friendly Customer Service

At True Innovations, we are committed to excellence. Our primary objective is to provide fast, friendly customer service to all of our customers.

Call us first!!


DO NOT RETURN TO STORE

For immediate help with assembly or product information call our toll- free number:


1-800-379-9773

Mon-Fri 6:00 am to 5:00 pm Pacific Time

Included Parts

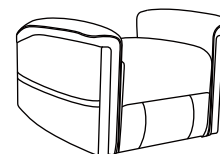

A Seat Cushion (1)

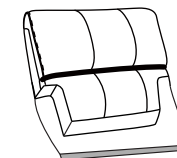


B Back Cushion (1)



C Power Cord (1)



D Remote Control (1)

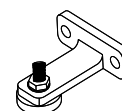

E Reinforcement foot (2)



F1 1/4"x1 3/4" L screws (4)



F2 butterfly nut (4)


NOTE:

Power Cord (C) and Remote Control (D) are located in the back of the seat cushion.

WARNING: Do not throw box away.




ELECTRIC POWER SUPPLY
ALIMENTATION ELECTRIQUE
NOTE: PLS DO NOT THROW AWAY

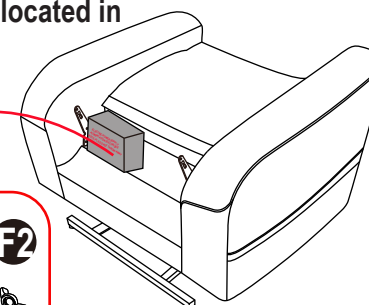




Step 1 Attach Base metal frame of seat cushion (A) to the Reinforcement foot (E) using $\frac{1}{4}$ "x $1\frac{3}{4}$ " L screws (F1), Fully tighten all screws with butterfly nut (F2).

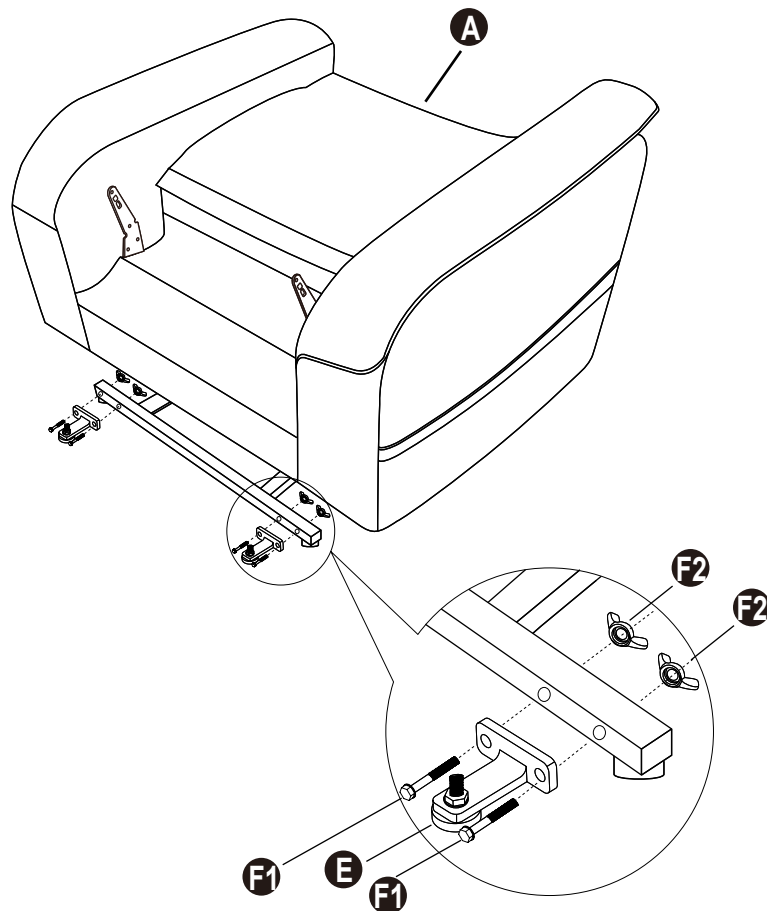


Step 2 Attach the Back Cushion (B) to the Seat Cushion (A) by sliding the receivers into the mechanism posts located on the Seat Cushion (A). From the front of the recliner, press down on both sides of the Back Cushion (B) to firmly secure. The Back Cushion (B) and Seat Cushion (A) will automatically lock into place. Then pull the back flap down tightly and attach the back flap to the bottom of the Seat Cushion (A).

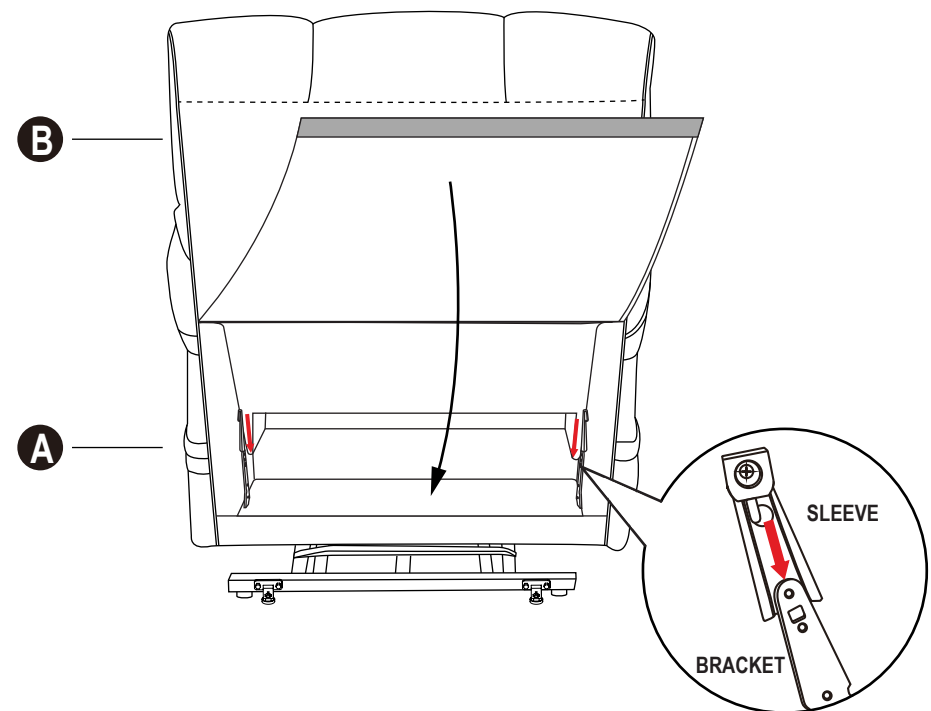




Step 3 Connect Female Plug and Male Plug of for power cord(C),connect Female Plug and Male Plug of for remote control (D).


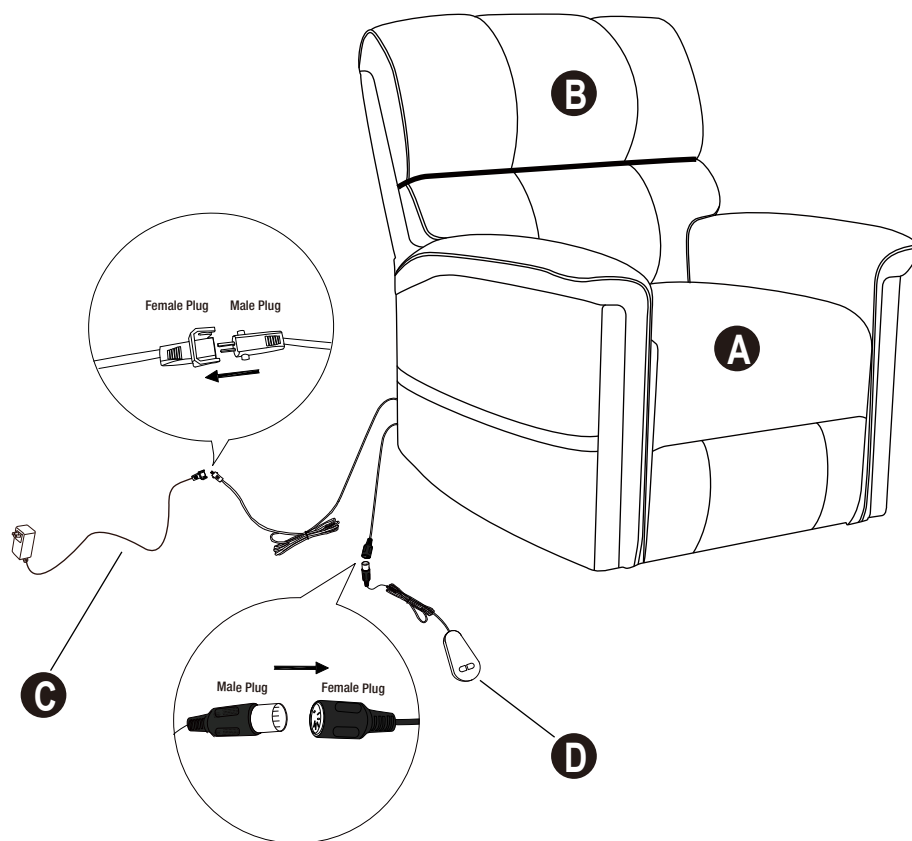


Step 4 While the cables connected as below, the lift recliner can work. When you're not using the remote control (D), you can put it in right side pocket.

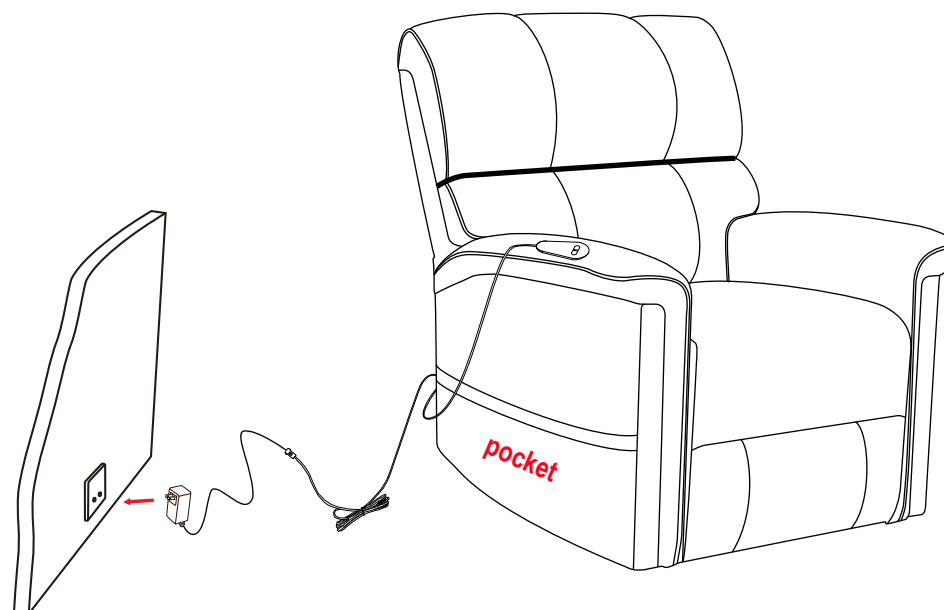





FUNCTION 1: RECLINE

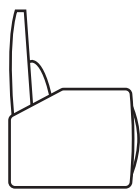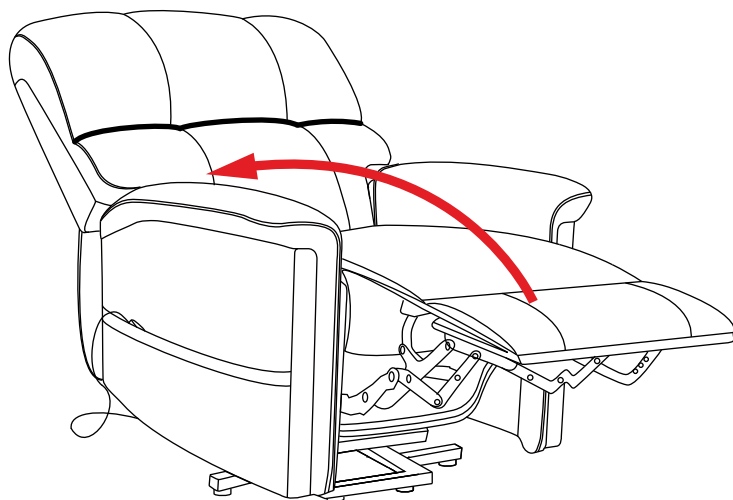

SIT

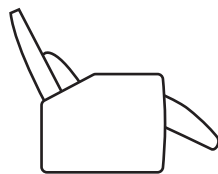

TV POSITION

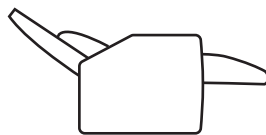


FULL RECLINE


⚠ WARNING

Do not get off from the recliner while in a reclined position to prevent fall.
Always put the recliner back to upright position and foot rest closed before getting off.

FUNCTION 2: LIFT

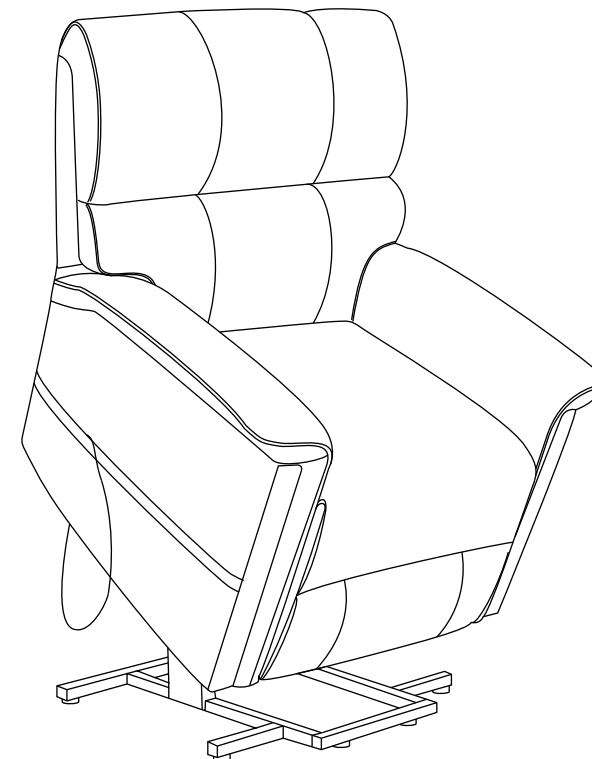**⚠ WARNING**

As the chair slopes forward slowly, get off it.

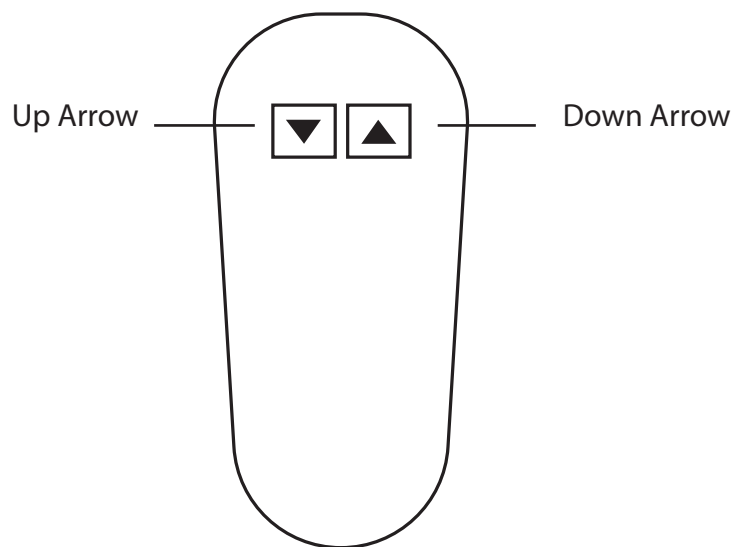

Up Arrow: Hold down up arrow to manually fold chair up and extend to maximum lifted position.

Down Arrow: Press and hold down arrow to unfold chair and recline chair to fully reclined position.

NOTES:


- Please read the instruction sheets completely before assembly.
- Remove all parts from carton and ensure all parts are included prior to assembly.
- Assemble the chair on the carton or carpet to prevent scratching the chair.
- Be careful not to cross-thread screws.
- Do not use power tools.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this futon only for its intended use as described in the operating instructions.
- Do not use outdoors. INDOOR USE ONLY.



⚠ WARNING

- Do not leave children unattended on or beside the recliner.
- Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
- Keep dry. Do not operate in a wet or moist environment.
- Do not stand on the seat cushion.
- Please do not place this product near sources of excessive heat.
- Please avoid use of sharp objects on or near this product.
- For persistent stains and overall cleaning, contact a professional upholstery cleaning company.
- Do not sit on the edge of the seat, armrests or on the backrest as this could make the recliner tip over.

- All warranties are limited to the original purchaser for normal use which is defined as a standard work week of forty hours by a person not exceeding weight capacity of 225 lbs (102kgs).

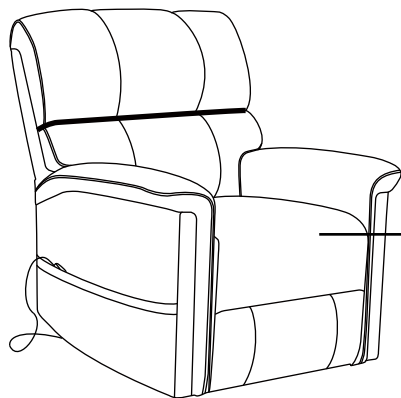

Maximum weight load
capacity is 225 lbs (102 kgs)