

Intelligent visual ceramic ear spoon V105 User manual

open

①

Scan the QR code and
download the app



②

Open the WIFI Visible
and click the 'start'



③

Hold the power button
to turn it on




Hold the power
button for 3sec



WiFi connect


①

Connect the
WiFi “WONSHAN-XXXX”

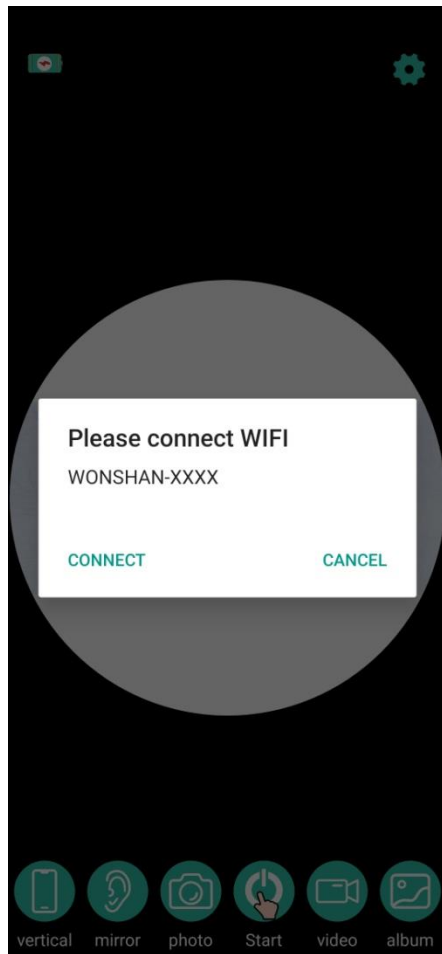

②

Return to the app and click
the _x0002_start _x0002_
button to access the app

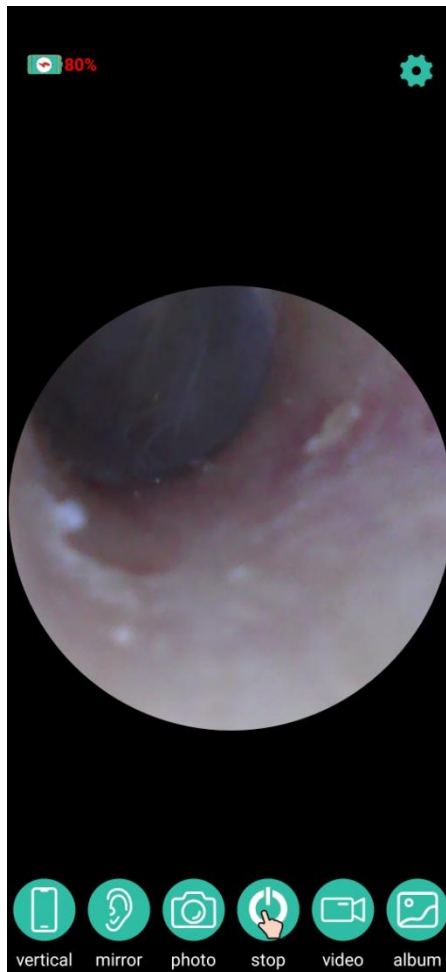

Replace the ear spoon head

Spiral type ear
spoon head

When
installing
the ear
spoon,
tighten it
clockwise

Be sure to
install it
before use to
prevent the
ear spoon
from falling
off



Remove the
ear spoon
by turning
it
anticlockwi
se


Accessories display

①

PC transparent ear spoon x 2



Product size

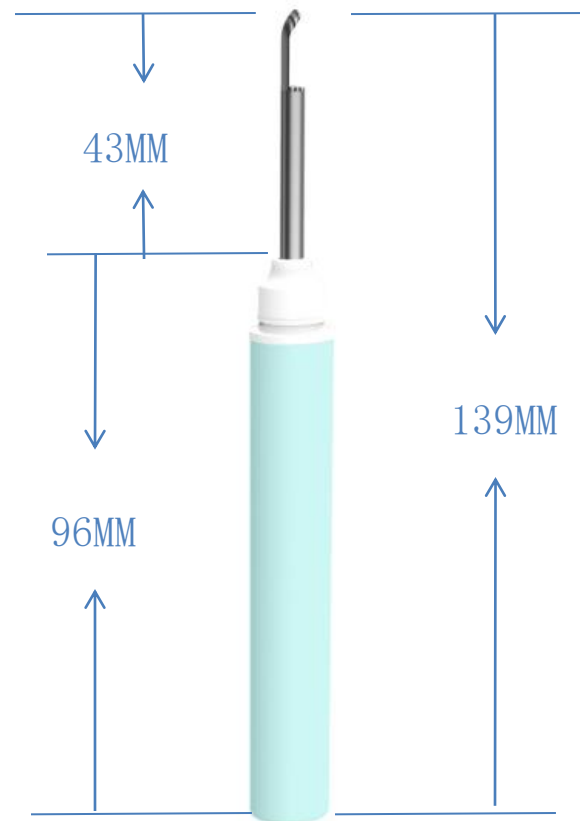

Trouble shooting

Failure phenomenon	Solution
Power on. No response	• Try Turning it on a few more times; • Please use your own charging device for charging, Charging time is 1 h.
The camera lights aren't working	• Please reboot the device; Check that the equipment is turned on.
app quits unexpectedly	• Reboot your phone; Or uninstall already installed APP, Scan the Qr code again to download the installation.
The equipment indicator light is red	• Please use your own charging device for charging, Charging time is 1 h.
After the device is connected to the appReal time image not displayed	• Disconnect the phone from the device WiFi, Reset the device, Reconnecting WiFi.
WiFi connection is unstable, Off and on	• Please restart the app and the device, Reconnecting device WiFi .
If the above operation fails to solve the problem, Please contact our customers. We will provide further solutions for you.	

Points to note

- 1、Cleaning the Lens, Professional Alcohol Swabs are recommended.
- 2、Please use this product in an environment where no one is running or fighting, To avoid impact.
- 3、This product is not suitable for children under 8 years of age.
- 4、Children are not allowed to use this product alone, So as not to inadvertently cause harm.
- 5、Do not use endoscopy in liquids, To prevent damage from water ingress.
- 6、This product has built-in rechargeable lithium battery, If not used frequently, Must also be recharged once a month.
- 7、In use, The temperature of the product will rise slightly (Up to 35 °C) , Please feel free to use.
- 8、To ensure the best use experience, When picking ears for children, Please keep the child in a sitting position.
- 9、For your ear canal clean and safe, The tip of the spoon is loose or used more than 30 times, Be sure to replace the ear spoon.



Product parameters

Brand

Product name	Intelligent visible ear pick stick
Model	E8
Host weight	19g
Network Standard	IEEE 802.11b/g
Antenna	Built-in antenna
Operating frequency	2.4Ghz
Image transmission rate	30fps
Image Sensor	CMOS
Working Environment temperature	-10-50 ℃
Battery	350mAh lithium battery
Battery life	About 120 minutes
Charging time	About 1 hour
Input current	DC 5V 500 mA
Lens diameter	4.5 mm
Optimum focal length	1.5 ~ 2 cm
pixel	3 million
Gravity sensor	3 axis

After sale warranty

User information	
User name:	contact number:
To address:	
Sales merchants:	contact number:



This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

- Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - -Reorient or relocate the receiving antenna.
 - -Increase the separation between the equipment and receiver.
 - -Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - -Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

