QUICK START GUIDE

BIG TalkingBrix™ speech device

WARRANTY

AbleNet manufactured products include a 2-year limited warranty. This warranty is against defects in materials and manufacturing for 2-years from the date of purchase. Full warranty details available at www.ablenetinc.com.

AbleNet, Inc.
2625 Patton Road
Roseville, MN 55113
United States of America
(651) 294-3101
ablecare@ablenetinc.com
www.ablenetinc.com

Made with recycled content

EN June 18, 2021 1:54 PM

Product Registration

Registering your product provides you access to ableCARE, product updates, and resources for your product. Scan the QR code below to register your product.

Getting Started

Scan the QR code below to watch a short getting started video or follow the instructions listed.

To get started:

1. On back of device, move the switch to REC.
2. Press and hold colored switch top.
3. Begin speaking for up to 10-seconds when light starts blinking.
4. Release colored switch top when finished.
5. On back of device, move switch from REC to ON to start using.
6. Your device is ready to use.

This device can do more! Full Instructions for Use are available at www.ablenetinc.com.


ableCARE™

Download the ableCARE app to access live support from the ableCARE Product Success team, an online knowledge base filled with videos and getting started information, and other resources.

Scan the QR code below to download the free ableCARE app to your phone or tablet.

3

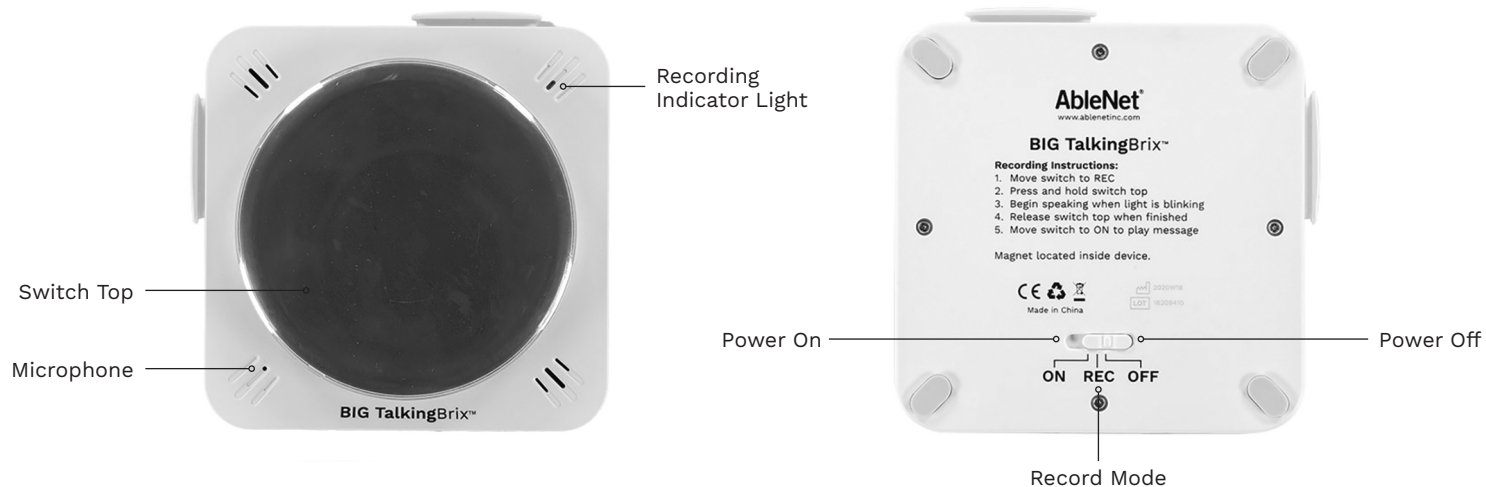Device
Overview

Download on the
App Store

Create printable picture symbol overlays for your device with the *FREE* AbleNet Symbol Overlay Maker app. Scan the QR code to download today!

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