

## XT Auto




**Notice!** Please remove the guide before giving the device to the patient. Only medical professionals could adjust pressure settings.

### Provider Setup Navigation

1. Press **"Menu"**  to any  Menu (See Operation Flow on the right) while in the Standby screen.
2. Hold down both the **"UP"**  & **"Down"**  buttons, and then simultaneously press the **"Menu"**  for 1 second. When the screen starts blinking, the Provider Mode would be unlocked.
3. Press **"UP"**  and **"Down"**  to adjust the settings.
4. Press **"Menu"**  to confirm the settings.
5. Press the **"Start/Standby"**  to go back to standby screen.


## Provider Mode Quick Guide

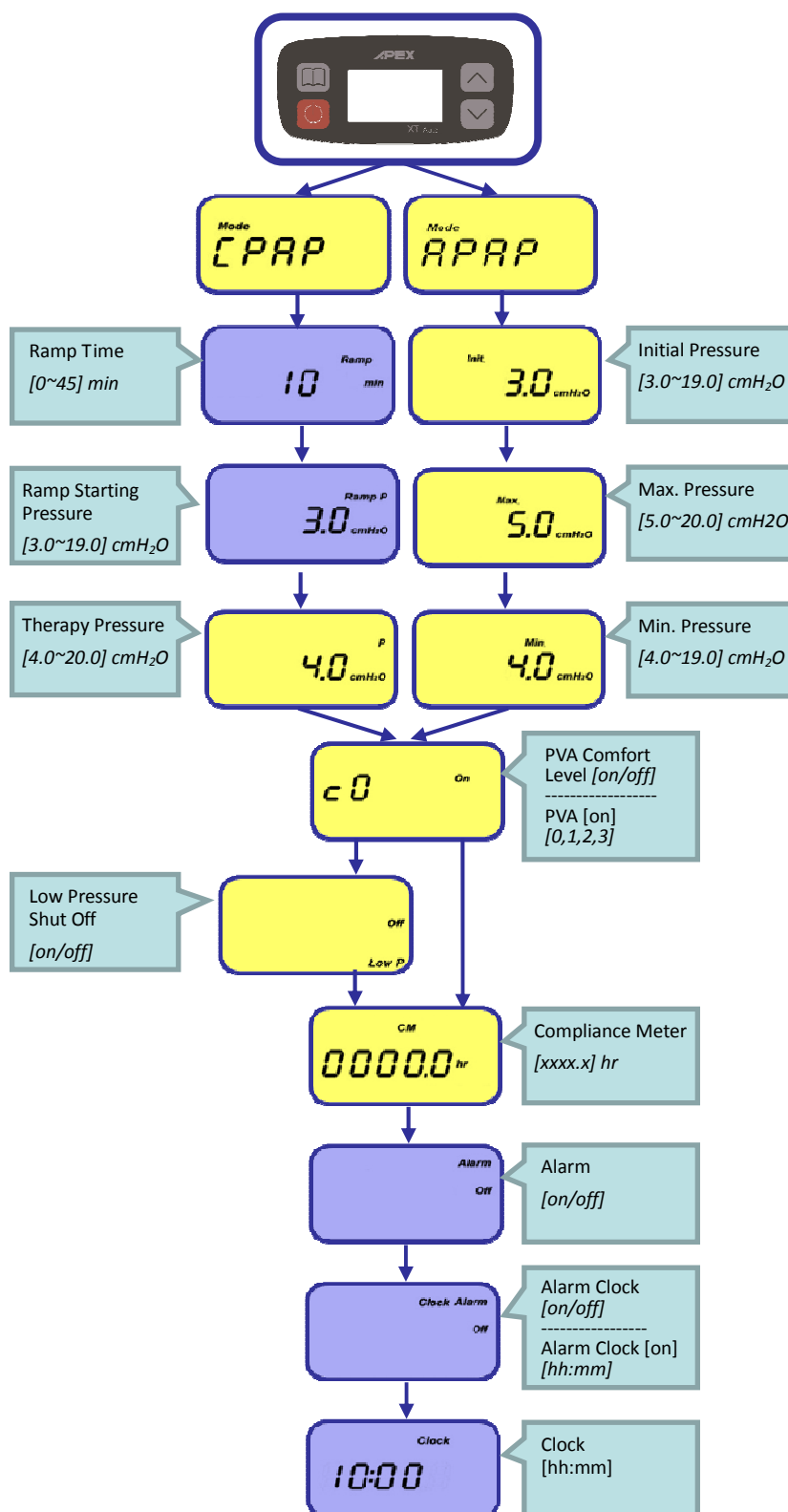

## Provider Mode Quick Guide

### Re-set the Compliance Meter

1. In the standby screen, press "Menu"  to select [CM xxxx.x hr] menu.
2. Accessing the Provider Mode first (Press   + ), then the screen starts blinking.
3. Press "UP"  or "Down" , then the screen will show [Clear].
4. Press "Menu"  to confirm the settings, and the Compliance Meter record will be erased. The screen will show [Clear OK] when finishing the process.

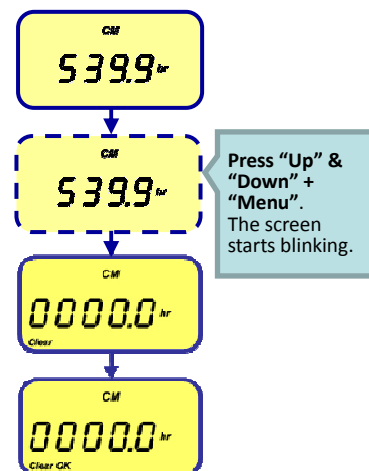

## Engineering Mode Quick Guide

### Accessing the Engineering Mode

1. Plug the power supply into the iCH device.
2. Before it shows [Standby] screen, hold down both the "Starts/Standby"  and "Menu"  buttons for 1 second, the Engineering Mode would be unlocked.

### Re-set the Total Meter

1. Accessing the Engineering Mode first, then the screen will show [TM xxxx.x hr] menu.
2. Hold down both the "UP"  & "Down"  buttons, and simultaneously press the "Menu" , and then the screen starts blinking.
3. Press "UP"  or "Down" , then the screen will show [Clear].
4. Press "Menu"  to confirm the settings, and the screen will show [TM Clr Clear], and the Total Meter record will be erased.
5. The screen will show [TM End Clear OK] when finishing the process.

