


Etac Crissy


Manual English

75369G 21-05-10

etac®

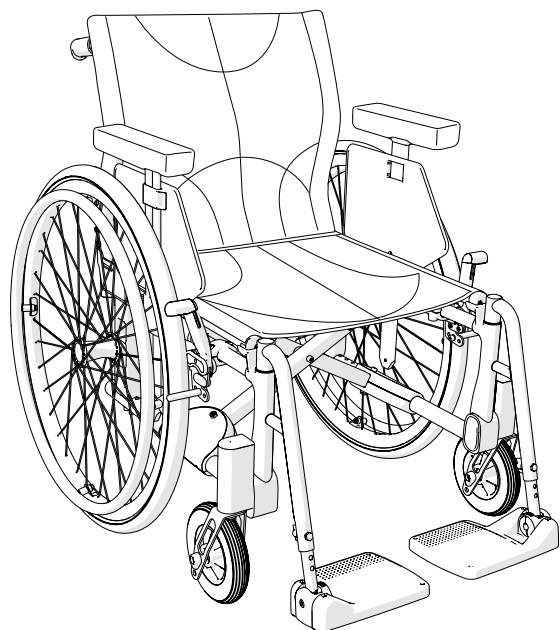

Etac Crissy Swing-Away

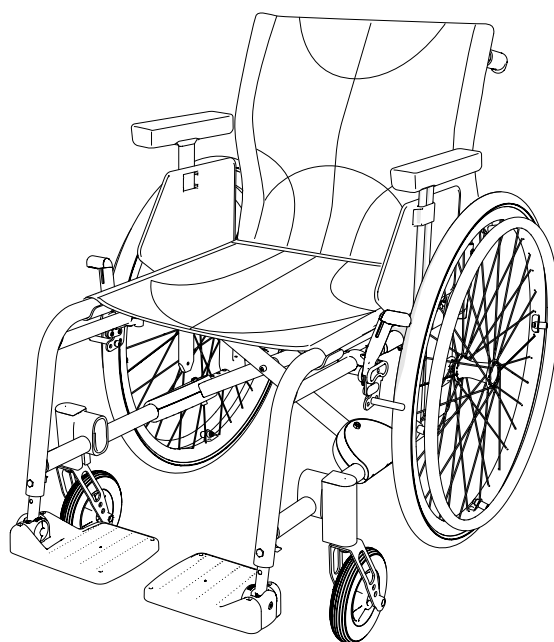

Etac Crissy Active

Contents

1. General – introduction		7. Fitting seat extender	
Intended use.....	3		34
2. Description		8. Settings Back support 2B	
.....	4		35
3. Chair options		9. Settings Back support 2C	
.....	5-7		36
4. Accessories		10. Readjustments	
.....	8-10	Brake, adjustment.....	37
5. Installation upon delivery		Adjustment of brake pin in certain rear wheel positions ...	38
.....	11	Variable angle-adjustable leg support (option).....	38
6. Settings		Amputee leg support (option).....	39
Front seat height.....	12	Anti-tip.....	39
Seat height adjustment using front fork attachment	12	Changing the camber angle.....	40
Seat height adjustment using the castor wheel's position in the fork	12	Toe-in	40
Seat height adjustment by replacing the front fork/castor wheel	13	Handrims	41
Front seat height table	13	Adjusting the distance of the handrims.....	41
Disassembly/assembly of front fork	14	Dismantling/mounting seat- and back support upholstery ..	42
Angle setting of front fork attachment.....	14	11. Accessories - adjustment, handling	
Rear seat height.....	15	Seat cushion and comfort wedge	43
Seat height adjustment using rear wheel position/rear wheel	15-16	Positioning belt	44
Rear seat height table	17	Head support.....	44
User manual		12. Technical data	
.....	19-30	Etac Crissy Swing-Away, Etac Crissy Active.....	45
Symbols and warnings	20	13. Service and maintenance	
Transportation in vehicles.....	21	Cleaning	46
How to use	22-23	Disinfection.....	46
Seating comfort settings	24-25	Rear wheels/castor wheels	46
Driving technique, use	26-29	Brakes.....	46
Care and maintenance.....	30	Touch-up paint	46
Troubleshooting chart	30	Miscellaneous	46
6. Settings (continued)		14. Tests and guarantees	47
Seat angle.....	31		
Seat depth	31		
Seating comfort	32		
Adjustment of the centre of balance	32-33		



1. General

Thank you for choosing a wheelchair from Etac. To avoid injury during handling and use, you should read this manual and keep it for future reference. You can also find it at www.etac.com. Here there are also updates concerning the product.

In this manual the User is the person sitting in the chair. The Carer is the person helping the user.

Explanation of symbols and general warnings can be found on page 20 in this manual.

Etac continuously makes improvements to its products. We therefore reserve the right to modify the products without prior notice. The dimensions given on drawings or other material are solely for guidance.

We cannot be held responsible for any printing errors or omissions.



This symbol can appear alongside text in the manual. It highlights actions where the wheelchair's, the user's or the carer's safety can be put at risk.



This symbol can appear alongside text in the manual. Here you can find advice and tips worth considering.



To avoid injury during handling and use, this manual must always be supplied with the wheelchair.

Explanation of the list of items


- Available options
- ⊕ Included as standard

Etac Crissy Active

Etac Crissy Swing-Away

Swing-Away
Active

In the middle of this manual, you will find a section containing user information.

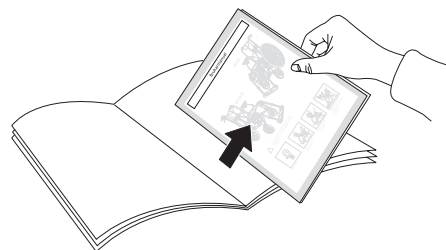

Tools needed (available as accessories)



Allen key
3, 4, 5, 6 mm



Ring spanner
8, 24 mm



Socket
spanner
19 mm

Intended Use:

Etac Crissy Swing-Away and Etac Crissy Active are manual, cross-folded, all-round wheelchairs with multiple functions, and are intended for use both outdoors and indoors. The wheelchairs can be used by people with restricted mobility or by people with walking difficulties. The chair seats are adjustable in height, angle and depth. The back supports are adjustable in height and shape. The foot supports are adjustable in height, angle and depth.

Contraindications:

Users with dementia or cognitive disabilities should not move in and out of the wheelchair without assistance.

In case of an adverse event occurred in relation to use of the device, it should be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to the manufacturer.

Etac Crissy Swing-Away and Etac Crissy Active can be supplemented and adjusted if needs change. A large range of options and accessories are available: various types of back support and leg support, anti-tips and various types of handrims and brakes, etc.

Etac Crissy Swing-Away and Etac Crissy Active can be reconditioned and re-used. Reconditioning instruction is available at www.etac.com.


Etac Crissy Swing-Away and Etac Crissy Active have the best prerequisites to create comfort, functionality and good manoeuvrability.



This product can be adjusted beyond safe use. The prescription of Crissy must be done by medical personnel and tested by a specialist.

2. Description/Standard version

Etac Crissy Swing-Away

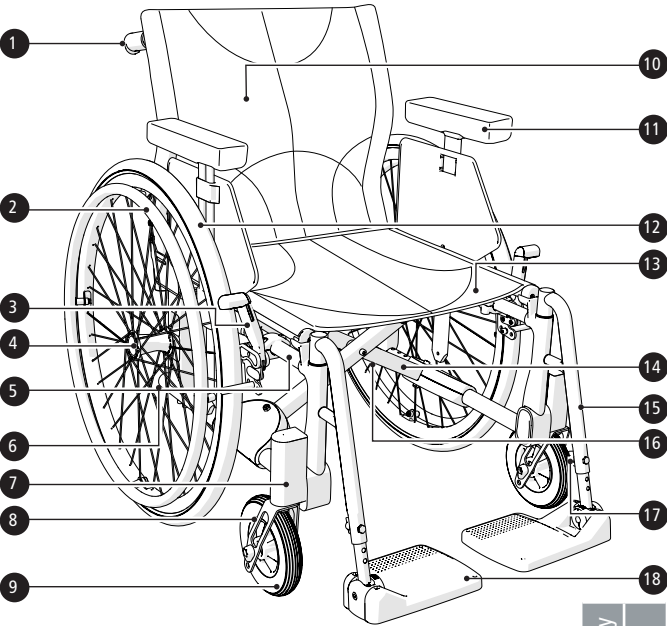

		Swing-Away	Active
1	Push handle	●	●
2	Handrim	●	●
3	User brake	●	●
4	Quick release hub	●	●
5	Side frame	●	●
6	Tilter	●	●
7	Front fork attachment	●	●
8	Front fork	●	●
9	Castor wheel	●	●
10	Back support upholstery	●	●
11	Arm support	●	●

Etac Crissy Active







		Swing-Away	Active
12	Rear wheel	●	●
13	Seat upholstery	●	●
14	Cross	●	●
15	Detachable leg support (on Crissy Active, fixed and integrated i ramen)	●	
16	Serial number	●	●
17	Knob footplate lock	●	
	Locking nut footplate lock		●
18	Footplate	●	●



3. Chair options









In the "Chair options" section you can see all the adjustments that can be made to the wheelchair. In addition, there are so-called "Non-standard combinations", e.g. for rear wheels, castor wheels and back support models. Contact customer service.

For more information about chair options, see www.etac.com



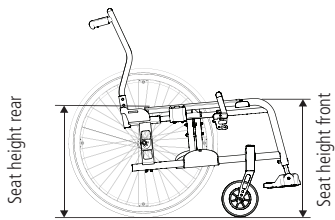
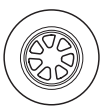
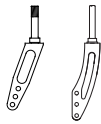
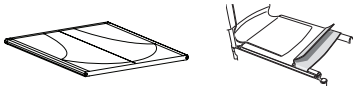
				Swing-Away	Active
Seat height		Front	40.5-55 cm	●	●
		Rear	38.5-51 cm	●	●
Seat angle		Max. forward angle: 2 cm lower front than back (3°) Max. backward angle: 6 cm higher front than back (9°)		●	●
Castor wheel		5" solid		●	●
		6" semi-solid		●	●
		8" semi-solid		●	●
Front fork		Short (three height settings for castor wheels)		●	●
		Medium (three height settings for castor wheels)		●	●
Rear wheel		22"		●	●
		24"		●	●
		24" Carbon (not with attendant brake)		●	●
Tyres		High pressure pneumatic 22" or 24" x 1 1/4		●	●
		High pressure pneumatic 22" or 24" x 1		●	●
		Solid		●	●
		Tyre with solid inner tube (only with attendant brake)		●	●
Handrims		Aluminium		●	●
		Stainless steel		●	●
		Plastic		●	●
		Foam rubber (not with attendant brake)		●	●
		Titanium (not with attendant brake)		●	●
Camber angle		2°		●	●
		0°		●	●
Seat		Padded seat, nylon, 10 mm foam		●	●
		Adjustable seat upholstery		●	●
Fixed back support with fixed push handles, 2B		Height of low back support 35-40 cm Height of high back support 42.5-47.5 cm Fixed angle +5°		●	●
Angle adjustable back support		Height of low back support 35-40 cm Height of high back support 42.5-47.5 cm Angle -5° – +5°		●	●
Back support 2C		Height 37-47 cm Angle +5°		●	●



3. Chair options

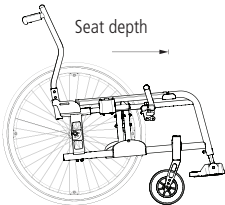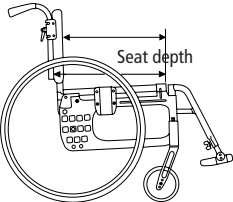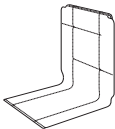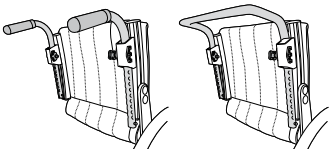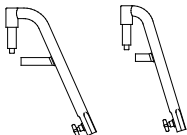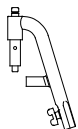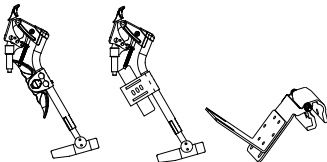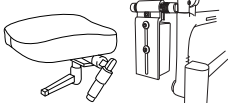				Swing-Away	Active
Seat depth		41/43/45 cm in frame (S, M, L) 37-45 cm with upholstery		●	●
Seat depth Back support 2C		41/43/45 cm in frame (S, M, L) 37-45 cm with upholstery		●	●
Back support upholstery		Adjustable upholstery		●	●
Back support cover		Airmesh (Back support 2B)		●	●
		Standard (Back support 2C)		●	●
Push handle/Push brace		Height-adjustable in 2 cm intervals. Detachable.	Push handle	●	●
			Push brace	●	●
Leg supports, detachable, swing-away		Normal knee angle		●	
		Narrow (only fits with 5" castor wheels)		●	
Leg support Lockable		Normal knee angle		●	
		Narrow		●	
Leg support Angle-adjustable		Angle adjustable with flip-up calf support from 40 cm seat width		●	
		Angle adjustable with cushion bracket and padded flip-up calf support		●	
		Angle adjustable with flip-up cushion bracket and padded calf support		●	
Leg supports, other		Amputee leg support		●	
		Amputee weight		●	




3. Chair options

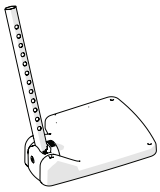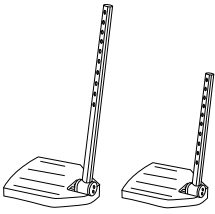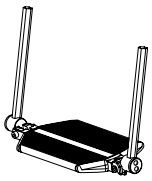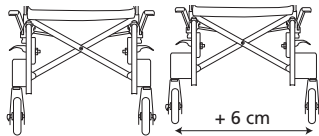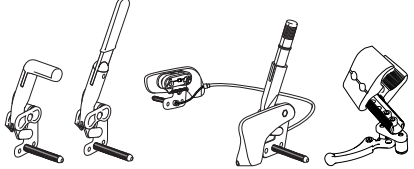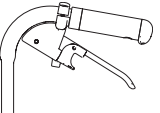				Swing-Away	Active
Foot support		15 cm deep Folding, angle and depth-settable. Attachment tubes adjustable in height	Standard attachment tubes Max. length: 51 cm. Min. length: 33 cm.	●	●
			Extended attachment tubes +10 cm	●	●
		13 cm deep. Folding, angle and depth-settable. Attachment tubes adjustable in height	Standard attachment tubes Max. length: 51 cm Min. length: 33 cm	●	●
			Extended attachment tubes +10 cm	●	●
		One-piece footplate. Attachment tubes adjustable in height		●	●
Front fork attachment		Standard		●	●
		Widened 6 cm		●	●
User brake		Angled handle		●	●
		Extended handle (foldable or fixed)		●	●
		Single hand brake		●	●
		Scissor brake With 22" drive wheels, not in positions A5–F5 and A4–B4. With 24" drive wheels, not in positions A5–B5.		●	●
Attendant brake		Brake handle 100 mm grip (only with 2C back support)		●	●
Frame colour	01 = Black 74 = Silver-white 75 = Ruby-red 76 = Champagne 77 = Menthol			●	●
Colour of rear wheel attachments	01 = Black 75 = Ruby-red 76 = Champagne 77 = Menthol			●	●

4. Accessories










Assembly instructions are provided with all accessories when they are supplied by Etac.
You can find the instructions as well as more information about the accessories on our website www.etac.com.

Settable= Adjusted using tools. Adjustable= Adjusted without tools.

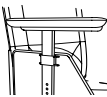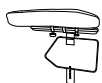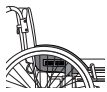			Swing-Away	Active
Arm supports	Height-adjustable 20–25 cm. Length 25 cm, solid top, black		●	●
Arm supports	Height-adjustable 20–25 cm. Length 38 cm, solid top, black		●	●
Arm support cover	Removable, 25 or 38 cm long, padded or gel, dark grey plush or black Dartex. Fits arm support with solid top.		●	●
Arm support cushion Hemi	Swing-away for fitting on long arm supports		●	●
Arm support locking Kit	A locked arm support is not detachable but its height can be set freely.		●	●
Side cover	Black		●	●
Mudguard	Right or left		●	●
Tray	Transparent, fits onto long arm supports		●	●





4. Accessories

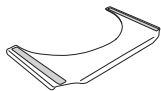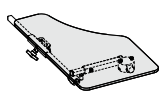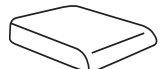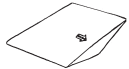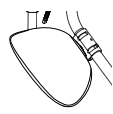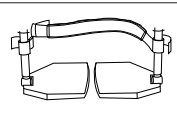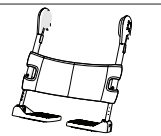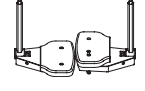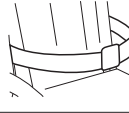			Swing-Away	Active
Hemi tray	Transparent, width 35-50 cm. Can be fitted onto short or long arm supports.		●	●
Anti-slip device for the Hemi tray			●	●
Hemi tray/Communication tray half	Right/left		●	●
Seat extender	With long seat depth		●	●
Seat cushion	Dark grey plush and black velour, 56 cm, cut according to seat depth set, washable cover		●	●
Seat cushion Tidy	Water-repellent, wipeable		●	●
Comfort wedge	For the cushion, evens out the front part of the seat		●	●
Cover for calf support	On angle-adjustable leg support, black Hygiene		●	●
Calf strap	Detachable, black nylon, length-adjustable		●	●
Padding	For calf strap		●	●
Calf pad	Plush. Fitted/adjusted with Velcro		●	●
Heel strap	Black nylon, adjustable length		●	●
Extended footplate	Fits onto the existing footplate, seat width 40-55 cm		●	●
Positioning belt	Two-piece with snap lock, fixing points on the wheelchair frame		●	●
Belt	With snaplock		●	●



4. Accessories

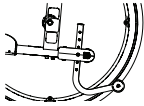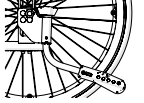			Swing-Away	Active
Head support, traditional	Adjustable in height, depth and angle. Only with Back support 2C		●	●
Fixed head support	For fitting onto push brace		●	●
Cross brace (for back support)	For Back support 2C. With snaplock		●	●
Back wedge	Three different sizes, fixed with Velcro. Creates support and shape in the back support. The wedges can be combined with other accessories and adapted to the desired size, and several wedges can be combined.		●	●
Back wedge	Dark-grey plush, 15 and 30 cm		●	●
Side cushion	Provides gentle lateral support. Used when the trunk support is too stabilising. (Combined with arm support and, if necessary, side stop)		●	●
Anti-tip	Foldable, settable in height		●	●
Anti-tip active	Settable in height, foldable from the users position. Not in combination with attendant brake		●	●
Cane holder	In two parts, one of which is an elasticated section fastened around the cane.		●	●
Spoke protectors	With grey or yellow print		●	●
Tool kit			●	●
Gloves	Information at www.etac.com		●	●




5. Installation upon delivery

Leg supports

Leg supports with foot plates are delivered unassembled and must be mounted upon delivery (see section "Handling").

Arm supports

Arm supports/side covers are delivered unassembled and must be mounted upon delivery.

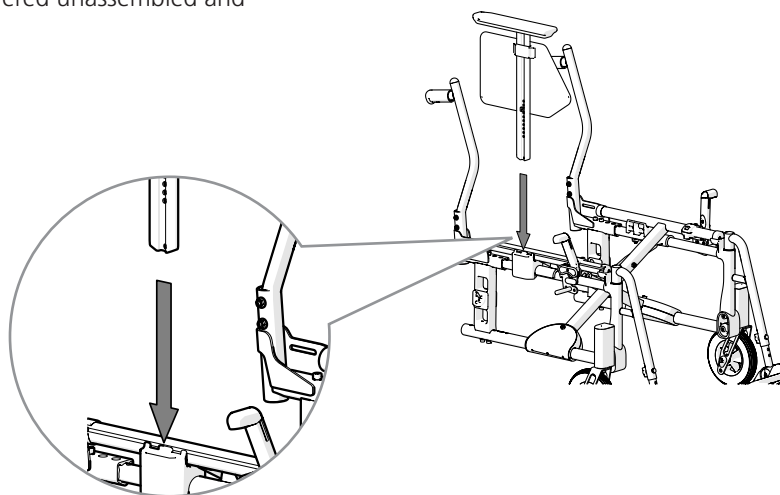

Swing-Away	Active
●	
●	●



6. Settings

Swing-Away

Active

Front seat height

Front seat height is adjustable using:

- Front fork attachment, adjustable height and angle
- Changing height in front fork.
- Changing front fork
- Changing castor wheels



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.


Seat height adjustment using front fork attachment

Loosen the screws (A and B) on both castor wheel attachments. Set the desired seat height by measuring from the floor to the top edge of the seat frame (press down the frame in the cup while measuring) and tighten the upper screw (B) on one front fork attachment. Adjust so that the attachment is perpendicular. (See 'Angle setting of front fork attachment' section.) Tighten the lower screw (A). Repeat the

procedure for the other attachment. Fine-tune the height of the front fork attachments by loosening the screws a little and carefully knocking on the front fork attachment. Make sure both castor wheels are touching the ground. Tighten the screws.



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

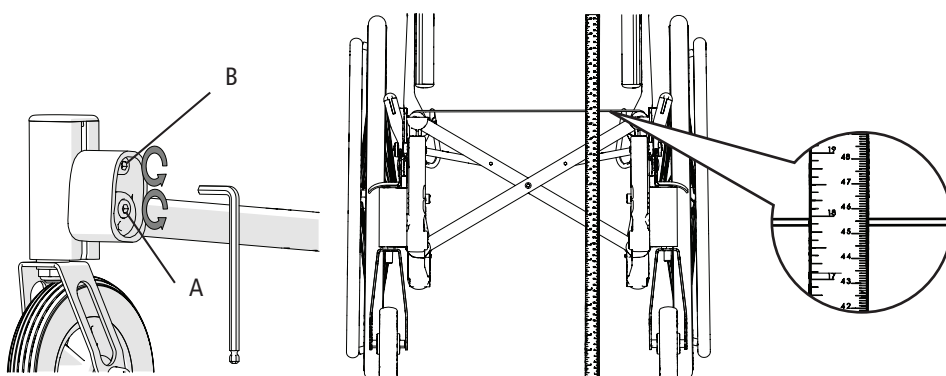

6 mm




Seat height adjustment using the castor wheels' position in the fork


See also section "Front seat height table".




Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

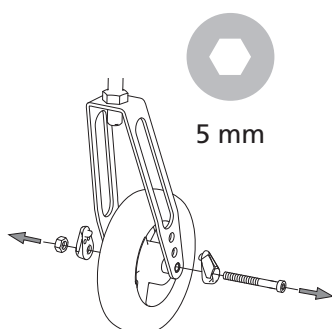

5 mm

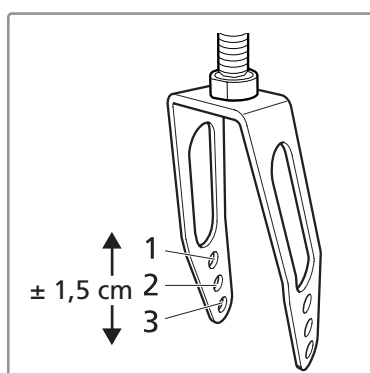


± 1,5 cm

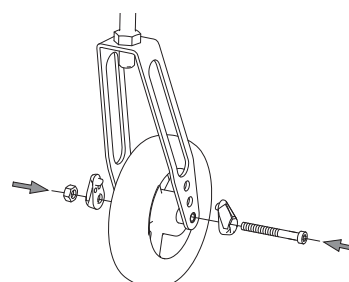

6. Settings

Disassembly/assembly of the front fork

- Remove the cover plug (A) and loosen the retaining nut (B). Pull out the front fork and transfer the bevelled washer (C) and ball bearing to the new fork.
- Mount the washers and ball bearings in the front fork attachment as illustrated. Place the washer (C) with the bevelled side towards the fork. Place the spring washer (D) uppermost under the retaining nut.
- Tighten the retaining nut until it comes to a stop and then release it by $\frac{1}{2}$ – 1 turn. The spring washer then has the correct tension, which reduces the risk of the castor starting to "wobble".
- Fit the cover plug.

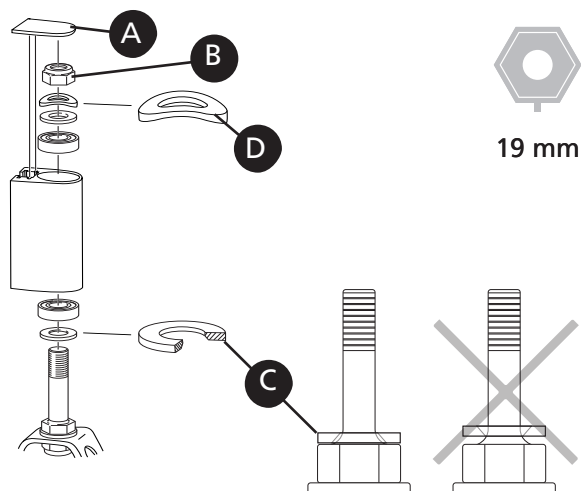

Angle setting of front fork attachment

The correct angle setting is important for the wheel-chair's manoeuvrability.

Loosen the attachment screw (A) about 2 turns so that the tag (C) moves freely inside the shaft.

Insert the Allen key in the hole (B) and turn until the attachment is 90° to the floor. Always start from a position that enables the wheel to move forward in the rotation (D).

Press in the shaft and hold it in the correct position while tightening the attachment screw.







Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.



Keep your eye on something vertical, e.g. a doorpost or a table leg when setting the angle.



6 mm

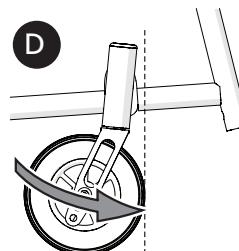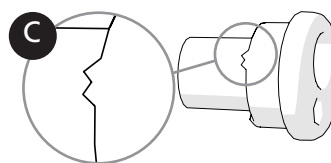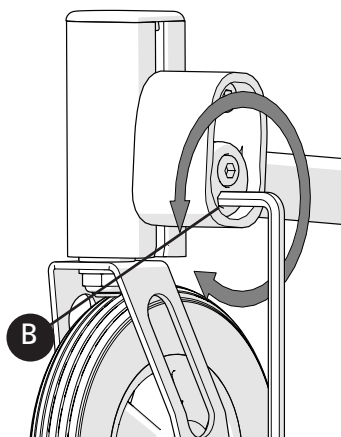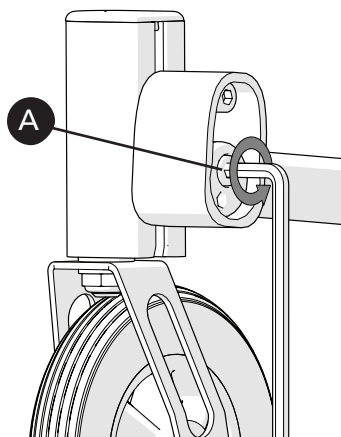

6. Settings


Rear seat height

The rear seat height can be adjusted using:

- The rear wheel's (camber washer's) position in the rear wheel attachment.
- The camber washer's position (rotation) in each rear wheel attachment opening.
- Changing the rear wheel



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

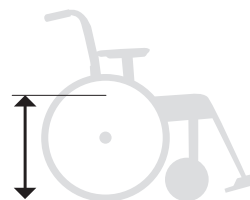


Seat height adjustment using the rear wheel position/rear wheel

1. Remove the rear wheel
2. Loosen the nut.
3. Open the kit with the washers and install at the required position (see table in section "Seat height table"). See also section "Changing camber angle").




Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

1

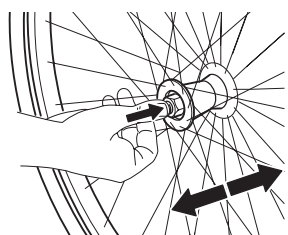


2

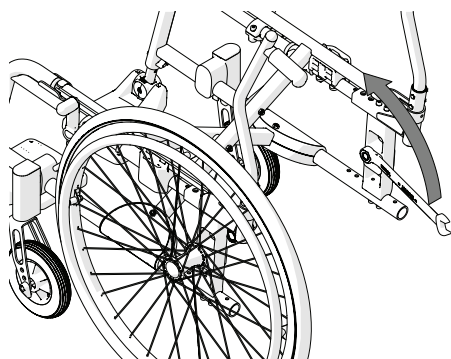


24 mm

3

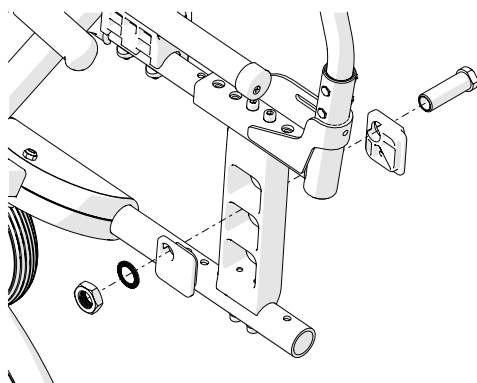

4

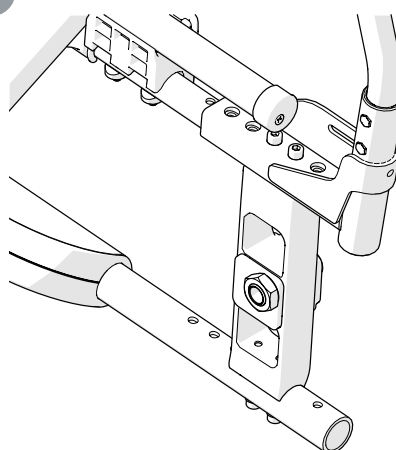

20 Nm

Swing-Away
Active

● ●





● ●

6. Settings


Swing-Away	Active
☐	☐



When fitting, check that the hub axle of the rear wheel is set correctly. The button in the hub should pop out to the original position.

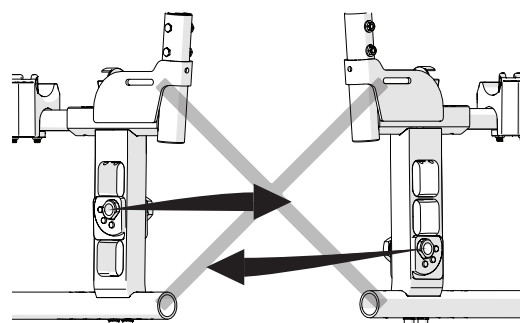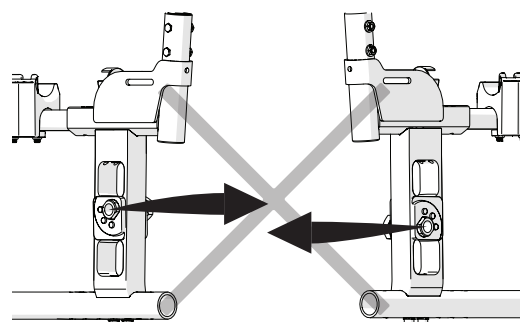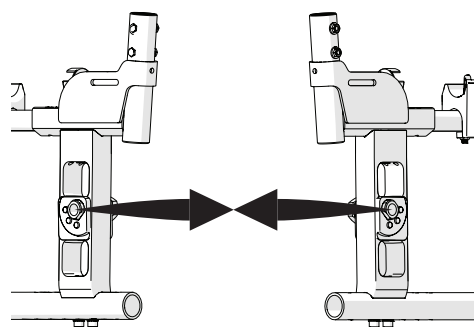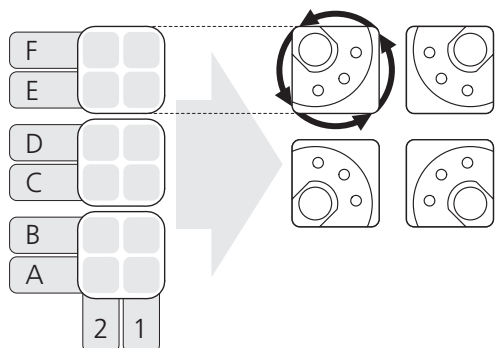


Fitting 0° Camber

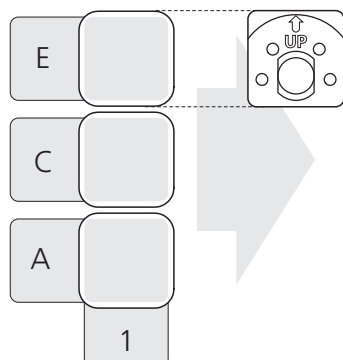


6. Settings

Rear seat height table


Rear wheel position height

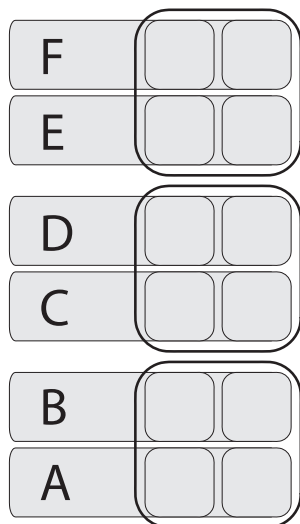

Seat height rear

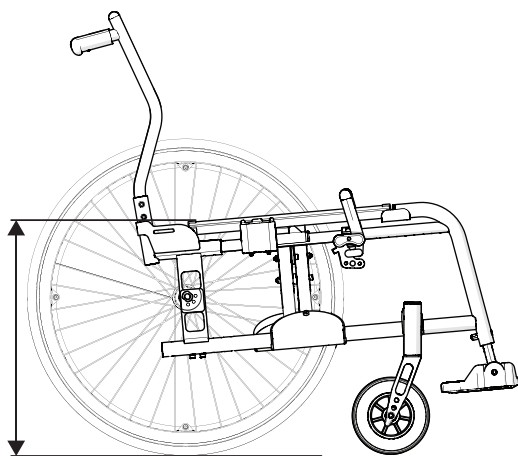

Distance rear wheel to seat

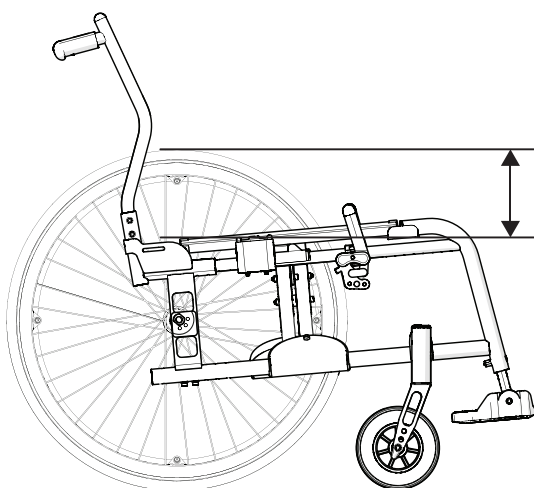














	Rear wheel	Seat height	Distance	Swing-Away	Active
Position A	22"	48.5 cm	7.0 cm	●	●
	24"	51 cm	9.5 cm	●	●
Position B	22"	46.5 cm	9.0 cm	●	●
	24"	49 cm	11.5 cm	●	●
Position C	22"	44.5 cm	11.0 cm	●	●
	24"	47 cm	13.5 cm	●	●
Position D	22"	42.5 cm	13.0 cm	●	●
	24"	45 cm	15.5 cm	●	●
Position E	22"	40.5 cm	15.0 cm	●	●
	24"	43 cm	17.5 cm	●	●
Position F	22"	38.5 cm	17.0 cm	●	●
	24"	41 cm	19.5 cm	●	●

6. Settings

Swing-Away


Active

Attendant brake

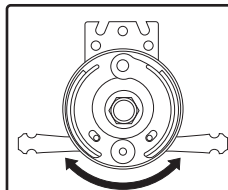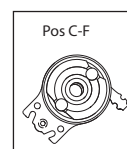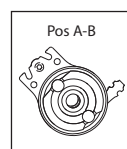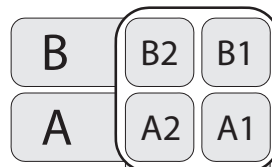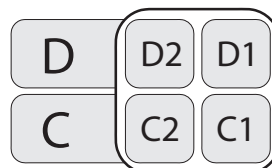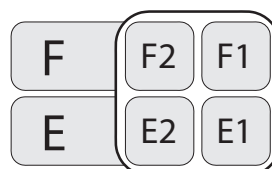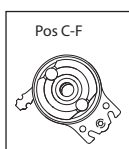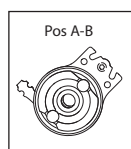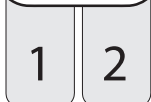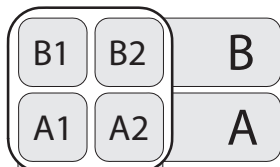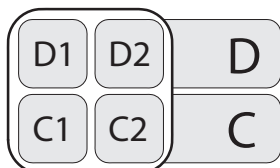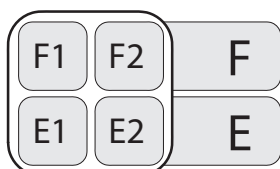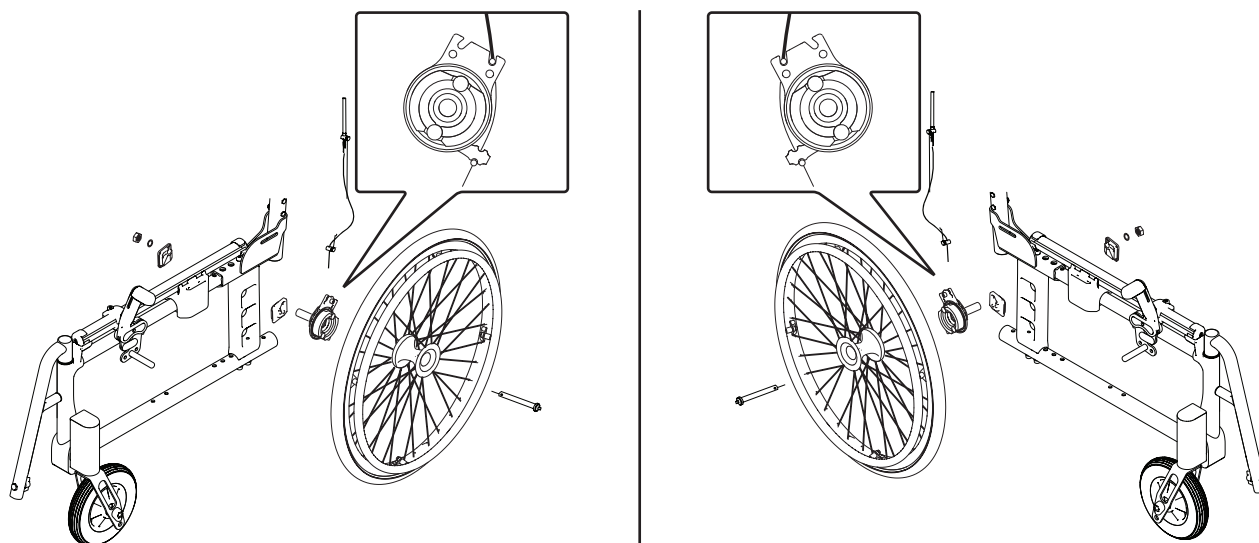

The brake arm can be turned.


When adjusting the brake loosen nut (A) and adjust the brake by screwing (B) up/down. When the brake has a reasonable effect tighten nut (A). This adjustment should be carried out when the brake is not functioning satisfactorily.

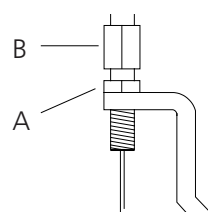


8 mm

User manual


75369G 21-05-10

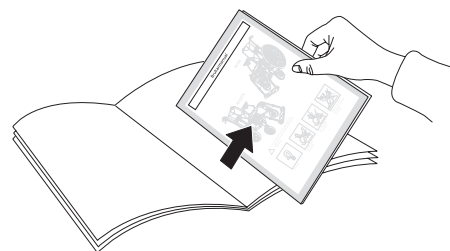

This part of the manual (User Manual) can be removed and supplied with the wheel-chair as it contains important user information.

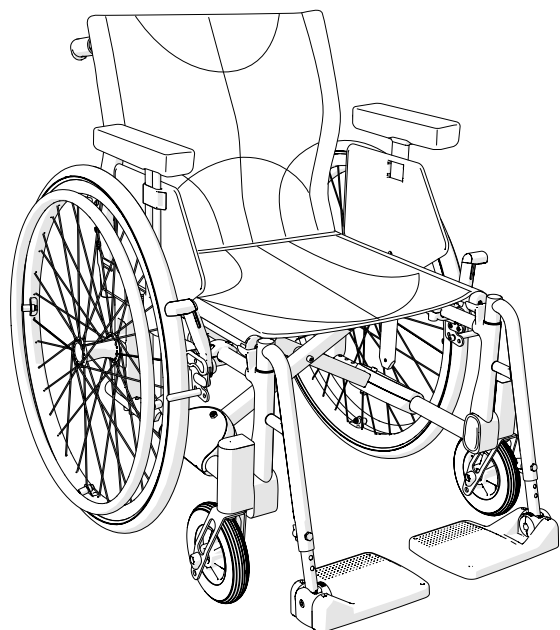

Etac Crissy Swing-Away

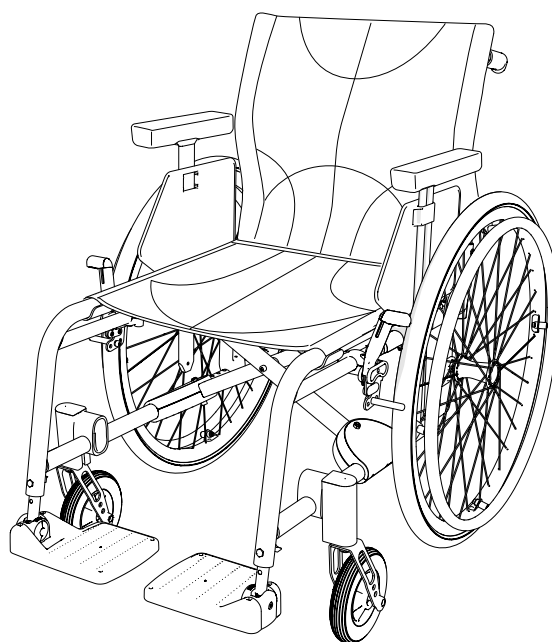

Etac Crissy Active



User manual

Table of Contents

User manual

Symbols and warnings	20	Driving technique, use	26-29
Transportation in vehicles	21	Care and maintenance.....	30
How to use	22-23	Troubleshooting chart.....	30
Seating comfort settings			
Back support cover	24		
Seat cushion	24		
Back support upholstery	24		
Back angling.....	24		
Footplates	25		
Calf/heel strap	25		
Arm supports	25		


Symbols and warnings



Removable part, do not lift!




Risk of tipping: Always have the anti-tips down when in motion.

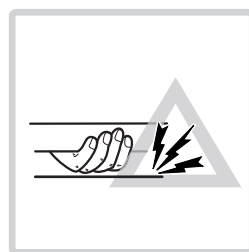

Risk of crushing

General warnings


The following must be taken into account when using the product.

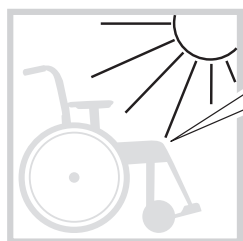

Bright sunlight could heat parts of the chair and cause burns



Do not lift the chair by removable parts.

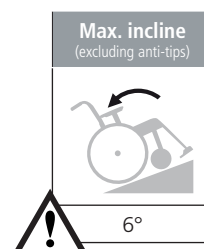

The chair may tip over if the backward tilt is greater than 6°. Use anti-tips!



After adjusting the seat height, centre of balance or back support angle, always ensure that you check the function of the anti-tips.

User manual


Transportation in vehicles

Etac recommends in the following order:

- 1) The user transfers to one of the seats in the vehicle and uses the vehicle's 3-point belt while travelling. The wheelchair is then placed in the boot or safely in the back seat so that it cannot overturn or roll.



- It is best to place the wheelchair in the car boot.
- If the wheelchair is placed on the back seat, ensure that it cannot overturn or roll. If possible secure the wheelchair with the car's safety belt.

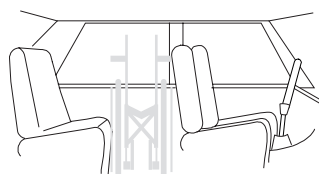

- 2) The wheelchair is secured facing forwards in the vehicle as per this manual; the user uses a separate 3-point belt secured in the vehicle. This is the way in which the wheelchair is tested and approved according to ISO 7176-19.







Securing

- Rear transport attachment must be used
- The restraining device must not be passed through the wheels or around the back tubes.



- The parking brake must be used
- The anti-tips must be activated
- Accessories that can be removed without tools must be removed and secured
- A correctly set head support must be used
- The back support should be level with or above the user's shoulders
- Handrims must be installed on the rear wheels

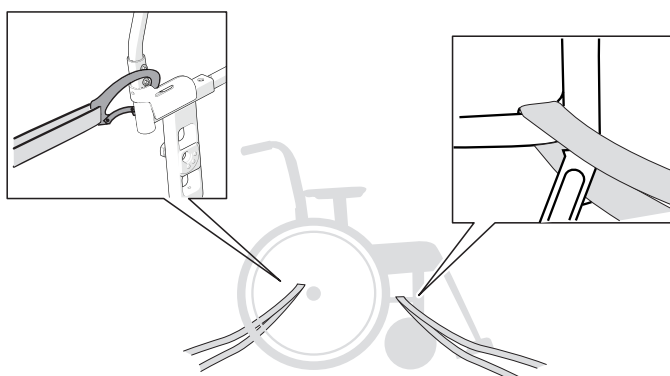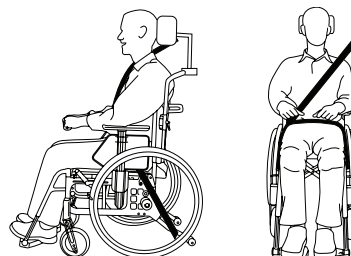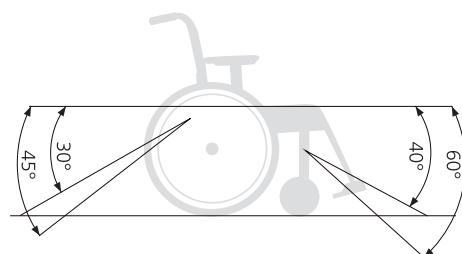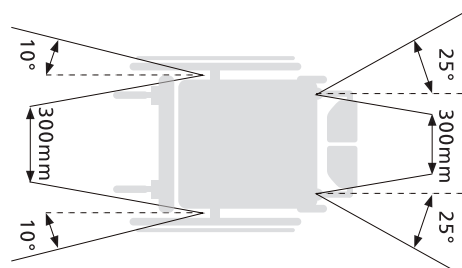

- 3) According to Directive 2001/85/EC, appendix VII, point 3.8.3, there may be specially marked wheelchair locations in vehicles that permit transport with a wheelchair facing in the direction of travel. If this means of travel is used, the user/carer must be aware while travelling, prepared for sudden movements and have the capacity to maintain a safe sitting position throughout the entire journey. The user's disabilities must not be of such an extent that he/she is not able to hold onto the handles fitted in the vehicle when there are changes of speed or direction.



- The user brake must be used
- The anti-tips must be activated
- Accessories that can be removed without tools must be removed and secured.
- A correctly set head support must be used
- The back support should be level with or above the user's shoulders
- A positioning belt should be used


User manual

How to use




Unfolding

- Push down on one side of the seat frame, using the whole of the flat of the hand.
- Do not hold the seat frame tube while unfolding as there is a risk of pinching your fingers.
- Flip down the foot supports.

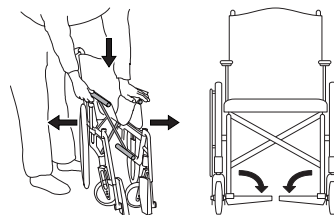

Folding


- If a push brace and/or cross brace are fitted they should be removed.
- Flip up the foot supports.
- Lift the seat upwards.

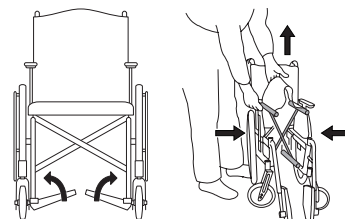

Rear wheel with quick release hub




During assembly, check that the quick release button pops out completely. Test by pulling the wheel.

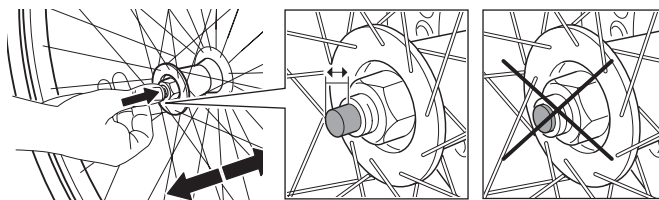


Anti-tip


- The anti-tips can be folded away/deactivated.

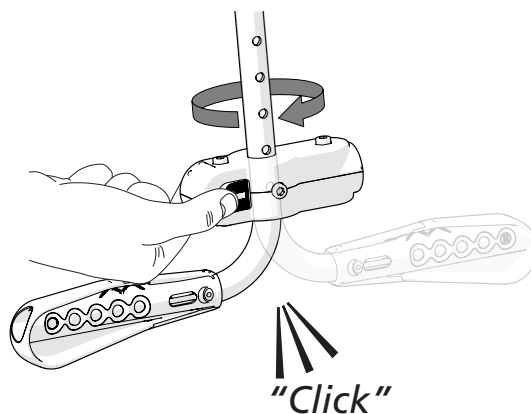

- When you activate the anti-tip, ensure that the locking pins are properly engaged.

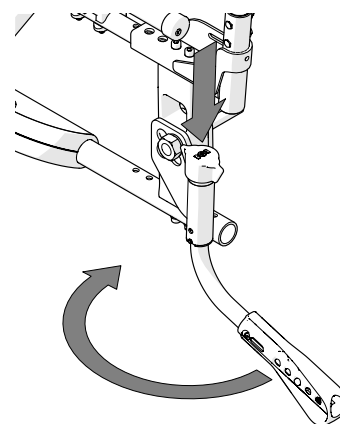

- After all adjustments to the chair, the anti-tip function must be checked. If the anti-tips need to be adjusted, contact your technical aids centre.

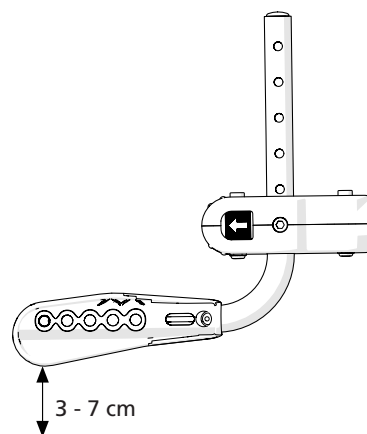

User manual


Leg support (Etac Crissy Swing-Away)

The leg supports can be lifted/swung to the side as needed. lockable leg supports by pressing the red button (A).

To attach the lockable leg supports, press down on the red button (A), put the leg support into place and release the button. The leg support is now locked (B).



Check that the leg supports have locked before lifting.

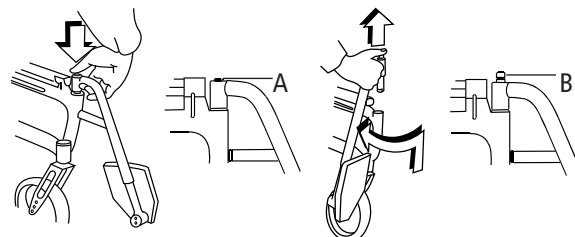

Push handle/push brace, height-adjustable

Loosen knob (A). Press in the red wedge and set the desired height. Ensure that the push handle is in a "snap lock" position, and tighten the knob.





Make sure knobs are securely tightened. This is especially important if the chair is to be lifted with the user sitting in it.



Ensure that the push handle's button lock has popped out below the push handle attachment.




The push brace must be removed before the chair can be folded.

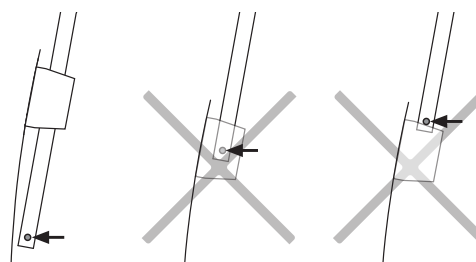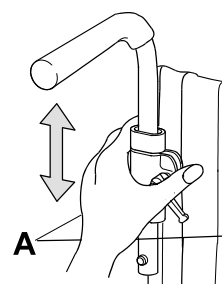


Attendant brake (motion brake)

- Brake while in motion by pulling the brake lever upwards.

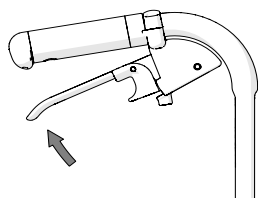

Attendant brake (parking brake)

- Activate the parking brake by pulling the brake lever upwards and pressing in button A.
- Release by pulling the lever upwards.

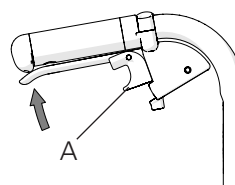

User manual

Seating comfort settings

Seating comfort

Good seating comfort depends on the needs and possibilities available to each individual. The following general guidelines here show what to consider when settings are made.


Select a seat cushion with care. Aside from comfort, it has an impact on stability and thereby freedom of movement. Seat cushions also have

different pressure distribution qualities.

The back support setting is a combination of angles, height and shape. This is why the angle settings and height for the back support may need to be adjusted when the back support cover is being adjusted.


Back support cover

- Loosen all back straps, but make sure that the Velcro is still stuck together.
- Place the back support cover so that it covers the upper edge of the back rails and forms a crease between the seat and the back support, allowing the user to properly sit in on the seat.

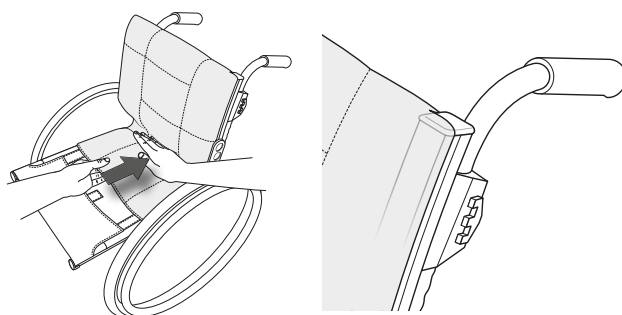

Seat cushion


- Make sure that the correct seat depth is set.
- Position the cushion as far back on the seat so that the rounded corners end up between the back rails. Make sure the cushion sticks to the Velcro.
- Measure a suitable seat cushion length and cut it according to supplied instructions.
- A more contoured seat cushion can be obtained by using a comfort wedge which is placed in the cushion cover's pocket under the seat cushion.

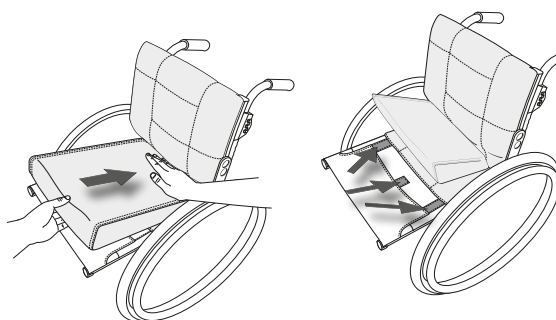

Back support upholstery

The back support upholstery is adjusted with the user sitting in the wheelchair.

- Make sure the user is sitting as far back in the seat as possible.
- Support the pelvis by tightening the strap found just below the lumbar region.
- Then adjust the upper strap to provide the user with upper body support and balance.
- The other straps are adjusted to provide space for the user's posterior and follow the back's natural curvature.

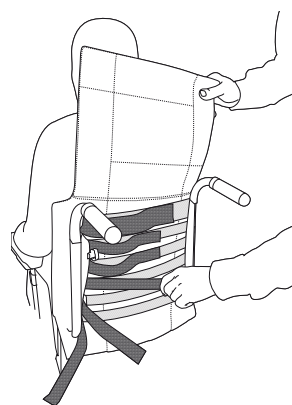

Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.



Do not over-tighten the upper Velcro straps as this can prevent the cross-brace from unfolding properly, i.e. the back tubes do not sit flush in the side-frames.



User manual

Seating comfort settings


Footplates

Adjust the height so that you feel the feet are supported and the thighs rest against the cushion.

For outdoor use the footplates should be raised 4 to 5 cm above the ground.




Tipping risk!
Never stand on the footplates

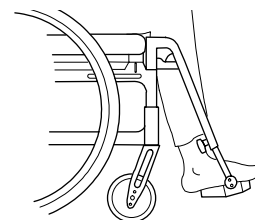

Angle adjustment





Adjust the angle of the footplates so that the ankles are at a 90° angle.

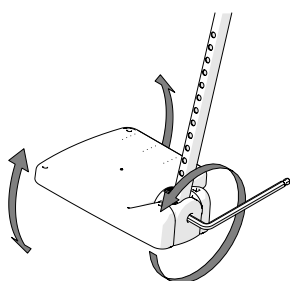

5 mm



Adjusting the height

Fully loosen the locking handle on the leg support. Remove the screw. Set the height. Replace the screw and handle. Tighten securely.

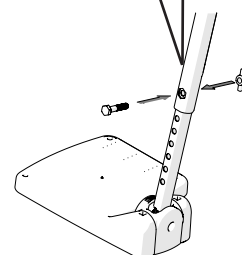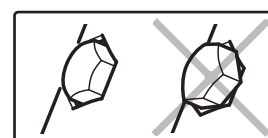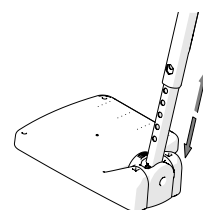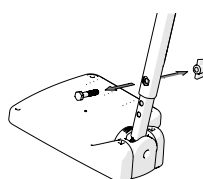


Changing the depth of the foot supports

Adjust the depth of the foot supports by switching the right and left foot support. Do not forget to adjust the angle.

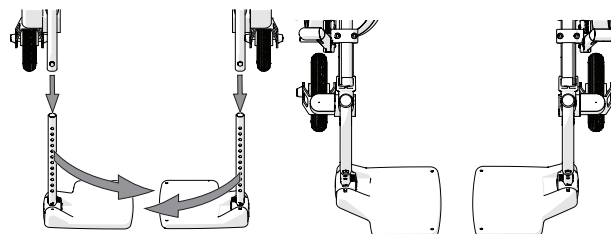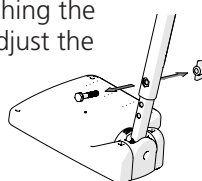

Calf/Heel strap

- Adjust the length so that the feet are centred on the footplates.

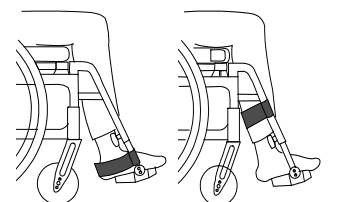


Arm support (adjustment)

- Loosen the screws on the inside of the arm support tray.
- Move the arm support up or down to the desired height and tighten the screws.
- Correctly adjusted arm supports provide upper body support. The arm support plate should be approximately 5 mm higher than the elbow when bent at an angle of 90°.



3 mm

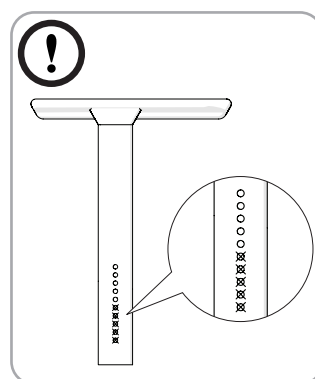

Only the top five holes may be used (applies to arm support height-adjustable 20–25 cm).

User manual

Driving technique, use




Parking

Turn the castor wheels forwards and activate the brakes.



Do not park the wheelchair on a hill when it has a user in it.

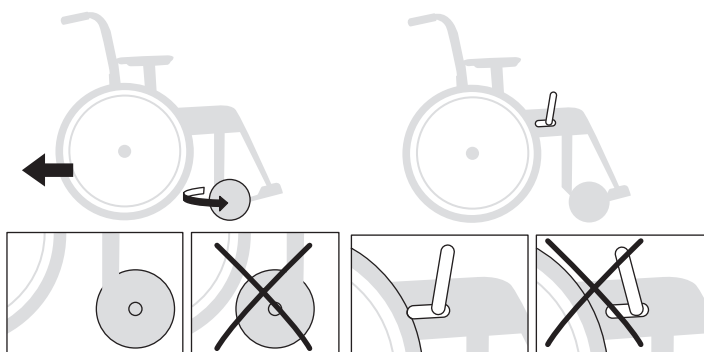


Transferring to/from the wheelchair



Turn the castor wheels forwards, activate the brakes, flip up the foot supports and twist the leg supports to one side.




Never stand on the footplates as there is a risk of tipping!

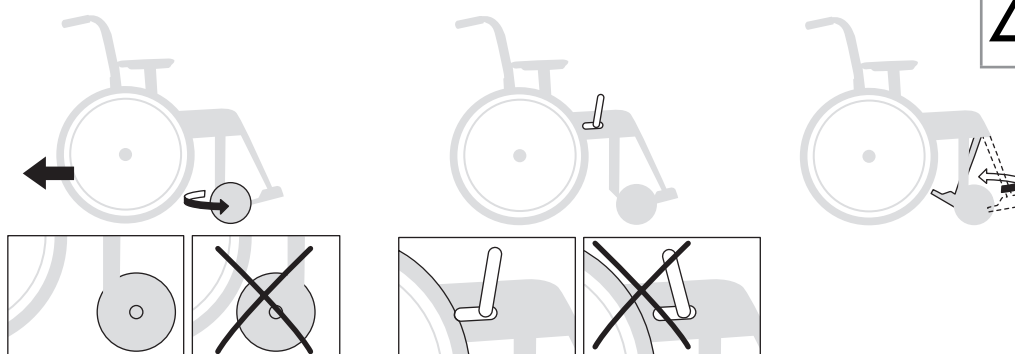



Sideways

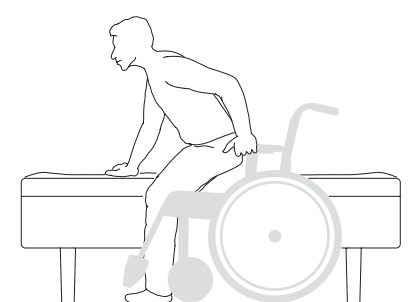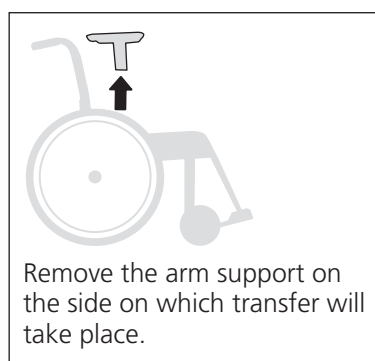


From the front

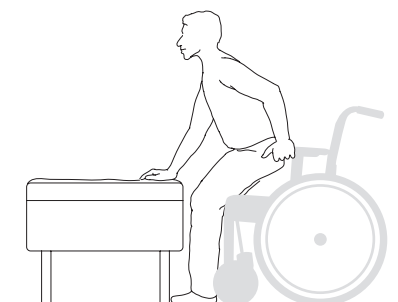

Lifting the wheelchair

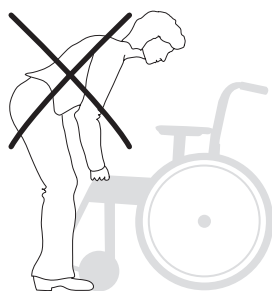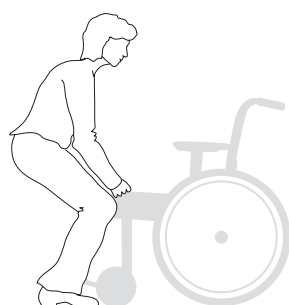

- Ensure that vertically adjustable push handles are properly tightened and that lockable leg supports are properly locked.
- Swing non-lockable leg supports in under the seat and lift the wheelchair by the front upper part of the frame.

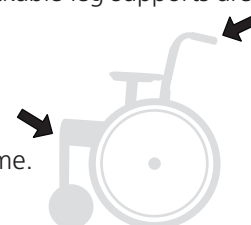

User manual


 The following illustrations show principles of driving technique and obstacle clearance.



User, up forwards

This technique is recommended only for experienced wheelchair users.

– Ensure that the anti-tips are disabled.

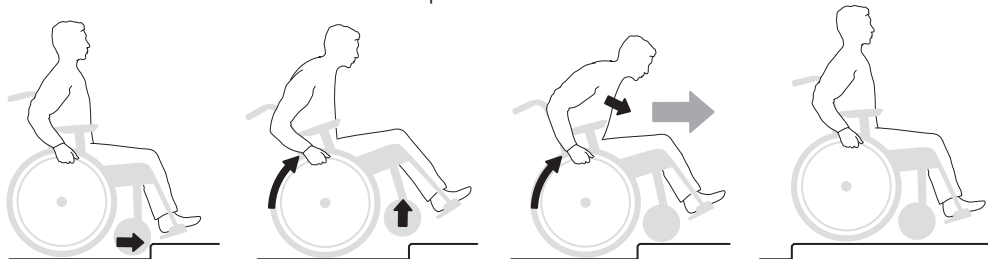


Enable the anti-tips afterwards.



User, up backwards

This technique only works if there is a low kerb/threshold, relative to the installed height of the foot supports.

– Ensure that the anti-tips are disabled.

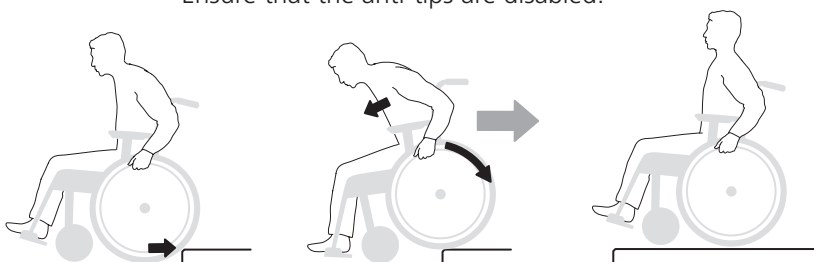


Enable the anti-tips afterwards.



User, down forwards

This technique is recommended only for experienced wheelchair users.

– Ensure that the anti-tips are disabled.

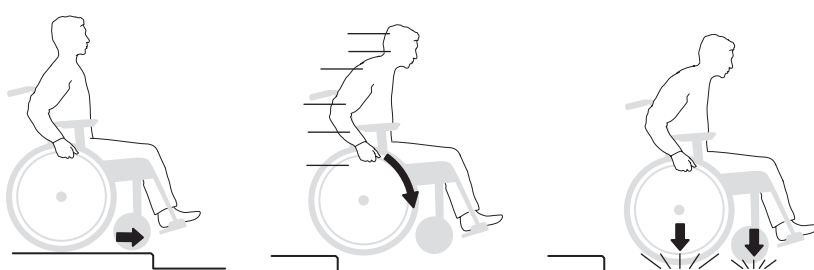


Enable the anti-tips afterwards.




User, down backwards

This technique is recommended only for experienced wheelchair users.

– Ensure that the anti-tips are disabled.

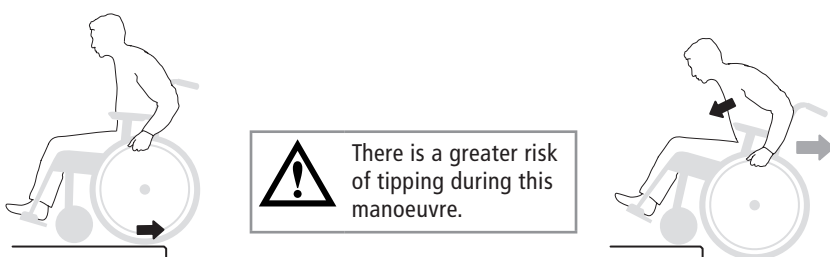

There is a greater risk of tipping during this manoeuvre.

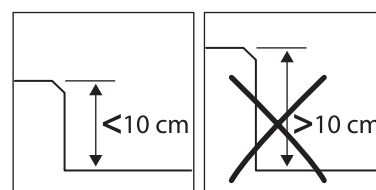


Enable the anti-tips afterwards.

User manual



Carer, up forwards

– Ensure that the anti-tips are disabled.

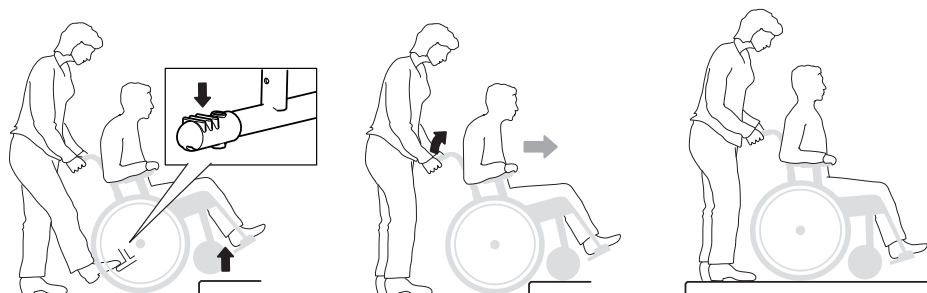


Enable the anti-tips afterwards.



Carer, up backwards

– Ensure that the anti-tips are disabled.

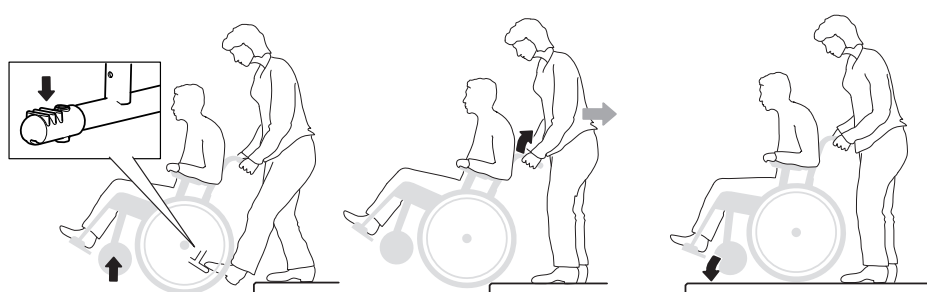


Enable the anti-tips afterwards.



Carer, down forwards

– Ensure that the anti-tips are disabled.

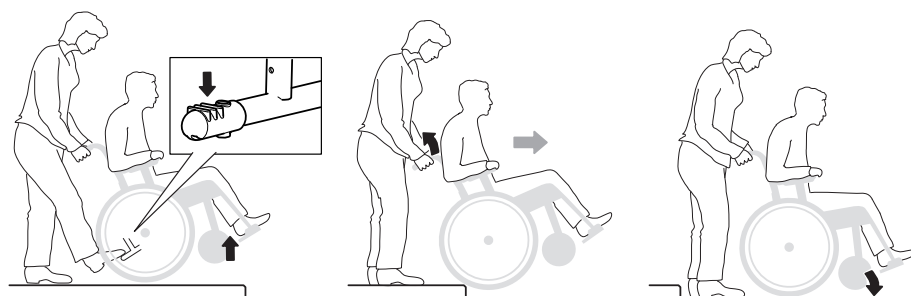


Enable the anti-tips afterwards.



Carer, down backwards

– Ensure that the anti-tips are disabled.

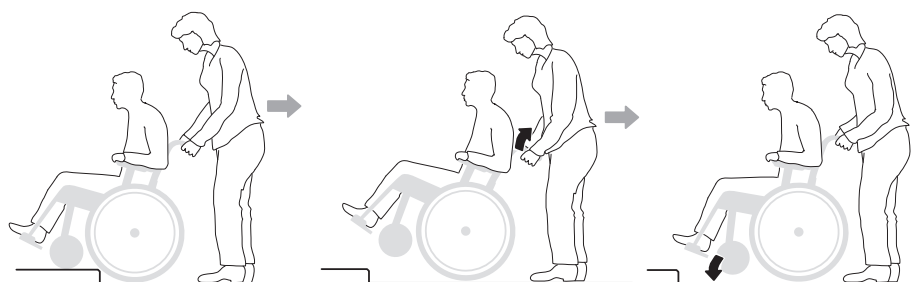

Enable the anti-tips afterwards.






User manual



Stairs, up

– Ensure that the anti-tips are disabled. Never use an escalator, even if a carer is available. Ensure that vertically adjustable handles are tightened.

We recommend that there are always two carers for this transfer: one walking behind and holding the push handles, and one walking in front, holding the frame (or the leg supports if these are lockable).

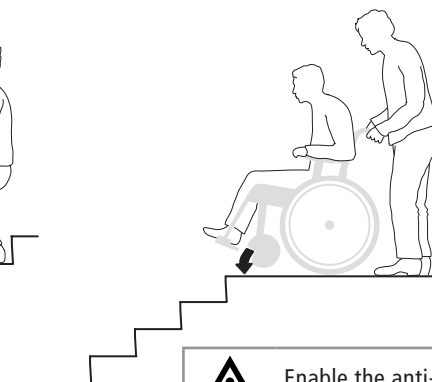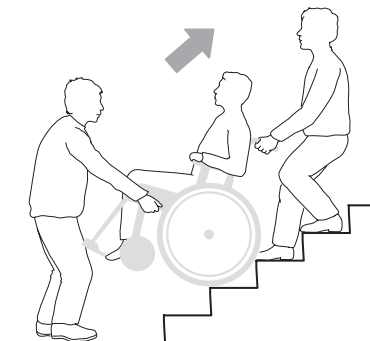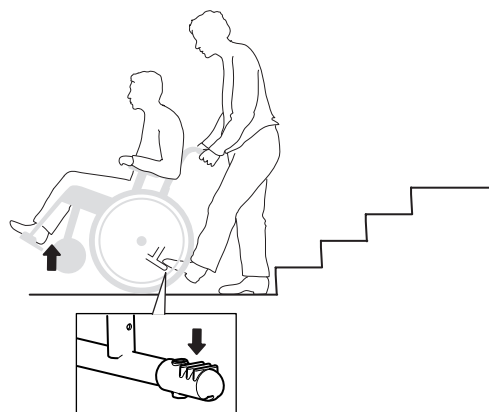




Enable the anti-tips afterwards.



Stairs, down

– Ensure that the anti-tips are disabled. Never use an escalator, even if a carer is available. Ensure that vertically adjustable handles are tightened.

We recommend that there are always two carers for this transfer: one walking behind and holding the push handles, and one walking in front, holding the frame (or the leg supports if these are lockable).

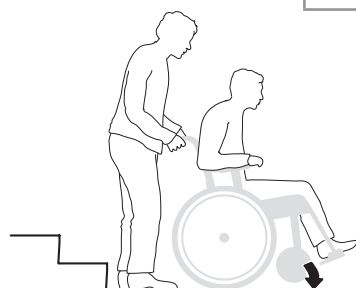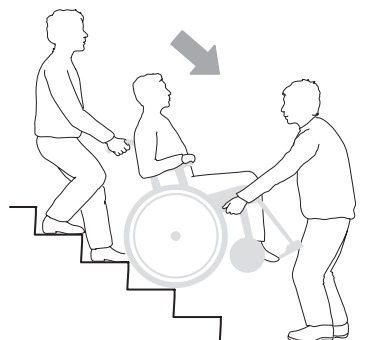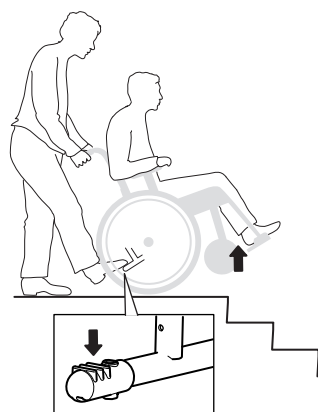

Enable the anti-tips afterwards.




Ascending/descending hills

Control your speed using the handrims, not the brakes!

Avoid turning round in the middle of a hill. Always drive as straight up/down as possible. If uncertain, ask for assistance.


Uphill

Lean forwards to correct your centre of balance.

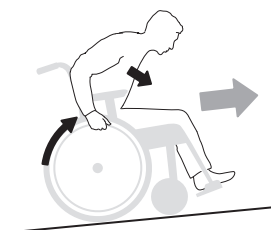

Downhill

Lean backwards to correct your centre of balance.

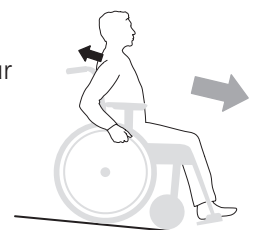

Enable the anti-tips afterwards.

User manual

Care and maintenance

Rear Wheels: Clean the quick release axles when necessary.

Castor wheels: Clean the castor wheel axles when necessary.

Frame: Clean the frame with a non-abrasive cleaning agent with a pH level between 5 and 9, or with a 70% disinfectant solution. Rinse and dry.

Upholstery: Wash, see label on upholstery.

For more severe problems, contact the Technical Aids Centre.



When necessary lubricate moving parts/joints with bicycle oil or similar.

Troubleshooting chart

The wheelchair pulls to the side	• Inflate the tyres (air pressure, see section "Technical data") • Check and adjust the height and angle settings of the front fork attachments • Check that the rear wheel attachments are correctly mounted • The user is distributing weight unevenly in the wheelchair • More strength being used on one side than the other
The wheelchair is "heavy" to push	• Inflate the tyres (air pressure, see section "Technical data") • Check that the rear wheel attachments are correctly mounted • Clean the castor axles of hair and dirt • Too much weight over the castors. Adjust the wheelchair's centre of balance
The wheelchair is "heavy" to turn	• Inflate the tyres (air pressure, see section "Technical data") • Check that the front forks have not been tightened too much • Clean the castor axles of hair and dirt • Too much weight over the castors. Adjust the wheelchair's centre of balance
Brakes not effective	• Inflate the tyres (air pressure, see section "Technical data") • Adjust the distance between brake and tyre
Rear wheels are "loose"	• Adjust the length of the axle shaft
Rear wheels hard to remove/re-place	• Clean and lubricate the quick release using e.g. bicycle oil • Adjust the length of the axle shaft
The castor wheels "wobble"	• The front forks have not been tightened sufficiently • Check and adjust the height and angle settings of the front fork attachments • Too much weight over the castors. Adjust the wheelchair's centre of balance
The wheelchair is hard to fold/unfold	• The upholstery is too tight • Clean and lubricate the cross-brace under the seat
The wheelchair feels "awkward"	• Inflate the tyres (air pressure, see section "Technical data") • Check that screws and controls are properly tightened



6. Settings

Seat angle


The seat angle is achieved through the difference between the front and the rear seat heights.



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

When the seat angle is changed:

- Adjust the angle of the front fork attachment.

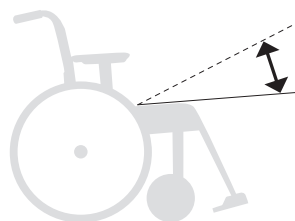

Seat depth




The seat depth can be adjusted by moving the attachment of the back support in the frame and lifting the front seat upholstery and sliding it backwards or forwards.

To move the back attachment, loosen the screws for back attachment/wheel bracket and move to any position.

NB. Maximum seat depth setting is 45 cm using the back attachment.



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

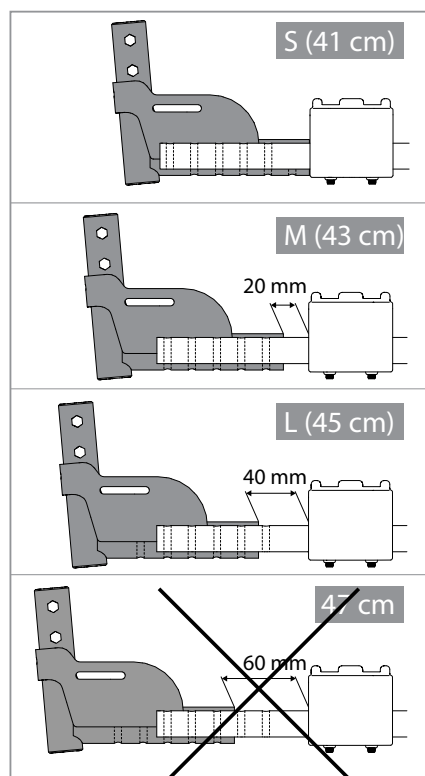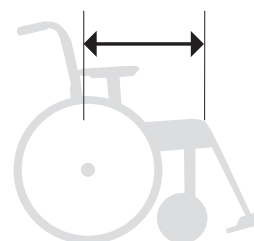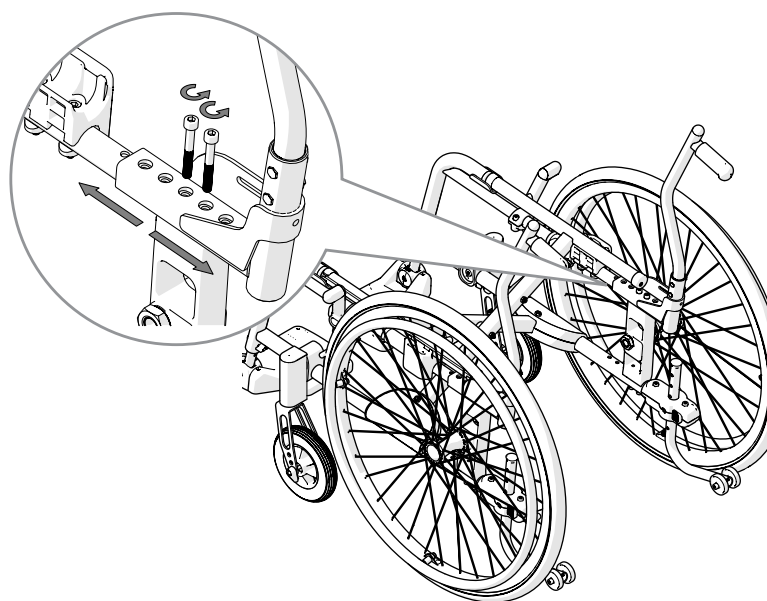

Swing-Away

Active

●

●

●

●


6. Settings

Seating comfort

After setting seat height, seat angle, seat depth and back angle, the back support upholstery should be adjusted to give support and comfort. See section “*Seating comfort settings* ” in the User Manual.

Adjustment of the centre of balance

Available positions for wheel mount at various seat depth settings.

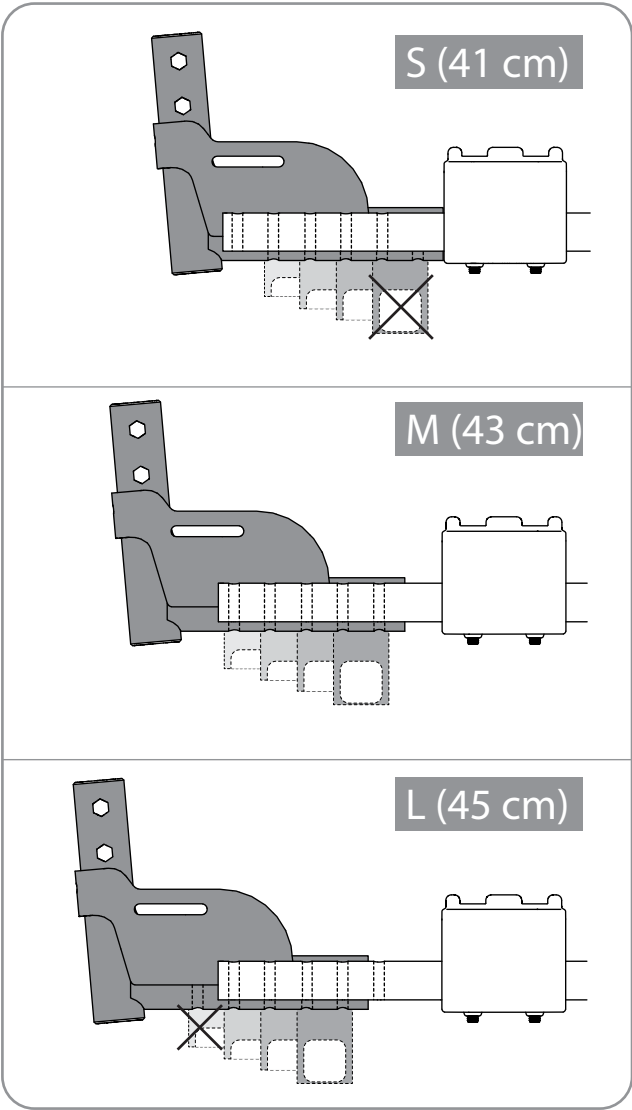

Swing-Away	Active
●	●
●	●



6. Settings

Adjustment of the centre of balance

The centre of balance can be adjusted by changing the position of the rear wheels.

- Move the rear wheel attachment forwards or backwards or rotate the camber washer.



Also check the brakes - see section "Brake, adjustment".



24 mm

When the position of the rear wheels is moved forwards the wheelchair becomes much more manoeuvrable, but the tendency to tip backwards increases. See section "Seat height adjustment using the rear wheel position/rear wheel" and section "Seat depth".



The centre of balance is also altered when the seat and/or back support angle are adjusted.









We recommend the use of anti-tips.



Ensure that rear wheels with quick release hubs are mounted securely.



12 Nm

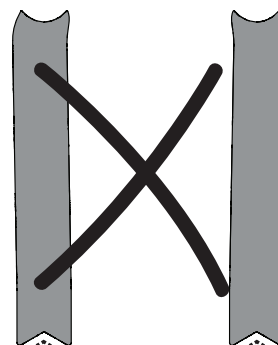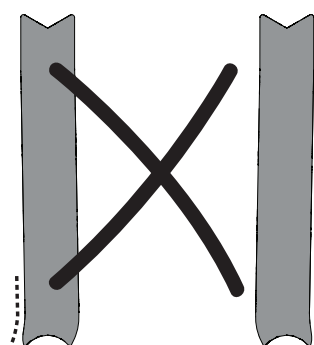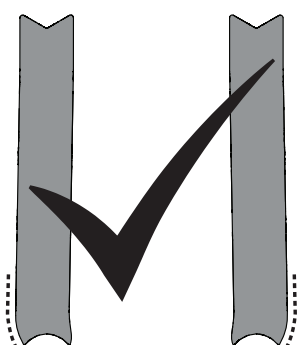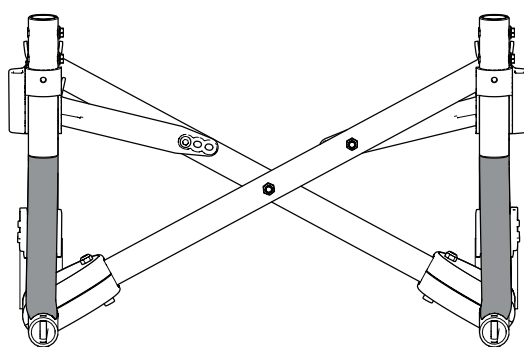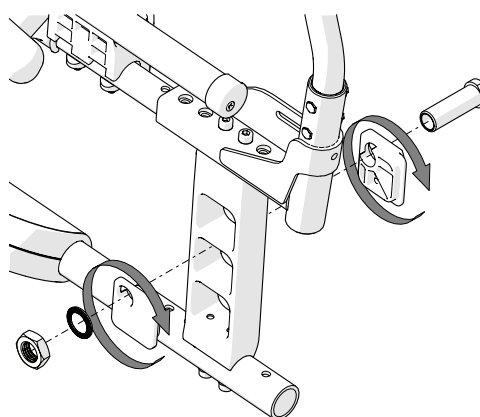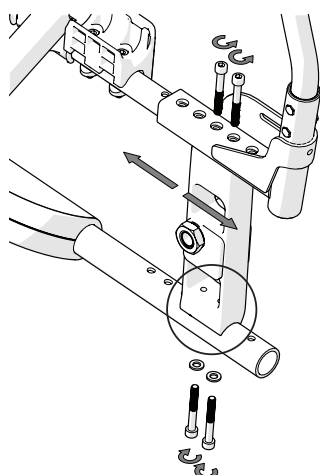


7. Fitting seat extender



Swing-Away

Active

Fitting seat extender

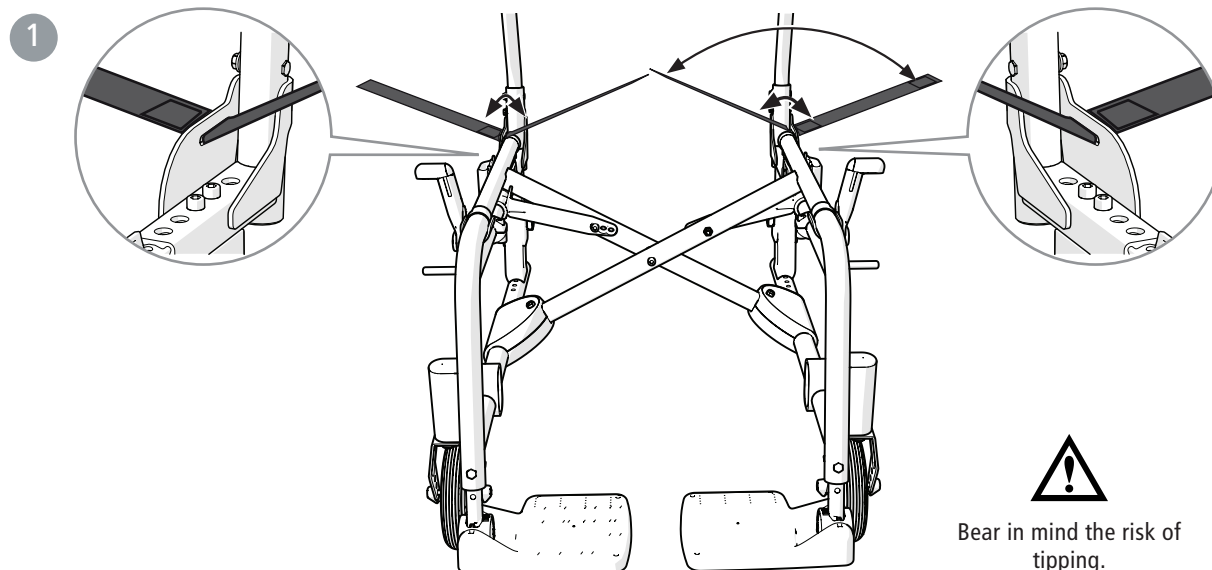

Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

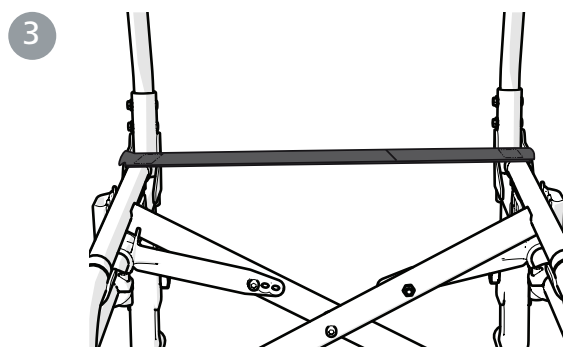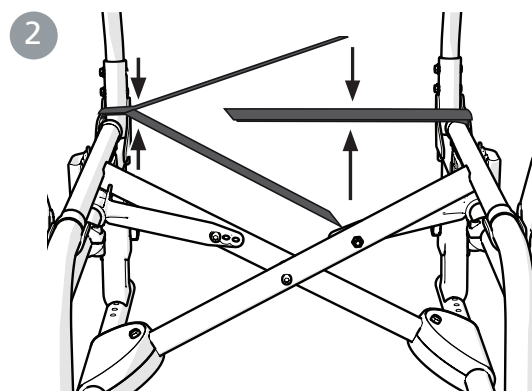


9. Settings Back support 2B

Adjusting back support height

The back support height is adjustable (See "*Chair options*"). Loosen the screws to adjust the height of the back support. Tighten the screws.



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

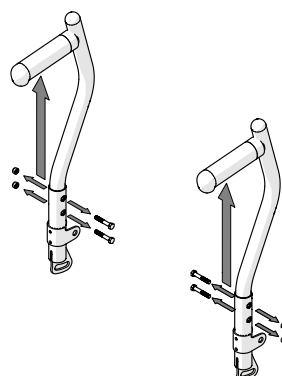

8 mm




3 Nm

Adjustment of back support upholstery

The contour of the back support upholstery is individually adjustable by using the five Velcro straps and the back support cover.

- Allow the cover plenty of room between seat and back support, so that it is possible to "sit in" against the back support.
- Loosen all the straps and ensure that the user is sitting as far back in the seat as possible.
- Tighten the straps so that they follow the contours of the back and give support to the lumbar region.

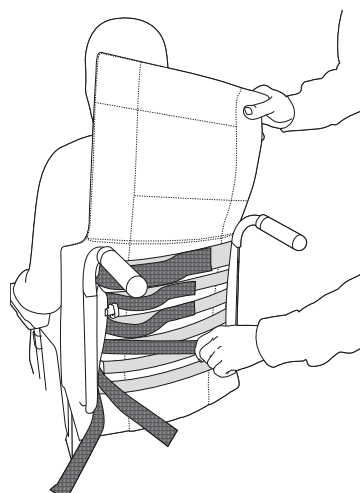

Do not over-tighten the upper Velcro straps as this can prevent the cross-brace from unfolding properly, i.e. the back tubes do not sit flush in the side-frames.




Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

Justering av ryggvinkel

Angle adjustable back support is optional (see "*The chairs possibilities*"). Adjust the angle by loosening the lower screw.




Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

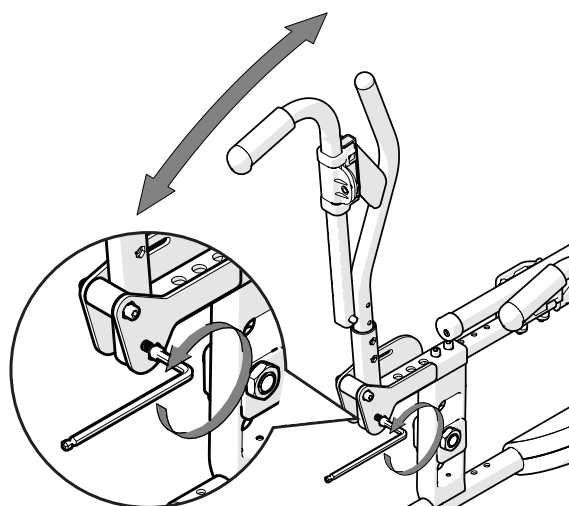

8. Settings Back support 2C

Adjusting back support height

The back support height is variably adjustable (See "*Chair options*"). Loosen the screws (A) 2-3 turns to set the height of the back support. Ensure that both sides are at the same height by checking against the ruler (B) on the back of the back tube. Tighten the screws securely.

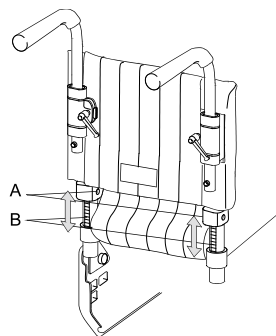

5 mm




Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.

Adjustment of back support upholstery

The contour of the back support upholstery is individually adjustable by using the five Velcro straps and the back support cover.

- Allow the cover plenty of room between seat and back support, so that it is possible to "sit in" against the back support.
- Loosen all the straps and ensure that the user is sitting as far back in the seat as possible.
- Tighten the straps so that they follow the contours of the back and give support to the lumbar region.

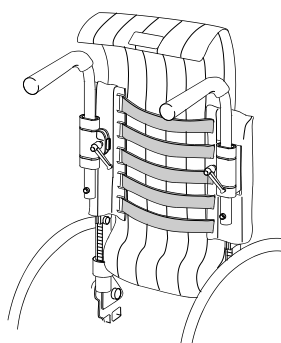

Do not over-tighten the upper Velcro straps as this can prevent the cross-brace from unfolding properly, i.e. the back tubes do not sit flush in the side-frames.



Bear in mind the risk of tipping.
Always check the positioning of the anti-tips.



10. Readjustments

Brake, adjustment


Swing-Away	Active
●	●

The brakes are variably adjustable.

1. Loosen the screws about one turn.
2. Adjust the brake position.
3. The brake block must be approximately 15-20 mm from the tyre when the brake is not being applied.
4. Re-tighten the screws. Test the brakes!




The braking effect is dependent on the air pressure in the tyres. The brakes are parking brakes and must not be used as motion brakes.

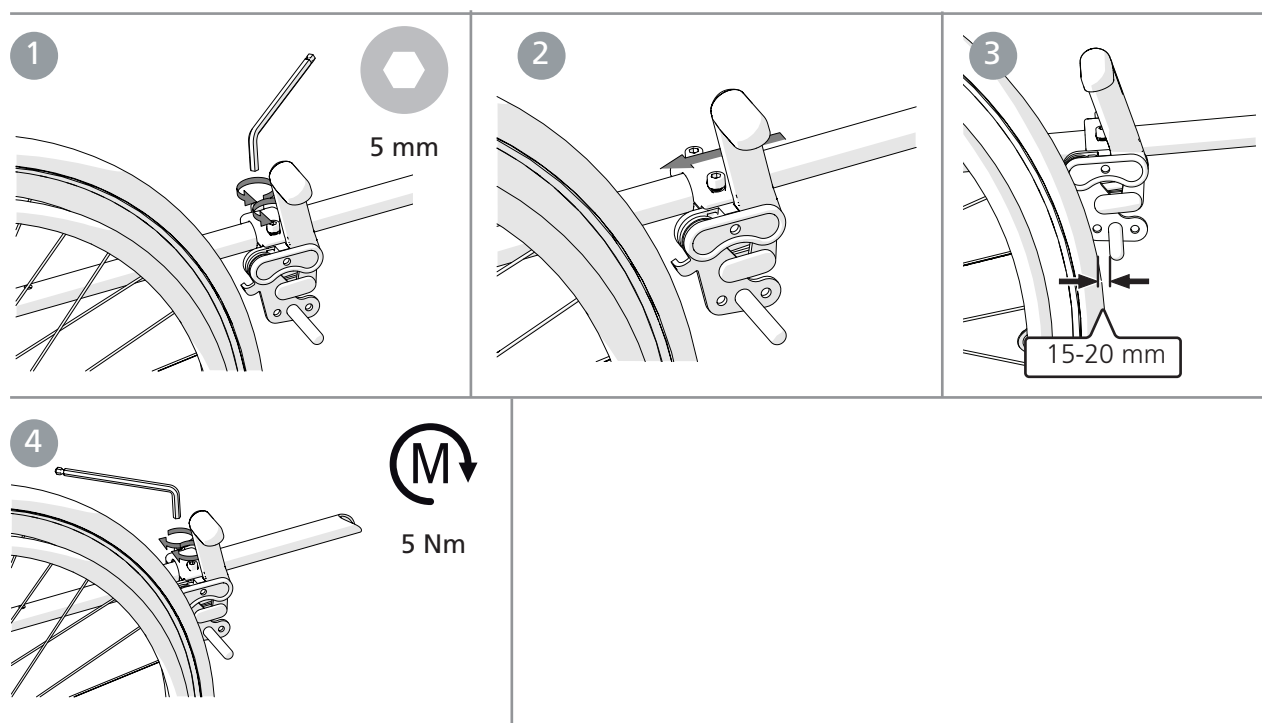


If the rear wheel is positioned in one of the forward settings it is possible for a leg support, when swung to the side, to unlock an applied brake. To avoid this, the brake handle should be adjusted so that you attain a "neutral" position. In this way the brake will only unlock when the handle is pushed close to the rear wheel.

1. Loosen the screw on the inside of the brake handle.
2. Remove the oval plate and turn it 180°. (The pin on the inside of the plate must sit in the upper hole.)
3. Screw it into place again.

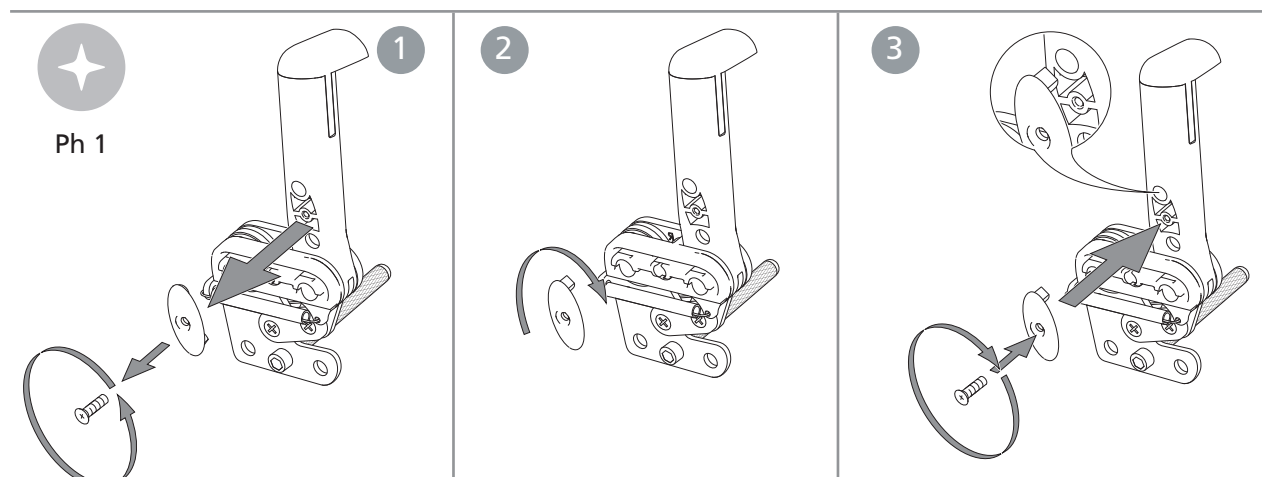

10. Readjustments

Adjustment of brake pin in certain rear wheel positions

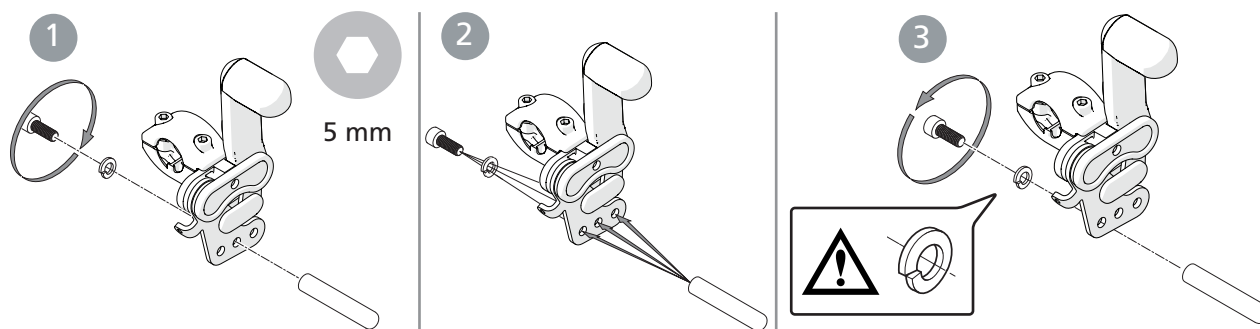






Leg supports, variable angle adjustment



Never stand on the footplates as there is a risk of tipping! When lifting a wheelchair, remove the leg supports and lift using the frame (lower or upper part).



Always make sure that the screw and footplate knob are installed even if the footplate is not used.

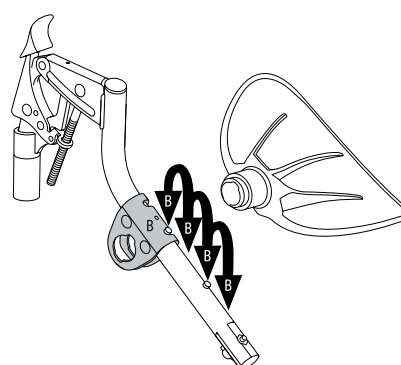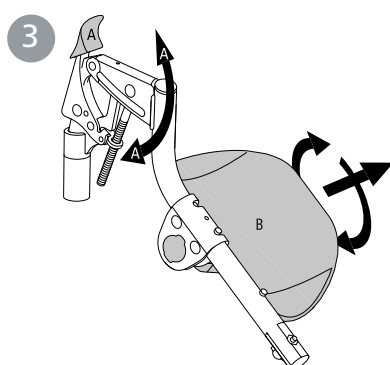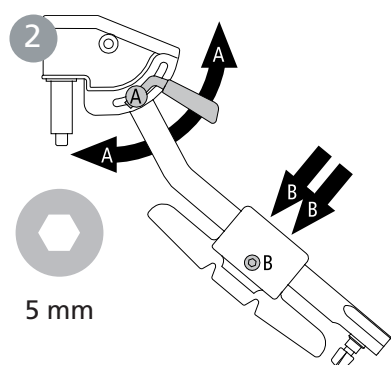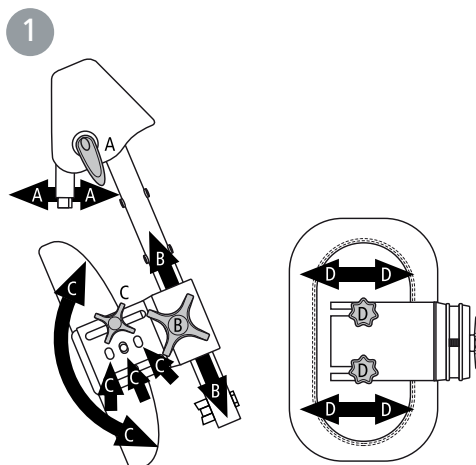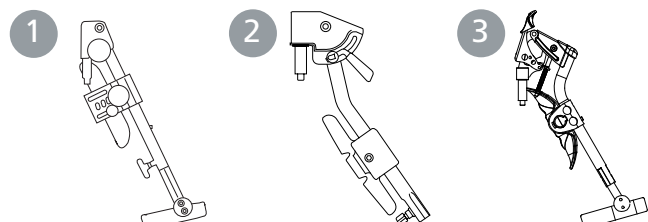

Swing-Away

Active



10. Readjustments

Amputee leg support




Never sit on the amputee leg support. When lifting the wheelchair, remove the amputee leg support, swing away the other leg support (or remove it) and lift using the frame (lower or upper part).



4 mm



6 mm

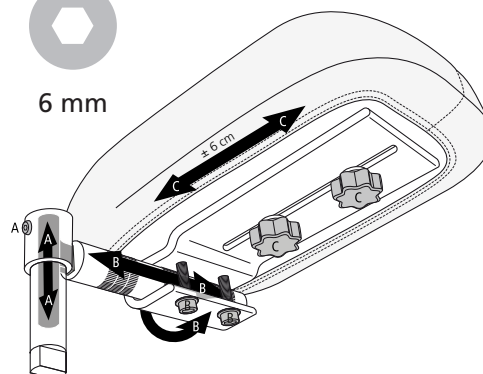

Anti-tip

The anti-tips are height-adjustable.



After adjusting the seat height, centre of balance or back support angle, always ensure that you check the function of the anti-tips.

The anti-tips are mounted on the tilters.



Anti-tip active is mounted on the inside of the rear wheel bracket.



5 mm



5 Nm

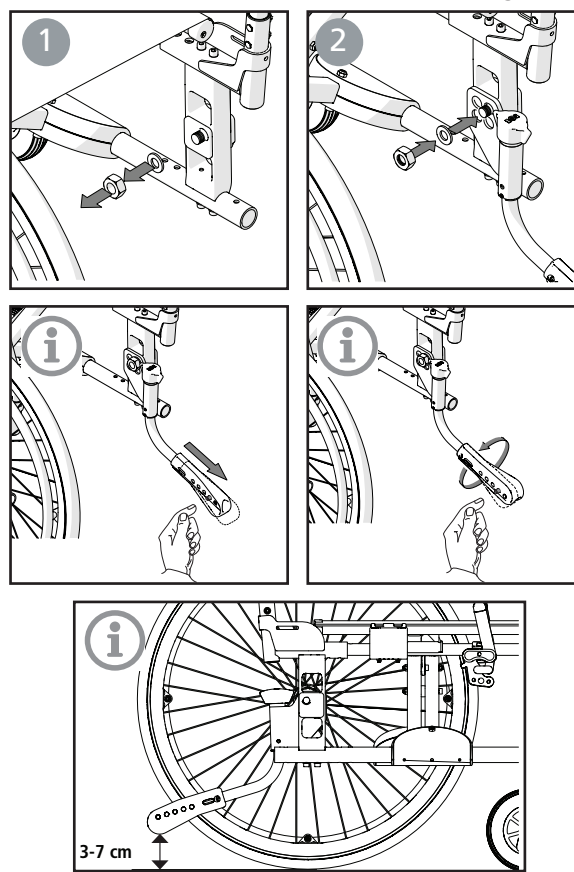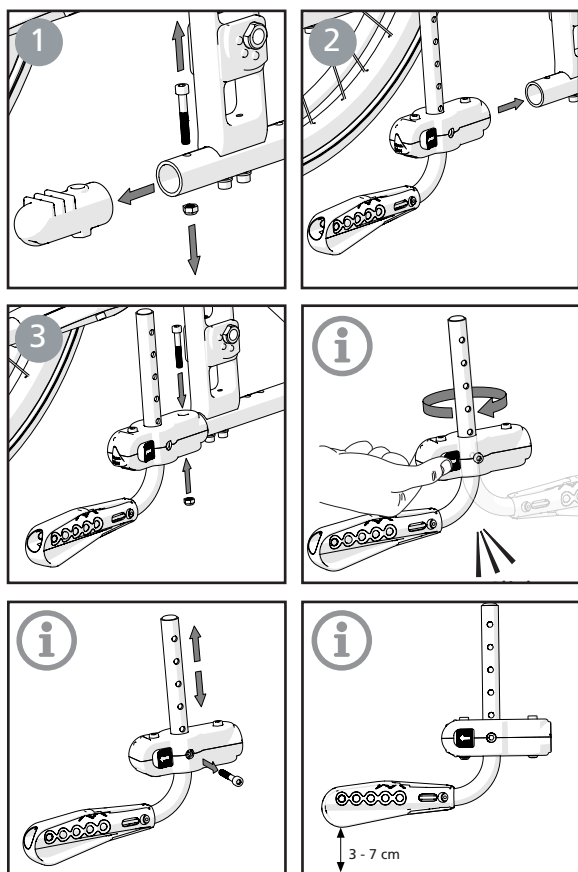

Swing-Away

Active

●

●

●

10. Readjustments


Swing-Away
Active

Changing camber angle


You can change the camber angle by changing the washers between the hub sleeve and the nut. (See also section "Seat height adjustment using the rear wheel position/rear wheel").

 Check that the hub axle of the rear wheel is set correctly. The button in the hub should pop out to its original position.


1

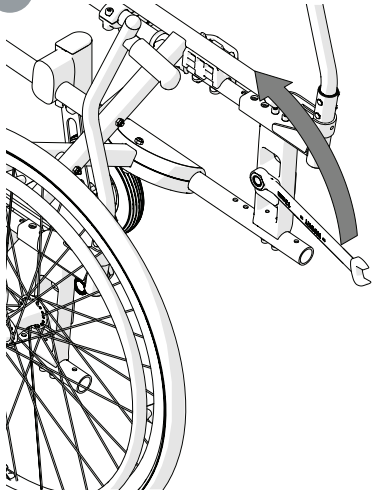



2

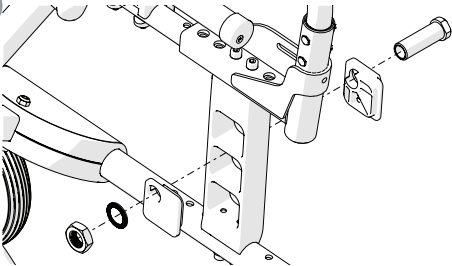





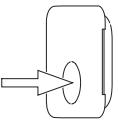

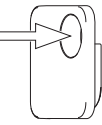


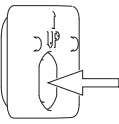

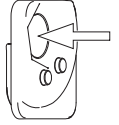


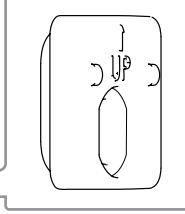

3

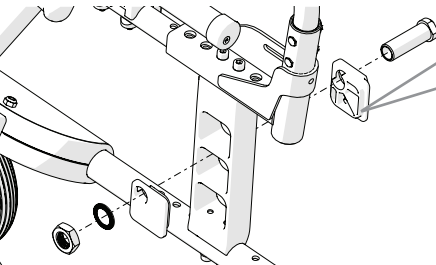

4

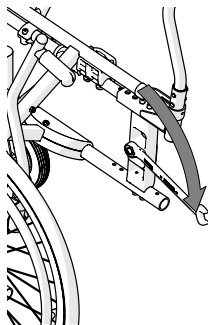

Toe-in, toe-out

Due to the design there is no need for a toe-in correction. Even when the rear wheel position is adjusted in height or centre of balance the rear wheel attachment makes sure that the adjustment of the toe-in is always correct.

10. Readjustments

Handrims

The way in which the user is able to grip the handrims is influenced by the handrim's material and its distance from the wheel. Foam rubber and plastic-coated handrims are available as options.



Foam rubber/plastic-coated hand rims give more friction.

Violent braking can cause mild friction burns (e.g. blisters).




Be aware that when passing through narrow spaces there is a risk of getting your fingers caught.



There is a risk of the user's fingers getting caught in the spokes. We recommend the use of spoke guards.

Adjusting the distance of the handrim

Aluminium, stainless steel, plastic-coated and foam rubber handrims for 24" and 22" rear wheels: The distance between the wheel and the handrim can be adjusted by moving or exchanging the two spacers.

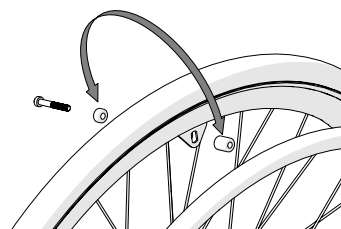



4 mm

Drive wheels with quick-release axle

The distance between the hub sleeve and the locking balls can be adjusted by adjusting the nut on the outside of the hub. The distance between the hub sleeve and the locking balls should be 0.2 to 1.2 mm.



When assembling, always make sure that the quick-release button in the hub springs out properly. Test by pulling the wheel.





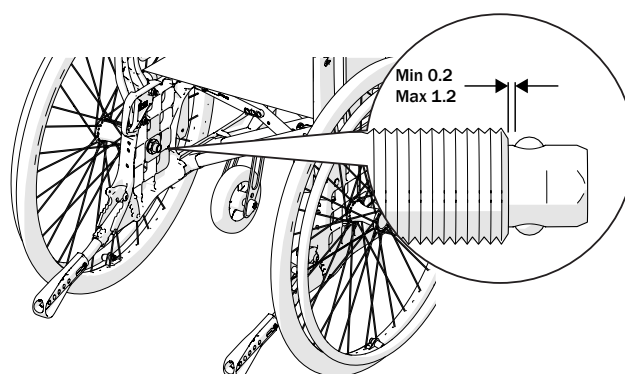
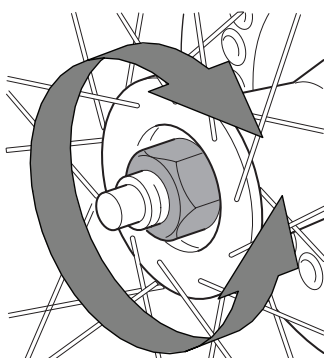


10. Readjustments

Mounting/dismantling upholstery

The seat upholstery is fastened lengthways to the seat frame, and can easily be removed from the frame by unscrewing the end caps.

The back support upholstery is removed by loosening the Velcro straps.

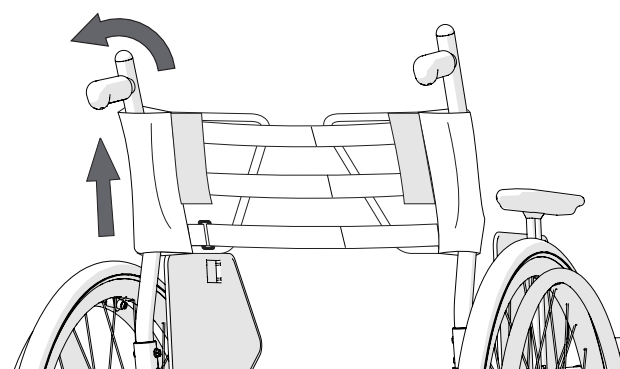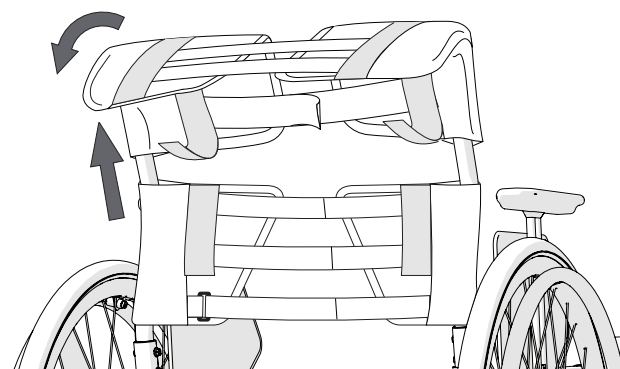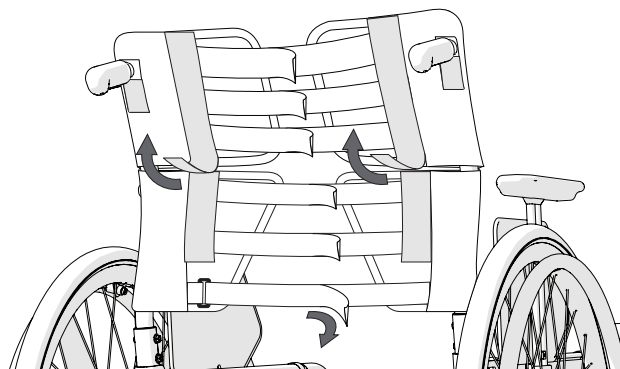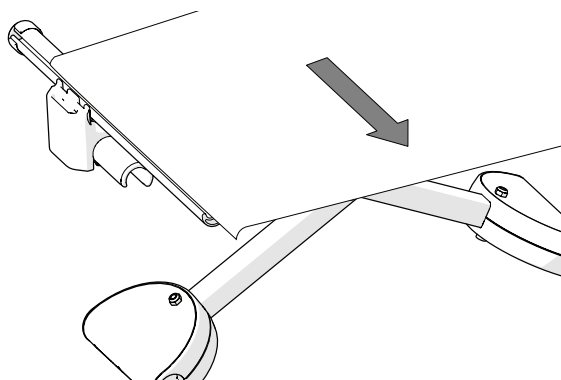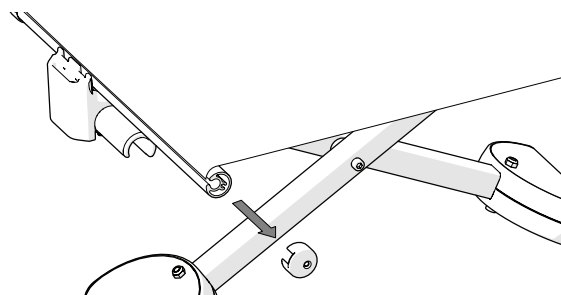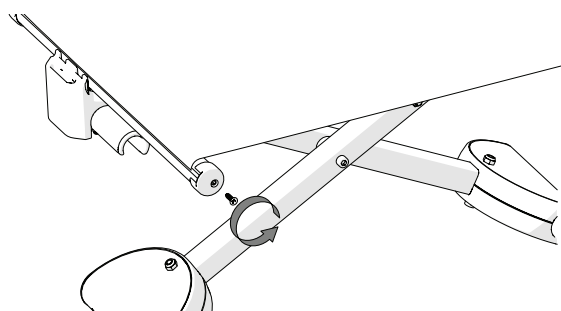

Swing-Away

Active



11. Accessories - adjustment, handling

Assembly instructions are always provided with accessories when they are delivered from Etac. Instructions are also available on our website www.etac.com


Seat cushion and comfort wedge

The cushion is cut to the required depth at the front or back edge. When measuring the length, ensure that the cushion is properly positioned between the back tubes with the rounded corners to the rear.




The cushion is a standard model and is not suitable for users with sitting sores.

A comfort wedge is a useful accessory for the seat cushion. It evens out the front edge of the seat, giving even greater seat comfort. Place the comfort wedge in the cushion-cover's inside pocket, under the cushion. The arrow on the comfort wedge should be on the top-side pointing forwards.

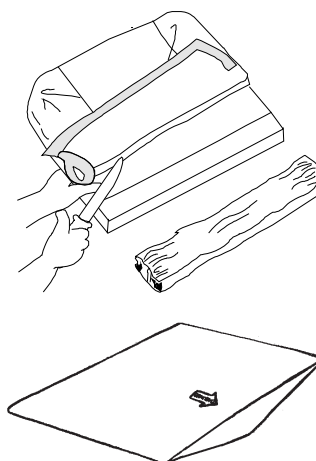

Sideguard/Splash guard

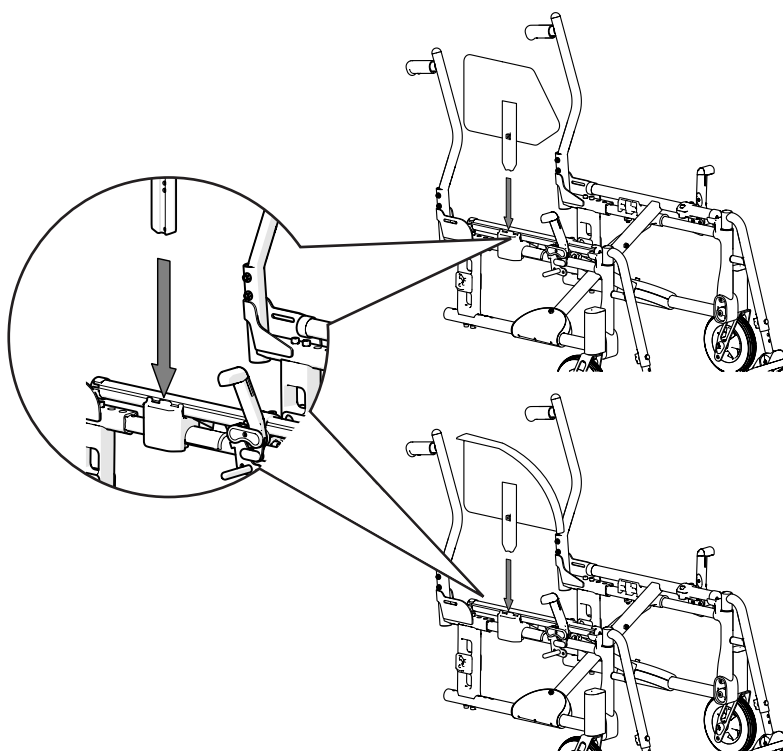

Swing-Away

Active



11. Accessories - adjustment, handling

Head support



The head support is mounted onto the push brace (option, Back support 2C). The head support is adjustable in height, depth and angle, and is detachable.



Before the head support is adjusted, ensure that the user has a good, secure sitting posture.



Check the balance of the wheelchair when the user leans against the head support. We recommend the use of anti-tips.

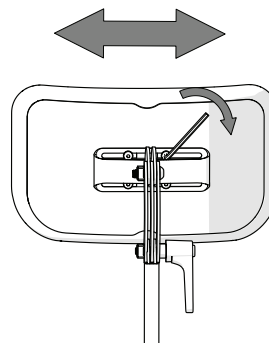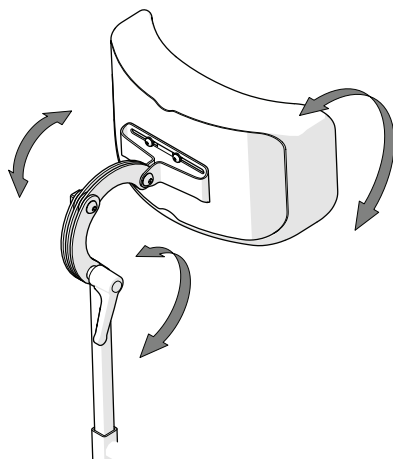

Swing-Away	Active
●	●



12. Technical data

	Type of chair	Art. no. Crissy	Total width cm	Seat depth from back support, cm	Seat width	Back support angle	Weight	Max user weight
Crissy Swing-Away	37.5 cm	13230103	57.5	37-45	37.5 cm	5°	13.5 (29.8 lbs)	125 kg
	40 cm	13230105	60.0	37-45	40 cm	5°	13.6 (29.9 lbs)	125 kg
	42.5 cm	13230107	62.5	37-45	42.5 cm	5°	13.7 (30.2 lbs)	125 kg
	45 cm	13230109	65.0	37-45	45 cm	5°	13.8 (30.4 lbs)	125 kg
	47.5 cm	13230111	67.5	37-45	47.5 cm	5°	13.9 (30.6 lbs)	125 kg
	50 cm	13230113	70.0	37-45	50 cm	5°	14.0 (30.8 lbs)	125 kg
Crissy Active	37.5 cm	13230203	57.5	37-45	37.5 cm	5°	13.0 (28.6 lbs)	125 kg
	40 cm	13230205	60.0	37-45	40 cm	5°	13.1 (28.9 lbs)	125 kg
	42.5 cm	13230207	62.5	37-45	42.5 cm	5°	13.2 (29.1 lbs)	125 kg
	45 cm	13230209	65.0	37-45	45 cm	5°	13.3 (29.3 lbs)	125 kg
	47.5 cm	13230211	67.5	37-45	47.5 cm	5°	13.4 (29.5 lbs)	125 kg
	50 cm	13230213	70.0	37-45	50 cm	5°	13.5 (29.8 lbs)	125 kg

	Seat height rear cm	Seat height front cm	Back sup- port height cm	Transport width cm	Total height cm	Total length cm	Leg sup- port length cm	Folded length cm	Leg support angle	Handrim size cm
Crissy Swing-away	38.5-51	40.5-55	35-47.5	32	73-98.5	94-102	30-64	71-79	108°-114°	56-60
Crissy Active	38.5-51	40.5-55	35-47.5	32	73-98.5	90-98	37-57	90-98	110°	56-60

Crissy Swing-Away	The weights given are for a chair with 24" rear wheels with quick-release hubs, back support 2B, high-pressure tyres, aluminium handrims, leg supports, foot plates, brakes, mudguards, 6" castor wheels and medium fork. The measurements given may vary by $\pm 2\%$.
Crissy Active	The weights given are for a chair with 24" rear wheels with quick-release hubs, back support 2B, high-pressure tyres, aluminium handrims, foot plates, brakes, mudguards, 5" castor wheels and medium fork. The measurements given may vary by $\pm 2\%$.

		Tyres							
Max. incline sideways	Max. incline down		High pressure 22"	High pressure 24"	Solid 22"	Solid 24"	Solid inner tube 22"	Solid inner tube 24"	High pressure 24" (Carbon)
		PSI/Bar	110/7.5	110/7.5	-	-	-	-	110/7.6
			550 mm	590 mm	550 mm	600 mm	550 mm	590 mm	592
>20°	>20°		25 mm	25 mm	32 mm	32 mm	25 mm	25 mm	23 mm

Fe/ALU/PA	ALU/PA	PES	Fe/CF/ALU	PP/TPR	ALU/PA/PE/PU

13. Service and maintenance

Cleaning

It is important to keep the wheelchair clean, both for your own comfort and the longevity of the chair. It is equipped with a drainage hole which ensures that it is easy to wash and dry.

Back support and seat upholstery, padded armsupport covers: See the washing instructions on the product.

Chassis, push handles and arm supports: Everyday maintenance: Clean the chassis with a commercially available detergent without abrasives having a pH between 5 - 9. Rinse and dry.

Disinfection

Disinfect the chair with 70% disinfectant spray. Air dry.

Advanced maintenance/Reconditioning: See reconditioning instruction at www.etac.com.

Please observe the application and safety instructions of the disinfectant manufacturer.

Rear wheel/castor

Tyre/inner tube: Check the tyre pressure (see side of tyre) at least once a month, also check the tread.

Spokes: Loose spokes can lead to wheel wobble. Consult a cycle dealer or your Technical Aids Centre if it is necessary to adjust the spokes.

Wheel axles: Clean the wheel axles from hair and dirt as necessary.

Ball bearings: Require no maintenance.

Handrims: If a handrim is damaged in such a way that it can cause hand injuries, it should be replaced.

Front fork attachment:
To achieve the best operating conditions, the attachments should be installed at 90°. Also check that the front forks are correctly tightened. See section "*Disassembly/assembly of the front fork*".

Brakes

The braking effect is dependent on the air pressure in the tyres. Heavy dirt can have a negative effect on the

brake mechanism. If adjustment is necessary, see section "*Brake, adjustment*".

Touch-up paint

Touch-up paint is available for minor scratches and chips in all the frame colours available for the wheelchairs.

Miscellaneous

If there is a fault in your wheelchair you should contact your dealer or Technical Aids Centre. A defective wheelchair must not be used.

If your chair needs reconditioning or repair, only original parts from Etac or components with equal quality should be used.

Etac will not be held responsible for damage or injury caused by use of non-original parts.



When necessary lubricate moving parts/joints with bicycle oil or similar.



14. Tests and guarantees

CE marking, tests:

The product conforms to the requirements of the Medical Devices Regulation (EU) 2017/745.

The product has been tested and complies with the requirements set out in EN 12183, which includes ISO 8191-2.

Combination agreements exist – see www.etac.com.

Etac's wheelchairs have been tested in accordance with ISO 7176-19. They have been crash-tested by the Technical Research Institute of Sweden. Testing was carried out using restraining device UNWIN_WWR/ATF/K/R and 3-point seat belt 907523, Klippan Safety AB.

Service life:

The anticipated durability and service life of the main product is at least five years when used as intended, and in accordance with the safety instructions, reconditioning manual and maintenance instructions available in the user manual. The main product comprises the seat and back frame. Other parts/accessories are managed in accordance with the manual and reconditioning manual. The actual service life may vary depending on how often and how intensively the product is used, but is at most 10 years. If the product is intended for use after the service life period specified by Etac, it is the responsibility of the product owner to ensure product functionality in accordance with the manual and reconditioning instructions. If this cannot be ensured then the product should be taken out of use.

The following methods of surface treatment have been used:

Lacquered surfaces = Polyester powder coating or ED lacquering

Non-lacquered aluminium parts = Anodised coating

Non-lacquered steel surfaces = Galvanised

Non-lacquered/galvanised steel surfaces = Stainless steel

The product is protected from corrosion.

Guarantee:

5-year guarantee against material and manufacturing defects. For terms and conditions, see www.etac.com.

Special adaptations:

comprise everything that falls outside the instructions and settings in the manual. Wheelchairs specially adapted by customers are not eligible for Etac's CE marking. Etac's guarantee no longer applies. If there are any doubts concerning the adaptation, please contact Etac for advice.

If the wheelchair is used in combination with another product, neither product may keep the CE marking unless a combination agreement is in place. Please contact Etac for up-to-date information.



To be sorted according to national regulations



Etac Sverige AB
Box 203
334 24 Anderstorp
Sweden

Tel 0371-58 73 00
Fax 0371-58 73 90
info@etac.se
www.etac.se

Etac R82 GmbH
Industriestraße 13
45699 Herten
Germany

Tel 02366-5006-0
Fax 02366-5006-200
info@etac.de
www.etac.de

Etac AB (export)
Box 203
334 24 Anderstorp
Sweden

Tel 46 371-58 73 30
Fax 46 371-58 73 90
info@etac.se
www.etac.com

Etac R82 B.V.
Tinweg 8
8445 PD Heerenveen,
Nederland

Tel +31 72 547 04 39
info.holland@etac.com
www.etac.com

Etac AS
Pb 249,
1501 Moss,
Norway

Tel 815 69 469
Fax 69 27 09 11
hovedkontor.norge@etac.com
www.etac.no

Etac Ltd.
Unit D4A, Coombswood Business
Park East
Coombswood Way, Halesowen
West Midlands B62 8BH
United Kingdom

Tel 0121 561 2222
Fax 0121 559 5437
enquiries@etac.uk.r82.com
www.etac.com

Etac A/S
Egeskovvej 12
8700 Horsens
Denmark

Tel 79 68 58 33
Fax 75 68 58 40
info@etac.dk
www.etac.dk

R82, Inc.
12801 E. Independence Boulevard
P.O. Box 1739
Matthews, NC 28106, USA

Tel 800 336 7684
Fax 704 882 0751
Information@snugseat.com
www.etac.com



Etac Supply Center AB
Långgatan 12
SE-334 33 Anderstorp



 **etac**®

