


Important: Remove this guide before handing the device to the patient. Pressures should only be adjusted by qualified personnel.

S8 AutoSet Spirit™ II System

Clinical Quick Setup Guide

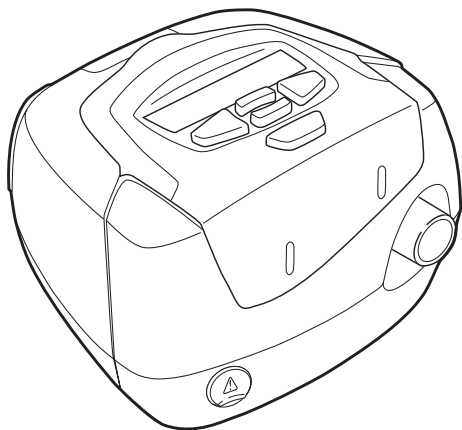

Setup is quick and easy . . . *Just Go!*

The following settings have been pre-set to accommodate typical patient treatment values.

CPAP mode

Start CPAP:	6 cm H ₂ O
CPAP:	10 cm H ₂ O
Ramp:	20 minutes
Max Ramp:	30 minutes
Humid:	None
SmartStart:	Off
Reminders:	Off
Mask:	Ultra
Leak Alert:	Off
Tube Length:	2 m
Language:	English
Smart Data:	On
Auto Appear:	On
Patient Menu:	Std
EPR:	Off
EPR Level:	0
EPR Inhale:	Med
EPR Access:	Off

AutoSet mode

Min Press:	4 cm H ₂ O
Max Press:	20 cm H ₂ O
Settling:	20 minutes
Max Settling:	30 minutes
Humid:	None
SmartStart:	Off
Reminders:	Off
Mask:	Ultra
Leak Alert:	Off
Tube Length:	2 m
Language:	English
Smart Data:	On
Auto Appear:	On
Patient Menu:	Std
EPR:	Off
EPR Level:	0
EPR Access:	Off

1

Accessing the clinical menu

- Turn on the device and wait until the standby (RAMP/SETTLING) screen appears on the LCD. At this stage you can adjust the ramp/settling time.
- Press  and  simultaneously for three seconds to access the clinical menu.

2

Accessing the quickview menu

- Press  and  simultaneously for three seconds to access the quickview menu.

3


Changing options and settings

- Press  to access the SETTINGS menu screen.
- Use  and  to scroll through the RESULTS, OPTIONS and SERVICING menu screens.
- From a menu screen, use  to access the parameters in that menu.
- Use  and  to scroll through menu parameters.
- Use  to view or change each parameter setting as required.
- Press  until you return to the standby (RAMP/SETTLING) screen.

See the Clinical Guide for full details.




S8 AutoSet Spirit II clinical menus

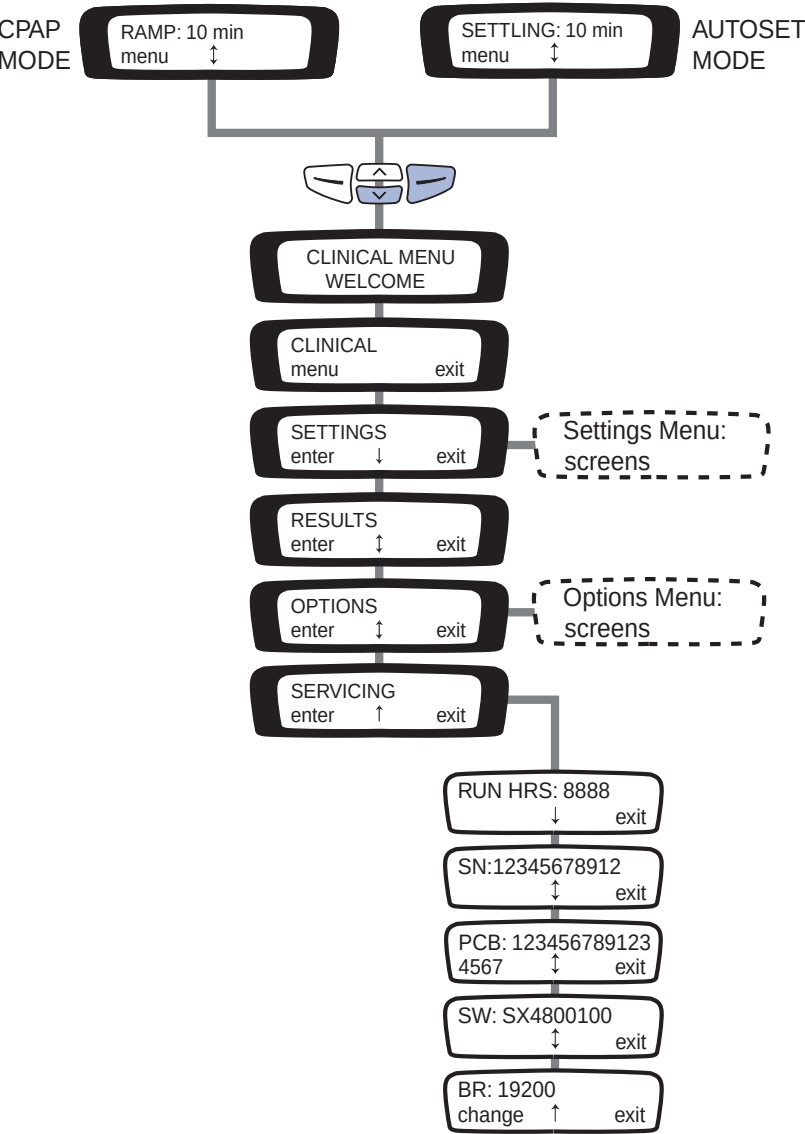


S8 AutoSet Spirit II clinical menus

Settings Menu: screens

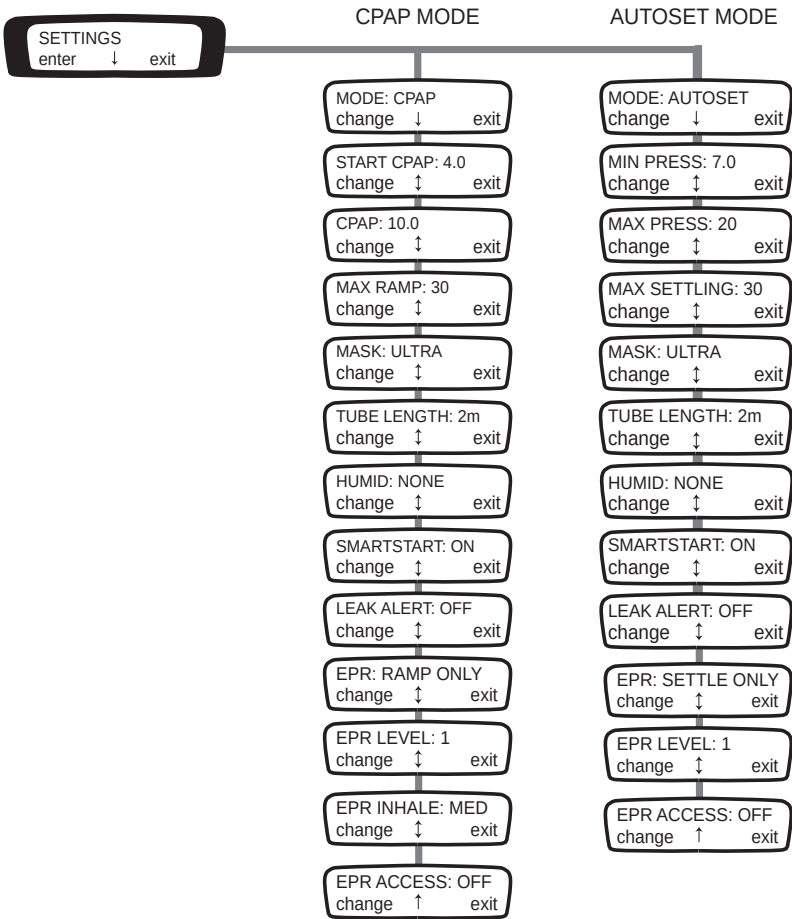

S8 AutoSet Spirit II clinical menus


Options Menu: screens

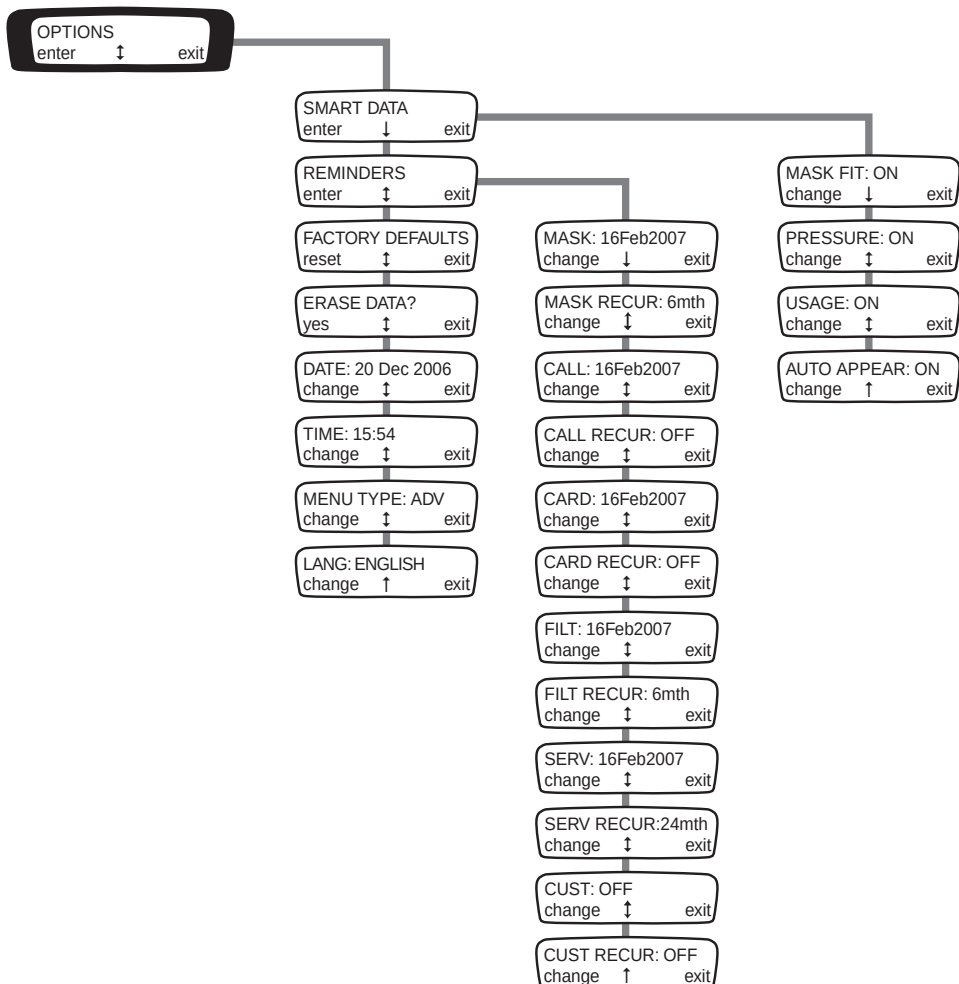

Note: For each reminder, you can set a time period for that reminder to keep re-occurring. For these recurring reminders to work, you must set the date for the first occurrence of the reminder as well as a recurrence time period.


For example, if you set MASK: 30Jun2008 and MASK RECUR: 6mth, the first reminder to replace the mask is displayed on 30 June 2008. The recurring reminder occurs 6 months from the date the reminder is acknowledged by the patient. The reminder will continue to recur at 6-month intervals unless it is set to OFF or the recurrence period is changed.

S8 AutoSet Spirit II quickview menu

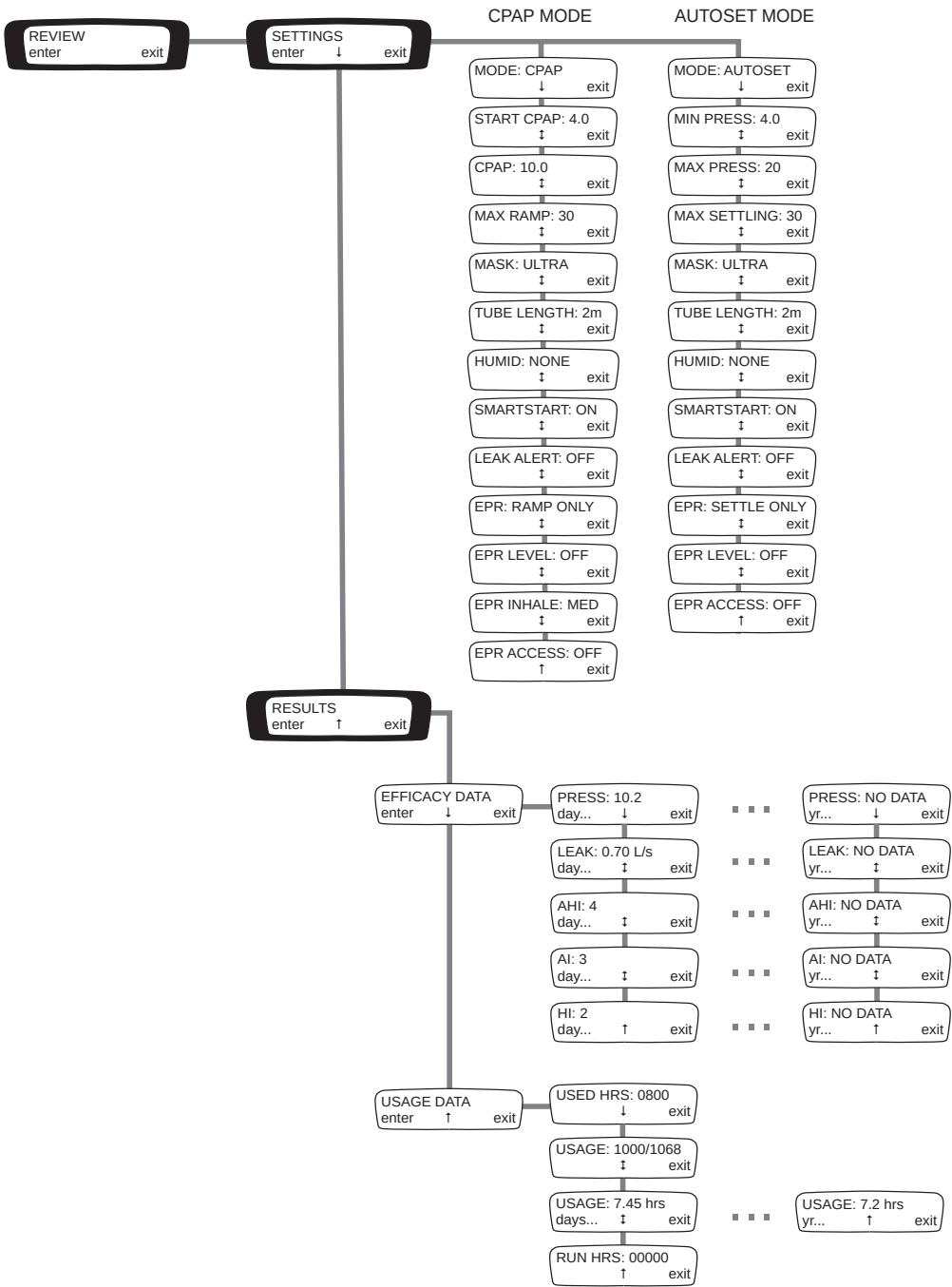

S8 AutoSet Spirit II patient menus


Detailed menus

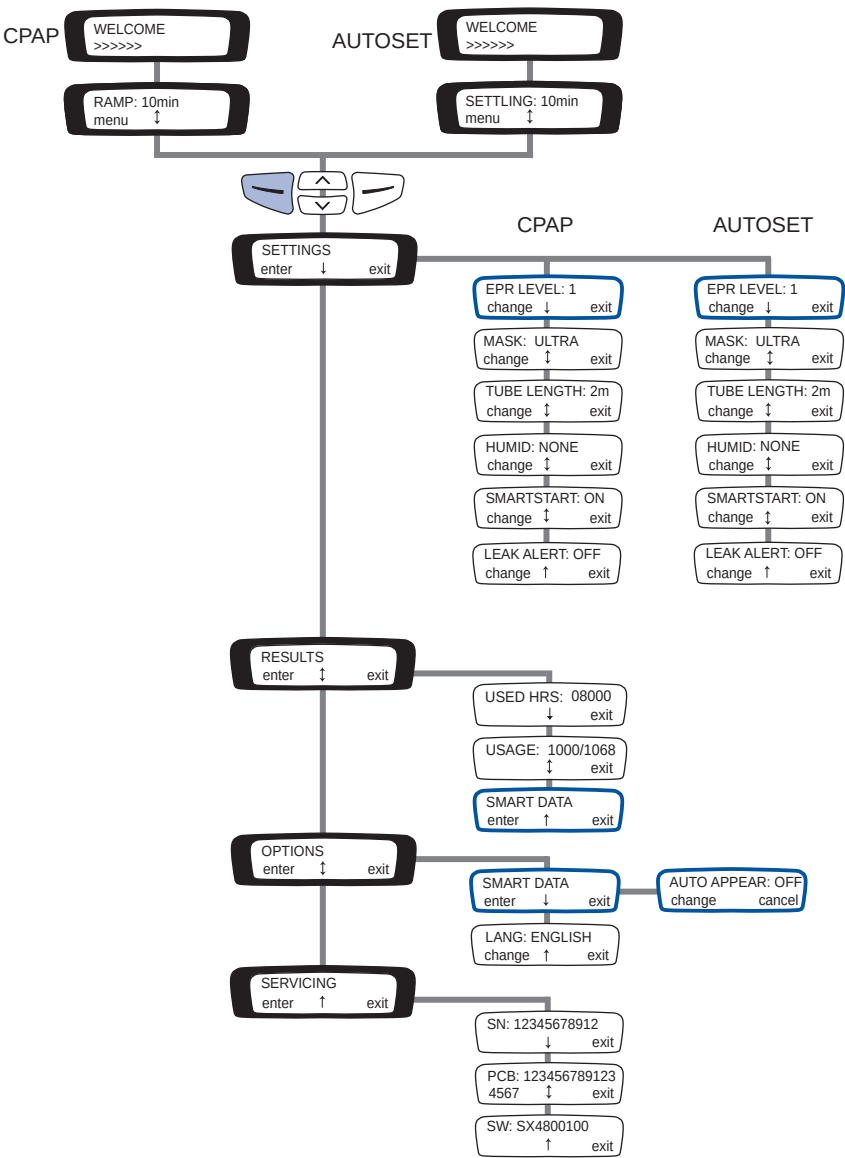

This screen only appears if access is set by clinician.

S8 AutoSet Spirit II patient menus

Standard menus

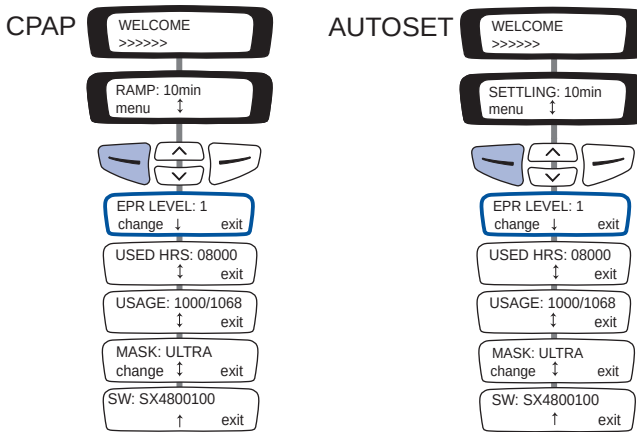

 This screen only appears if access is set by clinician.

ResMed Pty Ltd (Manufacturer) 1 Elizabeth Macarthur Drive Bella Vista NSW 2153 Australia • ResMed Corp (US Designated Agent) 14040 Danielson Street Poway CA 92064-6857 USA • ResMed (UK) Ltd (EU Authorised Representative) 96 Milton Park Abingdon Oxfordshire OX14 4RY UK • **ResMed Offices** Australia, Austria, Belgium, Brazil, China, Finland, France, Germany, Greece, Hong Kong, India, Ireland, Italy, Japan, Malaysia, Mexico, Netherlands, New Zealand, Norway, Portugal, Singapore, South Africa, Spain, Sweden, Switzerland, UK, USA (see www.resmed.com for contact details).

Protected by patents: AU 691200, AU 697652, AU 709279, AU 724589, AU 730844, AU 736723, AU 750095, AU 750761, AU 764761, AU 779327, AU 2002231456, CN 02804930.6, EP 0651971, EP 0661071, EP 0920845, EP 0934723, EP 1028769, EP 1126893, EP 1502618, EP 1144036, JP 3580776, JP 3778797, NZ 504595, US 5199424, US 5245995, US 5522382, US 5704345, US 6138675, US 6240921, US 6363933, US 6367474, US 6425395, US 6502572, US 6591834, US 6705315, US 6745768, US 6817361, US 6988498, US 7040317, US 7100608. Other patents pending.

Protected by design registrations: AU 302180, AU 302181, AU 302182, AU 302183, AU 302184, AU 302185, AU 302186, CN 200430121212.2, CN 200430121213.7, CN 200430121211.8, CN 200430121214.1, CN 200430121215.6, CN 200430121216.0, CN 200430121217.5, EU 269436, EU 815071, HK 0412901.9, JP 1257662, JP 1257663, JP 1266956, JP 1266957, JP 1267270, JP 1267271, NZ 405614, NZ 406923, NZ 406924, NZ 406925, NZ 406926, NZ 406927, NZ 406928, US D544598, US D553228, US D557406, US D557407, US D560795, US D561891, US D565160. Others pending.

S8 AutoSet Spirit and SmartStart are trademarks of ResMed Pty Ltd and AutoSet is registered in U.S. Patent and Trademark Office.