

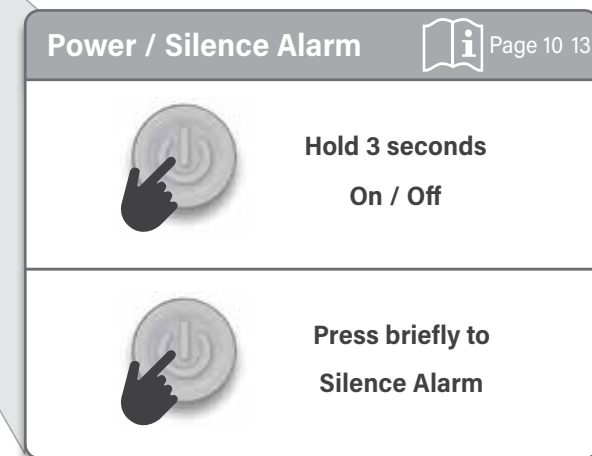
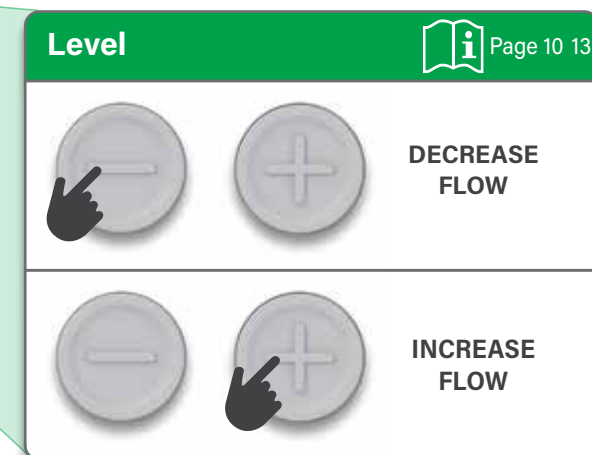
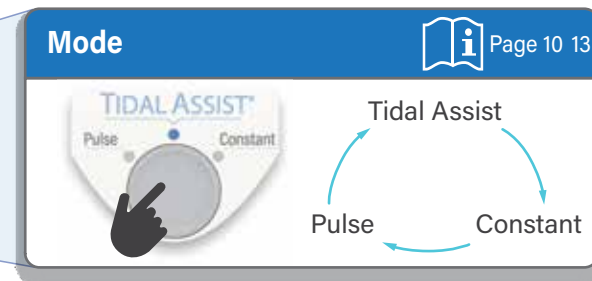


Quick Setup Guide *

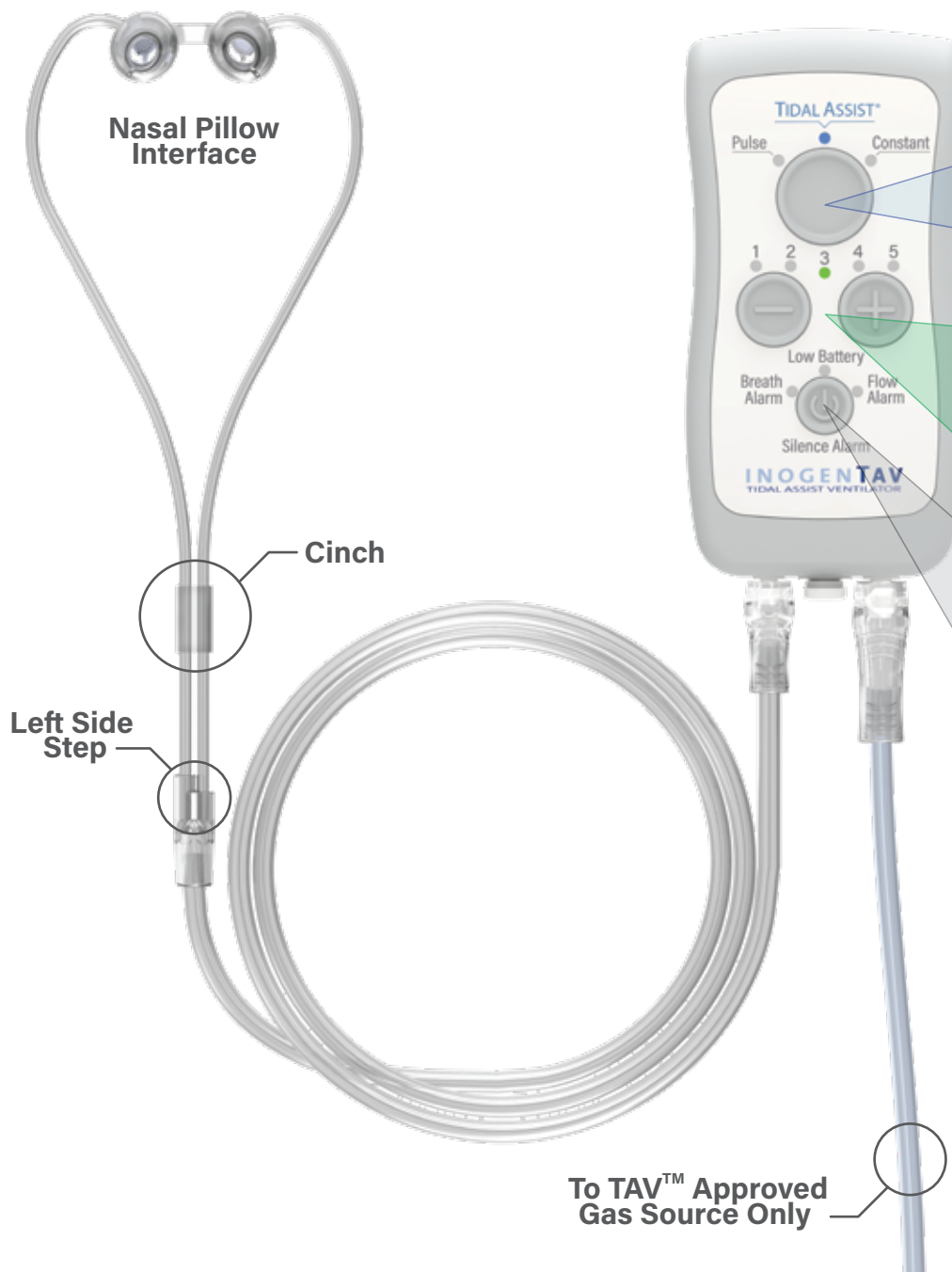

* Reference Instructions For Use for additional information
96-09783-00-01 C, TAV QUICK REFERENCE GUIDE

Quick Troubleshooting Guide *

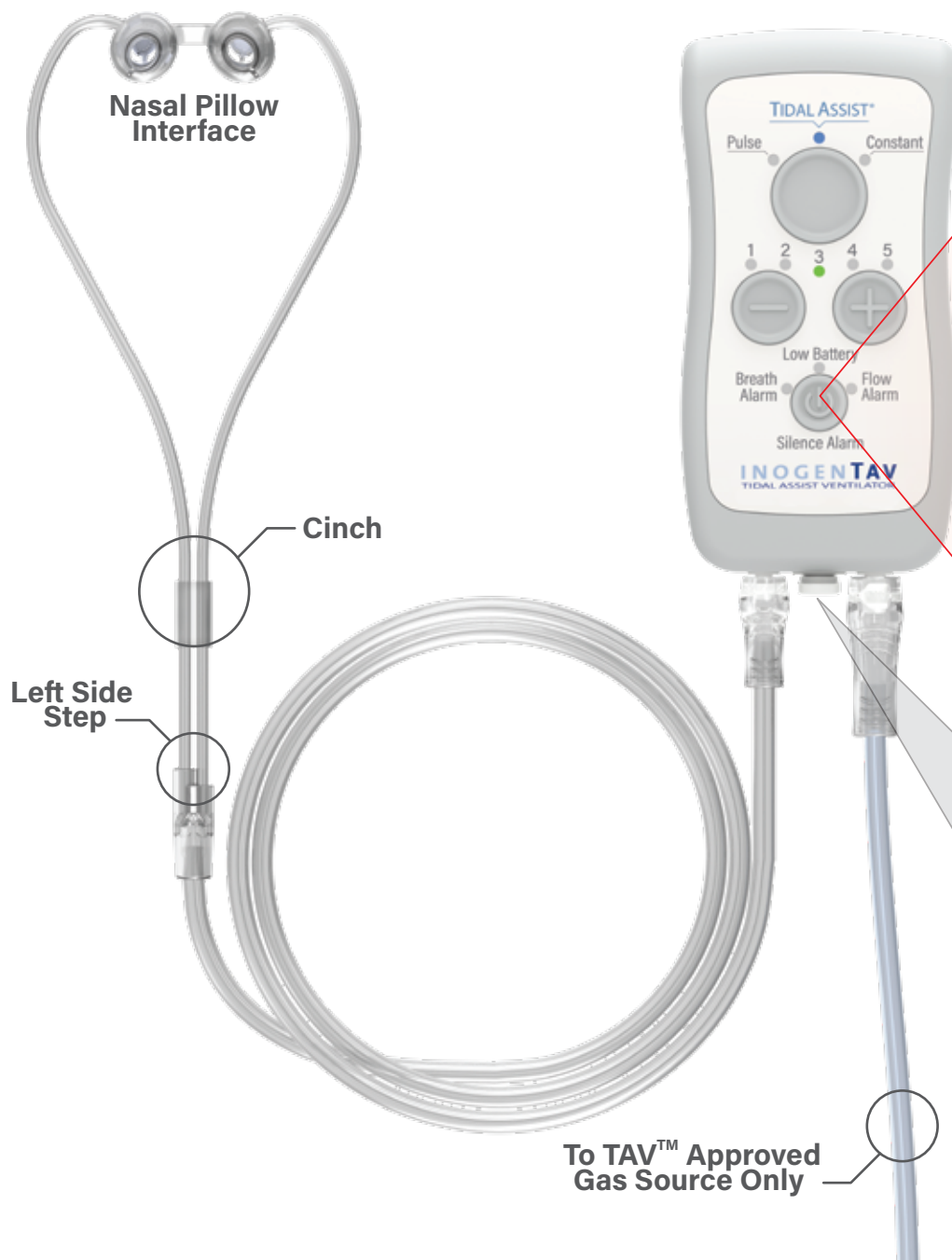

Alarms

Page 15



1. Check patient
2. Check nasal pillow interface fit
3. Check tubing for kinks and leaks
4. Check all connections
5. Turn Device off and then back on



Battery low - replace battery



1. Check tubing for kinks and leaks
2. Check connections
3. Check gas source
4. Check regulator adaptor
5. Turn Device off and then back on



1. Check Bypass Valve is closed
2. Turn Device off and then back on

Note: if alarms cannot be cleared, turn off device and call customer support. Use Bypass Valve if needed.

Bypass Valve (2LPM constant)

Page 14

Power off



Hold 3 seconds



Pull

* Reference Instructions For Use for additional information

96-09783-00-01 C, TAV QUICK REFERENCE GUIDE