

LPC01 Track Joining Kit - Fitting Instructions

This leaflet details how to use the LPC01 track joining kit to join a custom rail arrangement of the Thorlabs Curtain system. For the continuing safety of the operators of the equipment and for the protection of the equipment itself, the operator should take note of the **Warnings, Cautions and Notes** throughout this leaflet, and, where applicable, in any associated literature.

1.1 Assembly

1.1.1 Parts Included

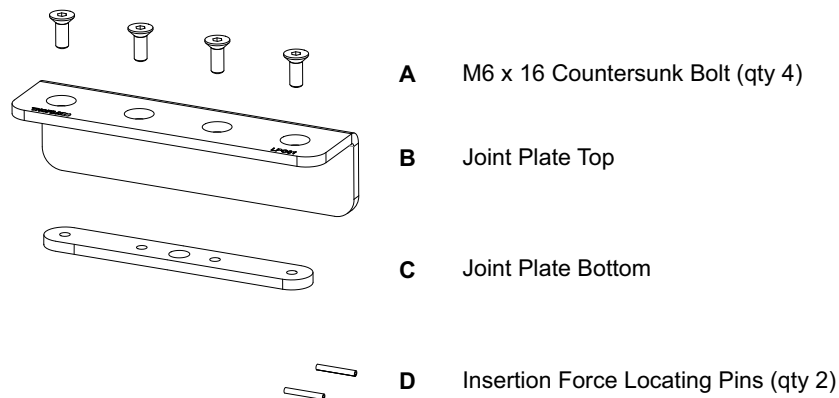


Fig. 1.1 Joining Kit Parts List

1.1.2 Joining the Curtain Rail

- 1) Lay the parts on the floor, or other suitable large, flat surface roughly in the positions shown in Fig. 1.2.
- 2) Insert the two spring pins (D) into the joint as shown below. Push in each pin up to the indent mark, such that about half its length is protruding.
- 3) Using the M6 x 16 countersunk bolts (A) and joint plates (B and C) provided, assemble the curtain rail as shown below. .

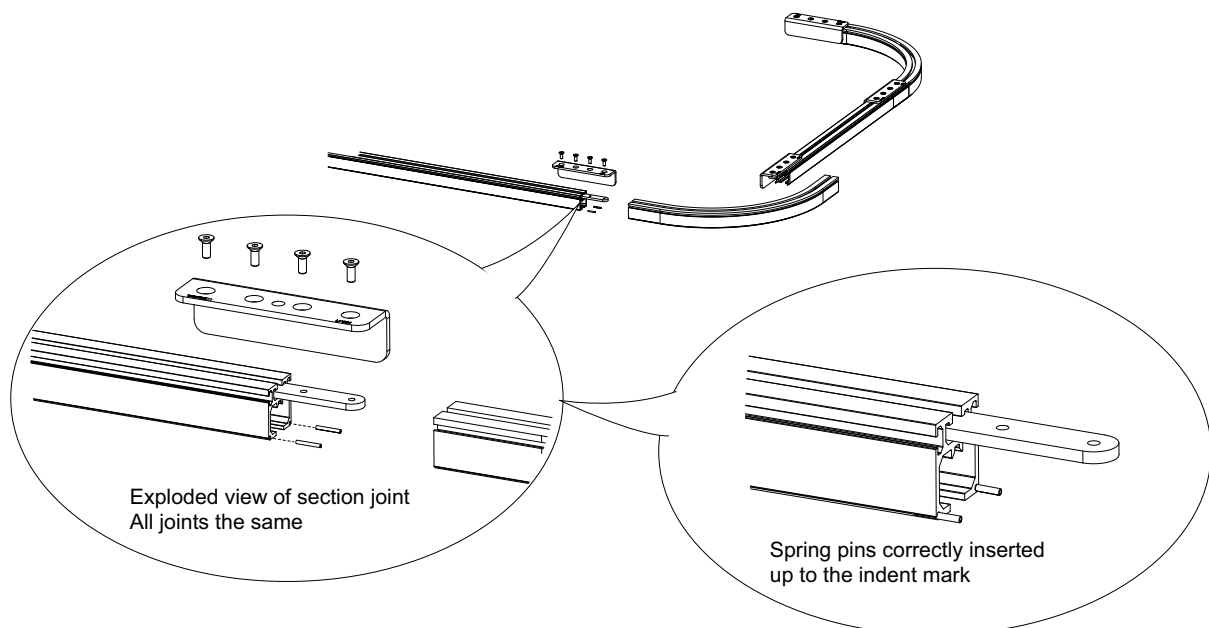


Fig. 1.2 Joining the Curtain Rail



Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



USA, Canada, and South America

Thorlabs, Inc.
sales@thorlabs.com
techsupport@thorlabs.com

Europe

Thorlabs GmbH
europe@thorlabs.com

France

Thorlabs SAS
sales.fr@thorlabs.com

Japan

Thorlabs Japan Inc.
sales@thorlabs.jp

UK and Ireland

Thorlabs Ltd.
sales@uk.thorlabs.com
techsupport.uk@thorlabs.com

Scandinavia


Thorlabs Sweden AB
scandinavia@thorlabs.com

Brazil

Thorlabs Vendas de Fotônicos Ltda.
brasil@thorlabs.com

China

Thorlabs China
chinasales@thorlabs.com

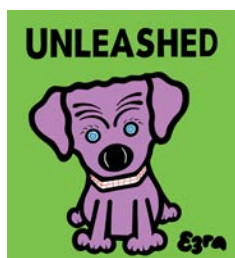

THORLABS

www.thorlabs.com

