

Sock / Stocking Aid, Deluxe Instruction & Safety Manual

Please read this manual before operation, for
maintenance instructions and safe usage



The enabling power
of applied knowledge



Introduction

The Sock / Stocking Aid is an assistive device for users with limited mobility, allowing them to put on socks or stockings without the need to bend down. The Terry cloth cover keeps the socks and stockings securely in place. Users can bend the board to adjust the fit for different socks.

Specification

Gutter Length	Strap Length	Weight	Material
230 mm	850 mm	88 g	PP(inner board) / Nylon Cloth / Terry Cloth / Cotton(straps)

Instructions for Use

1. Bend the board of the Sock / Stocking Aid to fit into the sock
2. Slide the sock to the end. Make sure the sock does not go over the edge of the Sock / Stocking Aid (figure 1)
3. Hold the straps and place the Sock / Stocking Aid in front of the foot
4. Insert the foot into the sock aid (figure 2)
5. Pull the straps until the sock aid is all the way pulled out from the sock

Figure 1

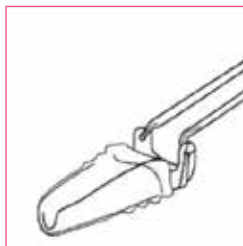


Figure 2



Attention

- Do not put in washing machines
- Periodically inspect for wear and damage.

Maintenance

Hand wash only. Do not put in washing machines.





PO Box 82 Bentley WA 6982
T 08 9351 7500 F 08 9350 5299
info@unicarehealth.com.au
www.unicarehealth.com.au

Authorised Distributor

