

SNOO Smart Sleeper

Hospital Quick Start Guide




What do I need to get started?

- SNOO bassinet (+ optional mobility cart)
- Reusable and/or disposable covers for SNOO mesh and mattress
- Clean SNOO Sleep Sack and fitted mattress sheet
- (optional) iPad with SNOO mobile app


How do I use the SNOO Sleep Sack swaddles?

1. Unzip the sack and spread open the inner bands. Lay the baby down on top of unzipped SNOO sack, placing shoulders 2-3” higher than the top of the sack.
2. Straighten the arms, pull the inner bands around, and attach them very snugly! The bottom edge of the bands should cover most of the hands.
3. Pull diaper flap up and secure over the band.
4. Pull the fabric at the top up and over each shoulder and zip the sack completely.

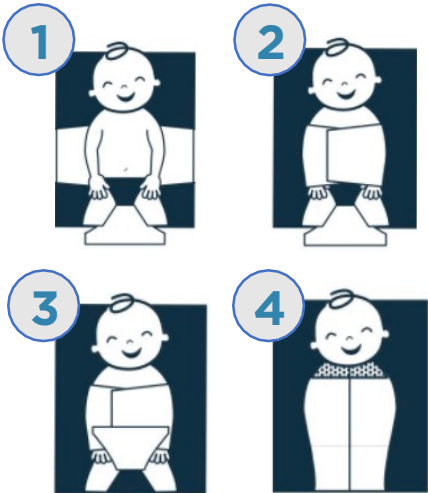

How do I place the baby in SNOO?

- Lay the swaddled baby inside the SNOO on his or her back.
- Slide the loop on the end of each swaddle wing over the safety clips on each side of the platform.

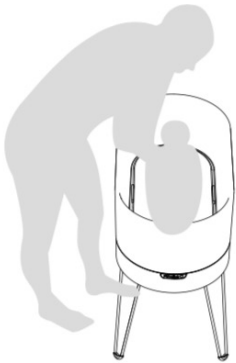

How do I turn SNOO on?

- Press the activity button on the front of SNOO. Once the white light stops blinking (usually takes <1 minute), push the button. The light will turn blue and SNOO will turn on.
- When the baby fusses, SNOO will wait a few seconds before responding. If the crying continues, the bed will move through 4 levels of higher sound/motion until the fussing stops or SNOO times out (after ~3 minutes of continued crying).
- To increase the level manually, push the activity button for several seconds.
- To stop SNOO, simply push the activity button once. To reduce the level of activity, you must use the SNOO Mobile App, or simply wait and SNOO will reduce its level – step by step – over the next 5-10 minutes.

- Baseline
- Level 1
- Level 2
- Level 3
- Level 4
- Action Stopped
- ⚙ Engage Safety Clips
- Ready Mode

Note: Light won't stop blinking? Check that the wing loops are slipped completely onto the safety clips.

Note: If the baby cries for longer than a minute, you can assume she needs a feeding, diaper change, cuddle, or other help. You do not need to wait for SNOO to stop itself.

How do I clean SNOO?

Prior to use:

- Place a waterproof white reusable mesh cover OR a blue disposable mesh cover to completely cover the bassinet.
- (If using blue disposables) Place the sealed mattress inside the removable blue waterproof cover.
- Place a clean cotton mattress sheet over the covered mattress.

Cleaning instructions:

- Use CaviWipes or a hydrogen peroxide-based disinfectant to clean the SNOO, mobility cart, reusable covers, and sealed mattress if visibly soiled, between patients, and every 7 days (if in use). The surface should remain visibly wet for 1 minute at room temperature.
- Launder in accordance with standard hospital laundering practices and protocols; a temperature of at least 160°F (71°C) for a minimum of 25 minutes is commonly recommended for hot-water washing.
- SNOO Sacks MUST BE cleaned between patients and when visibly soiled.

Note: For storage, cover cleaned SNOO with a plastic bag to designate it as clean and ready for use.

How do I turn on Motion Limiter and other optional features?

Manually

- Press the activity button on the front of SNOO. Hold for 30 seconds.

With optional mobile app (requires WIFI access)

- Check to make sure the SNOO is not actively moving while updating the settings.
- Open the SNOO app and head to Main Menu (3 white bars).
- Navigate to Settings, then My SNOO.
- Slide Motion Limiter ON. It will turn blue when activated.
- Head back to Start SNOO. You will now see an “M” at the top right-hand corner.

Optional Features Available with the SNOO Mobile App

- Motion Limiter:** For babies who like less jiggle
- Weaning:** To prepare for the transition to a crib
- Level Lock:** To keep sound and motion fixed
- Preemie:** Gentler rocking for tiny babies

