
Etac Swift Mobil-2


Manual

etac®


78691.13 23-01-18

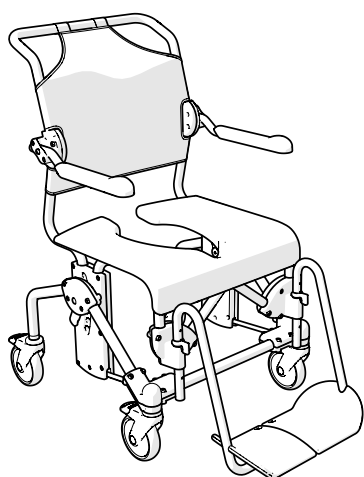Swift Mobil-2

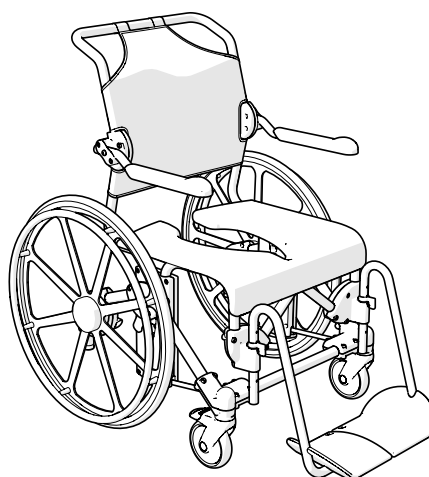Swift Mobil 24"-2

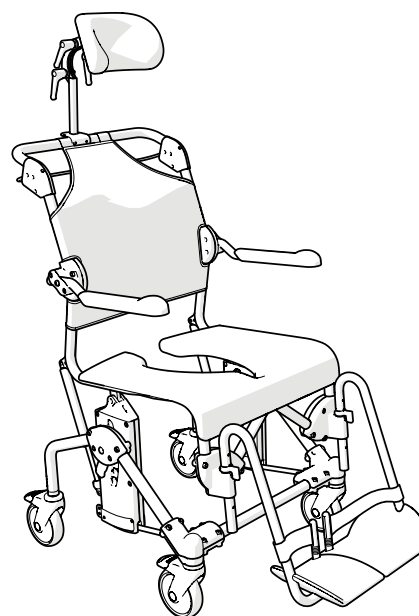Swift Mobil Tilt-2

Swift Mobil-2 XL

Swift Mobil Tilt-2 XL

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Swift Mobil-2	Swift Mobil-2 XL	Swift Mobil 24"-2	Swift Mobil Tilt-2	Swift Mobil Tilt-2 XL
●	●	●	●	●
○	○	○	●	●
●	●	●	●	●
●	●	●	●	●

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English

General

Thank you for choosing an Etac product. In order to avoid damage during assembly, handling and use, it is important to read this manual and save it for future reference. You can also find it at www.etac.com. You can select your language via the "International" and "Local websites" link. Here you will also find other product documentation, such as prescriber information, pre purchase guide, and reconditioning instructions.

In the manual, the user is the person sitting in the chair. The carer is the person helping the user.

Intended use:

The Swift Mobile 2 family is a range of mobile shower chairs and toilet chairs (hereafter called the product) intended for indoor use when performing hygiene tasks for people with restricted mobility. The product can be used in the shower, at a sink, or over the toilet, as well as for transfers to and from the hygiene room.

The product is intended for use in home environments and institutions.

The product is suitable for use in wet areas, but not pools or similar corrosive environments.

The product is intended to compensate for various types of injuries or disabilities that limit the mobility of the user.

The product is intended for adult users, generally persons who are at least 146 cm tall and who weigh at least 40 kg.

The product is not intended for any other use or environment beyond those stated above.

The product has been tested and complies with the requirements set out in EN ISO 10993-1 and ISO 17966.

The product conforms to the requirements of the Medical Devices Regulation (EU) 2017/745.

Contra indications:

There are no known contraindications.

This symbol appears in the manual along with text. It indicates instances where the safety of the user or the assistant may be at risk.

The images which are referred to on each page are found at the start of the manual (pages 3 to 11). See Figures A to J for the corresponding text.

This product should be recycled as per national regulations.

Item list: ● = Applies to, ○ = Accessories

See the item list on the right which displays the information applicable to each product.

Product description

See Figure A

A1. Bolt, **A2.** Leg support, **A3.** Seat, **A4.** Arm support, **A5.** Back support upholstery, **A6.** Item no label (serial no, production date)

A7. Head support

*The product's date of manufacture can be read from the barcode on the product. The number 11 is shown below the barcode in brackets. The number combination after these brackets is the date of manufacture.

Assembly

See Figure B

Technical data

See Figure C

Warranty / maintenance

See Figure D

5-year guarantee against material and manufacturing defects. For terms and conditions, see www.etac.com.

Service life 7 years. For complete information regarding the service life of the product, see www.etac.com.

D1. Clean the product with normal non-abrasive cleaning agents (pH value of between 5-9). Rinse and dry.
Disinfect the product with an alcohol-based disinfectant, 70%.
The product is otherwise maintenance-free.

D2. The product and its accessories can be safely cleaned at 85°C for 3 minutes (unless otherwise stated). This applies for reconditioning.

The product can be reconditioned. Refer to the reconditioning instructions at www.etac.com.

Constituent materials are resistant to common disinfectants.

The following surface treatment methods have been used to protect against corrosion:

Lacquered surfaces = Polyester powder coating

Non-lacquered steel surfaces = Fe/Zn-Ni

Reuse

The product is suitable for reuse. The number of times it can be reused depends on how often and in which way the product is used. Before reuse, thoroughly clean and disinfect the product and check its safety features. If necessary, take remedial action.

Storage

The product should be stored indoors in a dry place at a temperature above 5°C. If the product has been stored for a long time (more than four months), its function must be checked by an expert before use.

Swift Mobil-2	Swift Mobil-2 XL	Swift Mobil 24"-2	Swift Mobil Tilt-2	Swift Mobil Tilt-2 XL
●	●	●	●	●
○	○	○	●	●
●	●	●	●	●
●	●	●	●	●
●	●	●	●	●

See Figure E

- E1.** Make sure the castors are locked when users sit down on or get up from the chair.
 - An anti-tip attachment must be used if the user is an amputee.
- E2.** Never stand on the foot supports as there is a risk of tipping!
- E3.** The sitting position has an impact on stability.
- E4.** Chair stability is improved if the castors point outward.
- E5.** See the maximum user weight for your commode chair in the table (page 5).

E7. Do not use defective products.

Transfers

See Figure F

F4. Out-turned castors improve the chair's stability.

F5-F7. Fold up the arm support on the transfer side of the chair. Carry out the transfer.

See Figure G

G1-G7. The height is easily adjusted by inserting the bolt in the required position. See distances in the technical data section.

G8. The height of the arm supports and/or width between the arm supports can be adjusted. Loosen the screws, twist the support to the desired position and retighten the screws. Then loosen and refit the arm support plate in the correct position. Swift Mobil Tilt-2 XL is delivered with the arm supports mounted in the most inward position (seat width 48 cm).

G10. Fit the gas piston in position B to tilt the seat at angles from 0° to 35°

G11. When adjusting the brake, the brake pin should be set 18 to 22 mm from the rear wheel for good grip.

See Figure H

H2. Various height adjustment options are available by moving the locking pin. The leg supports have 6 different height settings.

See Figure 1

11. The assistant can adjust the sitting angle with ease by holding in the handle and tilting to the required position. Make sure that the wheels are locked if the tilt function is used when the user is sitting in the chair.

[illegible]

English

Accessories

See Figure J

In the section "Accessories", you can see all accessory options available for your product. Variations may occur in different markets. There are restrictions on accessory combinations, for more information, see www.etac.com.

See Figure J	Sw M.-2	Sw M.-2 XL	Sw M. 24"-2	Sw M. Tilt-2	Sw M. Tilt-2 XL	See Figure J	Sw M.-2	Sw M.-2 XL	Sw M. 24"-2	Sw M. Tilt-2	Sw M. Tilt-2 XL
J1. Leg and calf rest, angle-variable	☐	☐	☐	☐	☐	J23. Armrest cushion, pair	☐	☐	☐	☐	☐
J2. Footrest	☐	☐	☐	☐	☐	J24. Trunk support	☐	☐	☐	☐	☐
J3. Angle adj. footrest	☐	☐	☐	☐	☐	J25. Chest and hip positioning belts	☐	☐	☐	☐	☐
J4. Amputee legrest	☐	☐	☐	☐	☐	J26. Chest belt	☐	☐	☐	☐	☐
J5. Commode pan	☐	☐	☐	☐	☐	J27. Positioning belt	☐	☐	☐	☐	☐
J6. Pan-holder attachments	☐	☐	☐	☐	☐	J28. Crossbar 605 mm	☐	☐	☐	☐	☐
J7. Attachments for low pan-holder mounting	☐	☐	☐	☐	☐	J28. Crossbar 665 mm	☐	☐	☐	☐	☐
J8. Commode pan holders	☐	☐	☐	☐	☐	J29. Anti-tip stabilisers including bracket	☐	☐	☐	☐	☐
J9. Bucket holders	☐	☐	☐	☐	☐	J30. Anti-tip stabilisers excluding bracket	☐	☐	☐	☐	☐
J10. Pan	☐	☐	☐	☐	☐	J31. Tiliter	☐	☐	☐	☐	☐
J11. Lid with handle	☐	☐	☐	☐	☐	J32. Rear wheel kit	☐	☐	☒	☐	☐
J12. Splash guard	☐	☐	☐	☐	☐	J33. Headrest and bracket	☐	☐	☐	☒	☒
J13. Soft Comfort seat, 2 cm high, open at back	☐	☐	☐	☐	☐	J34. Heel straps	☐	☐	☐	☒	☒
J14. Soft Comfort seat, 4 cm high, open at back	☐	☐	☐	☐	☐	J35. Hemiplegic armrest	☐	☐	☐	☐	☐
J15. Comfort seat	☐	☐	☐	☐	☐	J36. Armrest lock	☐	☐	☐	☐	☐
J16. Soft Comfort seat	☐	☐	☐	☐	☐	J37. Extended back support	☐	☐	☐	☐	☐
J17. Soft Comfort seat - narrow opening	☐	☐	☐	☐	☐	J38. Short leg support	☐	☐	☐	☐	☐
J18. Soft cushion - standard opening	☐	☐	☐	☐	☐	J39. Widening kit for leg support	☐	☐	☐	☐	☐
J19. Soft cushion - oval opening	☐	☐	☐	☐	☐	J40. Anti tip front	☐	☐	☐	☐	☐
J20. Seat cover	☐	☐	☐	☐	☐	J41. Back Swift Mobil-2, upright	☐	☐	☐	☐	☐
J21. Soft back pad	☐	☐	☐	☐	☐		☐	☐	☐	☐	☐
J22. Comfort padding	☐	☐	☐	☐	☐		☐	☐	☐	☐	☐

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Etac Supply Center AB