

Instruction for use



IMM1064_IFU.Version7.2020-09-17



RescueSheet

immedia
by Etac




Index

Illustrations.....	3
User manual – English	4
Bruksanvisning – Svenska.....	9
Bruksanvisning – Norsk	14
Brugervejledning – Dansk.....	19
Käyttöohjeet – Suomi.....	24
Gebrauchsanweisung – Deutsch ...	29
Gebruiksaanwijzing – Nederlands ...	34
Instruction d'utilisation – Français...	39
Istruzioni per l'uso – Italiano	44
Instrucciones de uso – Español....	49


See www.etac.com



Illustrations

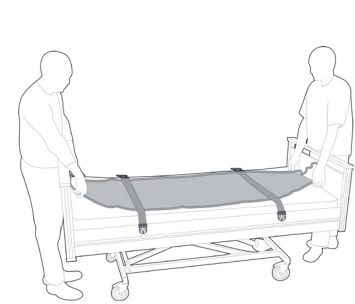


1

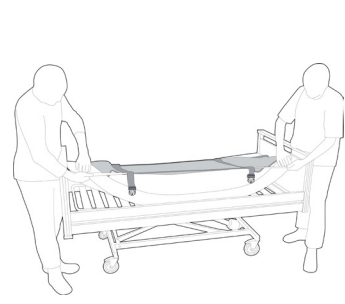


2



3

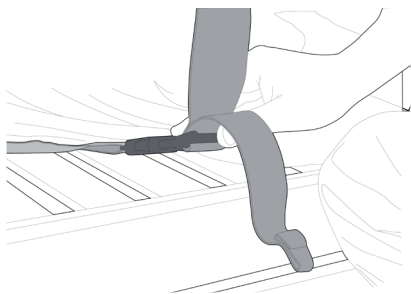

4

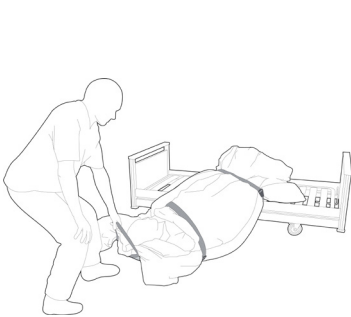

5



6

ILL.

GB

SV

NO

DA

FI

DE

NL

FR

IT

ES



General

Thank you for choosing an Etac product.

To avoid accidents and injury when moving and handling the products, please read this manual carefully.

The person described as “the user” in this “Instruction for Use” is the person who is lying or sitting on the product. The carers are the people who manoeuvre the product.



This symbol appears alongside the text in the manual. It draws the reader’s attention to points at which there may be a risk to the health and safety of the user or carer.

These products comply with the standards applicable for Class 1 products in the Medical Device Regulation (EU) 2017/745 on medical devices.

At Etac we strive to continually improve our products and therefore we reserve the right to make changes to products without prior warning. All measurements given on illustrations and similar material are for guidance only and Etac cannot be held liable for errors and defects.

The information given in this manual, including recommendations, combinations and sizing, does not apply to special orders and modifications. If the customer makes adjustments, repairs or uses combinations not predetermined by Etac, the Etac CE certification and Etac warranty will not be valid. If in doubt, please contact Etac.

Warranty: Two-year warranty on material and manufacturing defects, provided that the product is used correctly.

Expected lifetime: The device has an expected service life of 1 to 2 years under normal use. The service life of the device varies depending on usage frequency, loads and how often and how it is washed.

For further information on the Etac Patient handling, see www.etac.com.

In case of an adverse event occurred in relation to the device, incidents should be reported to your local dealer and the national competent authority in a timely manner. The local dealer will forward information to manufacturer.

Symbols in general

	Washing		Disinfect / Wipe off		Warning		Manufacturer
	Do not bleach		Never leave on the floor		Production year, month and day		CE -certification
	Tumbledry		Risk of sliding down		Batch no./Lot no.		Keep dry
	Iron		Keep out of the sun		User mass limit = maximum rated load		Non sterile
	Do not dryclean		Read the user manual		The product can be scrapped according to national regulations		Article number

MD Medical device

Intended use

RescueSheet can be placed permanently under the mattress to ensure quick and effective evacuation in case of fire.

We recommend frequent practice in the use of RescueSheet for all users.

Intended environment

Acute care, Long-term care, Home care



Practical handling

Positioning

RescueSheet fits permanently under the mattress, with the red straps down and fastened to the corners with elastic cords.

The loose ends of the red straps must be placed on the side not touching the wall.

Use

In an evacuation emergency, pull out the red straps and lay over the bed. Then fasten these straps to the straps on the other side, and tighten so that the mattress bends upwards on each side.

The user can now be transferred out of the room, and out to safety. The mattress and bed covering protect the user against blows, jolts and the cold.

Use the handle at the top of the bed and pull the user and mattress gently off the bed frame. Hold the handle high, giving the user maximum protection when the mattress touches the floor.

ILL.

GB

SV

NO

DA

FI

DE

NL

FR

IT

ES



Notice

Special features

Always check the product before use and after washing.

Never use a defective product. If the product shows signs of wear-and-tear, it must be scrapped.

Read these instructions carefully.

It is important that carers receive instruction in manual handling. Etac offers advice and training for carers. For further information, contact Etac.

Always use the correct manual handling techniques.

Encourage the user to assist where possible.

To ensure that the user feels safe and that every manual handling is smooth, always plan the manual handling in advance.



Never leave the product on the floor.



Always conduct a risk assessment, and ensure that the assistive product can be used for the individual user and, in combination with other devices, that it is safe for the user and carers.



During fire practice it is recommended to test the buckles under load. If the buckles are damaged or fail to function, replace the product.

If in any doubt - please contact Etac.

Material and cleaning

Material:

Fabric: Nylon
Straps: Polyester
Buckles: Acetal.

The product should be checked regularly, preferably each time it is used, and especially after washing.

Check that there is no damage to seams or fabric.



Never try to repair a product yourself.
Never use a defective product.

Washing instructions:

Do not use fabric softener.
This product can be cleaned with a solvent-free cleaning agent with a pH level between 5 and 9, or a 70% disinfection solution.



Products and accessories

Article no.	Product Description	Size (mm)
IM108088	RescueSheet	W750xL2000
IM108089	RescueSheet - extra wide	W850xL2000

ILL.

GB

SV

NO

DA

FI

DE

NL

FR

IT

ES



Service information

Inspection

First inspection:

- Is the packaging intact?
- Read the label on the packaging and check the article no. and product description
- Check that Short Instruction is enclosed – updated User Manual/Instruction for use can be downloaded from www.etac.com or contact your Etac Customer Service or the local dealer.
- Check the label on the product - does it include article no., product description, lot/batch no., cleaning instructions and supplier name?

Periodic Inspection:

- Make sure that your Manual transfer system always is in a perfect condition.
- After wash – always check materials, stitching, handles, buckles.
- If the product shows signs of wear and tear, it must be removed from service immediately.



Visual inspection/Check the product:
Ensure that the materials, seams, stitching, handles, buckles are intact/faultless.



Stability test:
Try to bend the product verify that it feels firm, stable and solid.



Mechanical load or stress/Test handles:
Draw hard in the handles in opposite directions and control the material, seams and stitching.
Test buckles:
Lock the buckle and draw in opposite directions, control the buckle and the stitching.



Test low/high friction:
Place the product on a firm surface or a bed, place your hands on the material and test the friction by pushing your hands down into the product.
Low friction – slides effortlessly
High friction – no sliding/moving

Check:

Manual transfer aids with handles, straps and buckles:



Manual transfer aids with for sitting and lying transfer and with high or low friction:



Production year, month and day



The product can be scrapped in accordance with national regulations



Read the user manual

