

LINAK Bluetooth® Adapter for bed



- Made for
- iOS devices with Bluetooth® 4.0 or newer
 - iOS devices with iOS 9.0 or newer
 - Android devices with Bluetooth® 4.0 or newer
 - Android devices with 5.0 or newer



Made for
 iPod  iPhone  iPad

iPod touch, iPhone, and iPad are trademarks of Apple Inc., registered in the U.S. and other countries.

“Made for iPod touch”, “Made for iPhone,” and “Made for iPad” mean that an electronic accessory has been designed to connect specifically to iPod touch, iPhone, or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPod touch, iPhone or iPad may affect wireless performance.


The Bluetooth word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Ascion, LLC is under license. Other trademarks and trade names are those of their respective owners.

Table of contents All manuals and user guides available at www.linak.com

The Adapter	3
Connect to bed	4
How to download the App for Android.....	4
How to download the App for iOS	4
How to put the adapter in pairing mode	5
Pairing Bluetooth	6
Buttons.....	7

The Adapter

The Bluetooth® Adapter allows you to adjust your LINAK product wirelessly from an app or Bluetooth® hand control. Simply plug the Bluetooth® Adapter in your HOMELINE® control box or HOMELINE product and download the app: Your application can now be adjusted via your smartphone or tablet. The Bluetooth® Adapter is a very compact unit requiring a minimum of space. No mounting is needed, just plug in the adapter and you are ready.

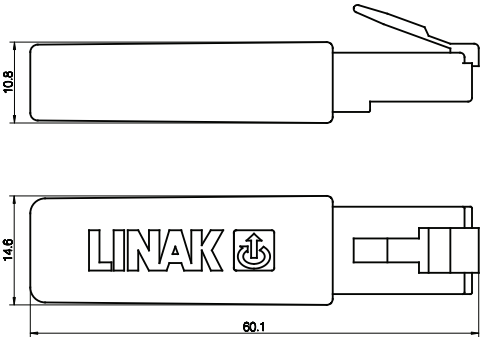

Connect to bed

Connect the Bluetooth® Adapter to the CBH Advanced or the TD3 simply by plugging it in. The Bluetooth® Adapter is a very compact unit requiring a minimum of space. No mounting is needed, just plug in the adapter and you are ready.

How to download the App for Android:

- Step 1:**
Please use your Android device to log into Google Play Store.
- Step 2:**
Type “Bed Control” in the search column.
- Step 3:**
Select the “Bed Control” app best suited for your product. If your product supports memory positions, select the “Bed Control Memory” app.
- Step 4:**
Click install to download the app to your device.


How to download the App for iOS:

- Step 1:**
Please use your Apple device to log into App store.
- Step 2:**
Type “Bed Control” in the search column.
- Step 3:**
Select the “Bed Control” app best suited for your product. If your product supports memory positions, select the “Bed Control Memory” app.
- Step 4:**
Click install to download the app to your device.




How to put the adapter in pairing mode

Plug the Bluetooth Adapter into a powered HOMELINE® control box or a TD3.

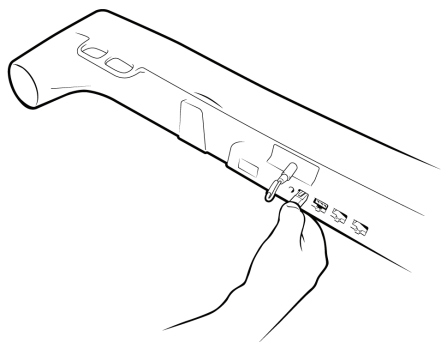

The pairing mode will last 3 minutes and is indicated with a blue flashing LED light that is visible through the transparent plastic. (It can take up to 10 seconds before pairing mode starts).

Download the software you want in App Store or Google Play store and follow the instructions in the programs.

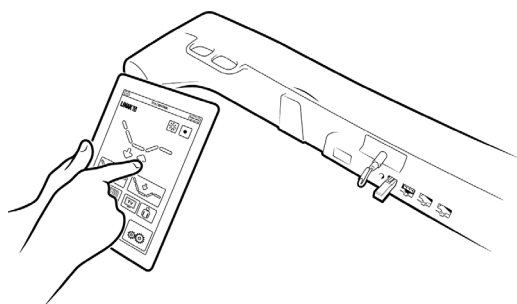

Pairing Bluetooth

All manuals and user guides at all-guides.com

When the Bluetooth® Adapter is plugged into a powered TD3 or CBH Advanced, it will be in pairing mode for 3 minutes. A blue blinking LED will indicate this.

You pair your device to the Bluetooth® Adapter in the app. Enable Bluetooth on your smartphone or tablet. Open the Bed Control app or Bed Control Memory app and a "search bed" page will be available.

When searching; All devices will be called BLE2LIN followed by a 4-digit number.

A list of Bluetooth Adapters and previously connected Bluetooth Adapters will appear, tap the adapter you want to connect to. Your device is now paired to the adapter and the app will display the operation page.

The Bluetooth Adapter automatically reestablish the connection to the paired devices once they are within range and turned on.

The only difference between Bed Control and Bed Control Memory is that there is no memory positions displayed in Bed Control.

Buttons

Select back or leg

Select which part to operate back or leg. Backrest (A) and leg rest (B)

UP (C)

Moves the selected part up

Down (D)

Moves the selected part down

Synchronous drive up and down (E), (F)

Moves both back and leg up or down at the same time – independent of the selected part

Light (G), (H)

Turns on or off the under bed light. The light box must be purchased as an accessory.

Storing positions (I)

4 stored positions can be directly accessed by pressing the corresponding button

Storing positions

Via the settings page 4 positions are possible to store. Simply press the button and the actual position is stored. Via the bed control page the bed moves to this position when the corresponding icon is pressed.

Switch bed

To select another connected bed, press the Bluetooth icon in the top right corner of the screen. Then select a connected bed from the list or scan to find new devices.

Renaming bed

The paired application can be renamed in the Android app on the settings screen. The iOS version does not support renaming the device, so the new name will only be saved in the app.

